INSTALLATION GUIDE

RephiDuo S Pack 1（EDI/RO feed）

RephiLe Catalogue Number:
RRSYNM0IX
Corresponding Millipore Part Numbers:
SIMPAK0R1


Suitable Millipore Systems:
Simplicity, Simplicity 185
Synergy, Synergy 185

This product is comparable to the corresponding Millipore cartridge. Please refer to Millipore's User Manual for a detailed instruction on installation.

A brief description of the installation process is described below:

1. Press "OPERATE/STANDBY" for 2 seconds to let the system into **Standby** mode.
2. Remove the final filter from the water dispense outlet.
3. Put a container underneath the dispense outlet, open the dispense valve for a few seconds to empty the system.
4. Locate the cartridge Lock button (shown in picture), push down the lock button, lift the cartridge security bar to the left.
5. Remove the exhausted cartridge from the system.
6. Open the new cartridge package, remove the two blue protective covers from the new cartridge, then wet the O-rings with pure water.
7. Snap-in the new cartridge by pushing it into its position, close the security bar and lock the cartridge to its position.
8. Degas the cartridge: Remove the final filter. Connect the degas connector to the water outlet and the sink. If cannot locate the degas connector, put a large container underneath the water dispense outlet.
9. Press "OPERATE/STANDBY" for 2 seconds to get the system into Operate mode. Open the water dispense valve to dispense 1-2L of water, then shut off the water dispense valve.
10. Remove the degas connector or beaker, install a final filter, then degas the final filter (see installation guide for the final filter).

The system is ready for use.

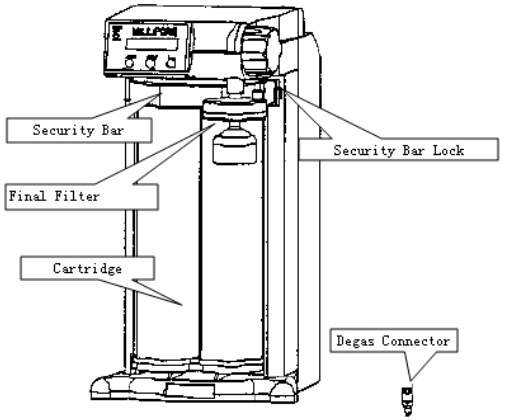

RephiSolo, RephiDuo and RephiQuatro are trademarks of RephiLe Bioscience, Ltd.
© 2010. RephiLe Bioscience, Ltd. All rights reserved.
All other trademarks are trademarks of their respective manufacturers.