

avantor[™]
delivered by **VWR**[™]

VWR.COM

BLICKMAN



IV Stands

Optimal continuous fluid delivery

Includes:

- IV Stands
- IV Stand Accessories
- IV Irrigation Safety Systems


THE BLICKMAN DIFFERENCE

IV Stands

- Our IV poles are optimal for maintaining continuous fluid delivery for every patient, available in a wide range of adjustable heights with a variety of easy locking options that can be operated with one hand.
- Our IV poles are all available in high quality chrome and stainless steel, and set on 4, 5, or 6 leg bases with casters, no projections, and great balance/tip resistance. Some bases are gravity weighted to help compensate for heavier loads. All bases have the ability to spin 360° in tight spaces.
- Most of our models have 2-, 4-, 6- and 8-(rake) hook adapters available. All come with colored secure grip tips.

Key Product Features

*See page 9 for more details on all optional features

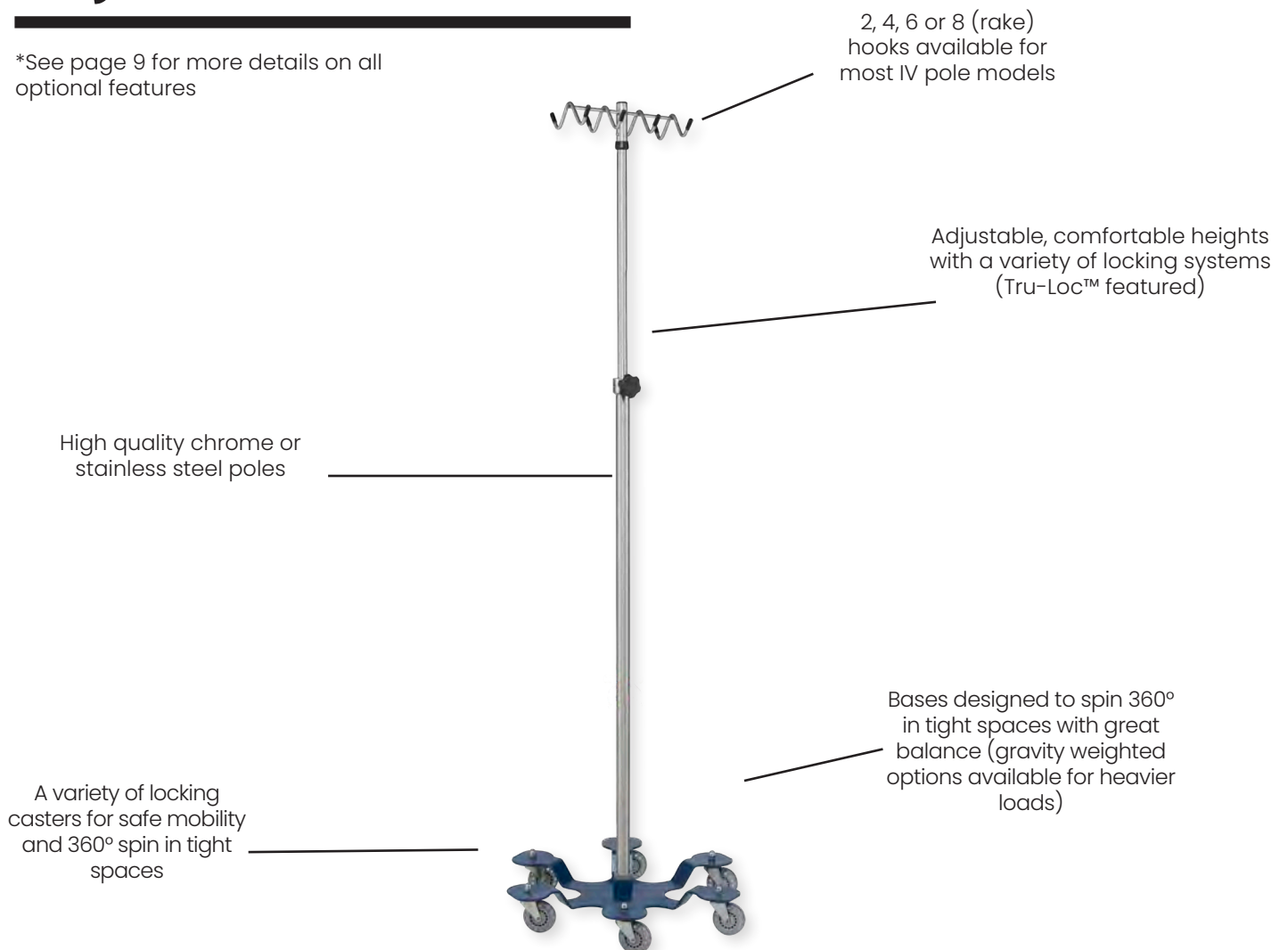


IV Stands

4 Leg

4 LEG CHROME IV STANDS

Description: We offer two 4-leg chrome IV stand models, one with a standard base, and one with a heavyweight base. Both base types include Wall Saver™ tips to protect against wall scuff marks. They both move on 2" plastic twin wheel hooded swivel casters and they both come standard with two hooks. A four hook removable adapter is available for model 1305.

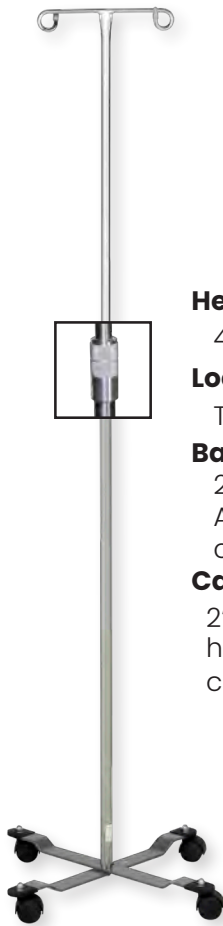

Height adjustment:

47 1/2" - 85 1/2"

Lock:

Twist

Base:

22" diameter
Apoxy coated
chrome

Casters

2" plastic twin wheel
hooded swivel
casters



Height adjustment

51 1/2" - 93"

Lock:

Tru-loc™

Base:

22" diameter
Chrome plated
heavy bar steel

Casters

2" plastic twin wheel
hooded swivel
casters


Name	Standard Base	Standard Base	Heavyweight Base
# of Hooks	2 Welded Hooks	4 Removable Hooks	2 Removable Hooks
Model #	1305	1305-4P	1310
VWR Catalog #	76545-294	76545-296	76545-298
Shipping Weight	13 lbs.	13 lbs.	14 lbs.

IV Stands

4 Leg



4 LEG STAINLESS STEEL IV STANDS

Description: We offer two types of 4-leg stainless steel IV stands, one with a standard base, and one with a heavyweight base. They come standard with either 2 or 4 removable hooks. Contact customer service for 6- or 8-(rake) hooks.

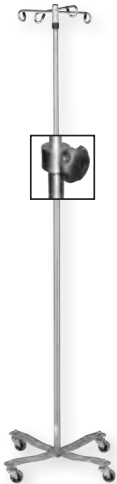

Name	Standard Base
# of Hooks	2 Hooks
Model #	1410SS
VWR Catalog #	76545-316
Shipping Weight	13 lbs.
# of Hooks	4 Hooks
Model #	1410SS-4
VWR Catalog #	76545-318
Shipping Weight	14 lbs.

Height Adjustment: 51 3/4" - 94"
 Mechanism: Tru-Loc™
 Base Material: Stainless steel with Wall Savers™, 22" diameter
 Casters: 2" hard rubber swivel casters

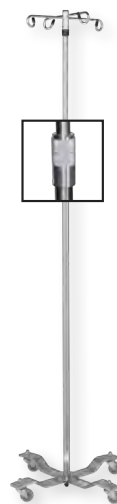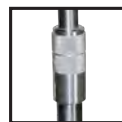

Name	Heavyweight Base	
Lock Type	Twist Lock	Thumb Control
# of Hooks	2 Hooks	
Model #	8890SS	8889SS
IVWR Catalog #	76545-256	76545-252
Shipping Weight	24 lbs.*	24 lbs.*
# of Hooks	4 Hooks	
Model #	8890SS-4	8889SS-4
VWR Catalog #	76545-258	76545-254
Shipping Weight	24 lbs.*	24 lbs.*

Height Adjustment: 66" - 100"
 Mechanism: Twist Lock, Thumb Control, Foot Operated
 Base Material: Heavyweight stainless steel flat bar, 22" diameter
 Casters: 2" roller blade swivel casters

IV Stands

5 Leg

5 LEG IV STANDS

Description: Our 5-leg IV stands come in either chrome (standard & heavy duty) or stainless steel poles and stand on a 24" diameter, polished die cast aluminum base with Wall Saver™ caps.

Chrome Pole

Standard (1")


# of Hooks	2 Hooks	4 Hooks
Caster Heights	2 1/2" Plastic	2 1/2" Plastic
Height Adjustments	52 1/2" - 93 1/2"	52 1/2" - 93 1/2"
Model #	1315-2	1315-4
VWR Catalog #	76545-300	76545-302
Shipping Weight	16 lbs.*	16 lbs.*

Heavy Duty (1 1/4")

# of Hooks	2 Hooks	4 Hooks	6 Hooks
Caster Heights	3" Rubber	3" Rubber	3" Rubber
Height Adjustments	56" - 100"	56" - 100"	56" - 100"
Model #	1350-2	1350-4	1350-6
VWR Catalog #	76545-304	76545-306	76545-308
Shipping Weight	17 lbs.*	17 lbs.*	17 lbs.*



Height Adjustment Mechanism: Tru-Loc™. Base Material: 5 Leg, Polished die cast aluminum, 24" diameter. Casters: 3" hard rubber OR 2 1/2" twin wheel swivel casters

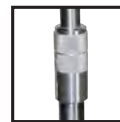

Stainless Steel Pole

Lock Type

Tru-Loc™

Twist Lock

Thumb Control

# of Hooks	2 Hooks	4 Hooks	6 Hooks
Caster Heights	2" Rubber	3" Rubber	3" Rubber
Height Adjustments	52" - 94"	74" - 110"	67" - 98"
Model #	1415SS	7792SS	7794SS
VWR Catalog #	76545-260	76545-264	76545-248
Shipping Weight	16 lbs.*	15 lbs.*	16 lbs.*

of Hooks

4 Hooks

Caster Heights	2" Rubber	3" Rubber	3" Rubber
Height Adjustments	52" - 94"	74" - 110"	67" - 98"
Model #	1415SS-4	7792SS-4	7794SS-4
VWR Catalog #	76545-262	76545-270	76545-250
Shipping Weight	16 lbs.*	16 lbs.*	17 lbs.*



Height Adjustment Mechanism: Tru-Loc™, Twist Lock, Thumb Control Base Material: 5 Leg, Polished die cast aluminum, 24" diameter. Casters: 2" or 3" hard rubber swivel casters


*Shipping weight with an asterisk is the combined shipping weight of separate parts

IV Stands

5 Leg

5 LEG WASHABLE IV STAND

Description: Our Hackensack 5-leg IV stand is cart washable as it was built specifically to withstand high temperatures and cleaning fluids with weep holes in the pole and in the base for fluid drainage. This model comes with a stainless steel pole with a twist lock height adjustment mechanism. Choose from 2-, 4-, and 8-(rake) removable hooks.

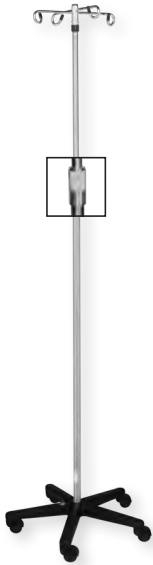

Features Include:

- ✓ 22 7/8" diameter, 5 leg black composite base
- ✓ 2 1/2" twin wheel hooded swivel casters
- ✓ 1" outside diameter base pole, and 3/4" outside diameter extension pole
- ✓ Height extends from 74" to 110" off the ground



Name	Washable IV Stand
Lock Type	Twist Lock
# of Hooks	2 Hooks
Model #	7795SS-2
VWR Catalog #	76545-272
Shipping Weight	16 lbs.*
# of Hooks	4 Hooks
Model #	7795SS-4
VWR Catalog #	76545-274
Shipping Weight	16 lbs.*
# of Hooks	8 Hooks
Model #	7795SS-RH
VWR Catalog #	76545-276
Shipping Weight	16 lbs.*



*Shipping weight with an asterisk is the combined shipping weight of separate parts

IV Stands

6 Leg

6 LEG IV STANDS



Description: Our 6-leg IV stands are designed for hospital applications requiring heavier loads by using low center of gravity bases. They come in either a chrome or a stainless steel pole with a Tru-Loc™ height adjustment mechanism that allows the pole to reach a minimum of 52" high off the ground and a maximum of 90" off the ground. Our 6 leg IV stands come standard with 4-, 6-, and 8-(rake) removable hooks.

Features Include:

- ✓ 23 3/8" diameter base with low center of gravity for improved stability
- ✓ Heavy duty 3" roller blade swivel casters for easy maneuverability
- ✓ Steel powder coated finished base for durability (many colors available)
- ✓ Accommodates up to 4 infusion pumps and up to 8 IV solution bags while preventing tubing line entanglement



Powder Coated Base

Stainless Steel Base

Name	Chrome Pole	SS Pole	SS Pole
# of Hooks	4 Hooks		
Model #	1370-4	1370SS-4	1371SS-4
VWR Catalog #	76545-310	76545-284	76545-288
Shipping Weight	29 lbs.*	29 lbs.*	26 lbs.*
# of Hooks	6 Hooks		
Model #	1370-6	1370SS-6	1371SS-6
VWR Catalog #	76545-312	76545-286	76545-290
Shipping Weight	29 lbs.*	29 lbs.*	26 lbs.*
# of Hooks	8 (rake) Hooks		
Model #	1370-RH	1370SS-RH	1371SS-RH
VWR Catalog #	76545-314	76545-282	76545-292
Shipping Weight	29 lbs.*	29 lbs.*	26 lbs.*



Shipping weight with an asterisk is the combined shipping weight of separate parts


MR Conditional

IV Stands

MR IV STANDS



Description: Our 5 leg MR conditional IV stand consists of a stainless steel pole on a 24" diameter, polished diecast aluminum base with Wall Saver™ caps and 2" MR casters. It is height adjustable from 74" to 110" via a twist lock adjustment mechanism. Choose from either 2- or 4-hooks or contact customer service for 6 or 8 (rake) hooks.

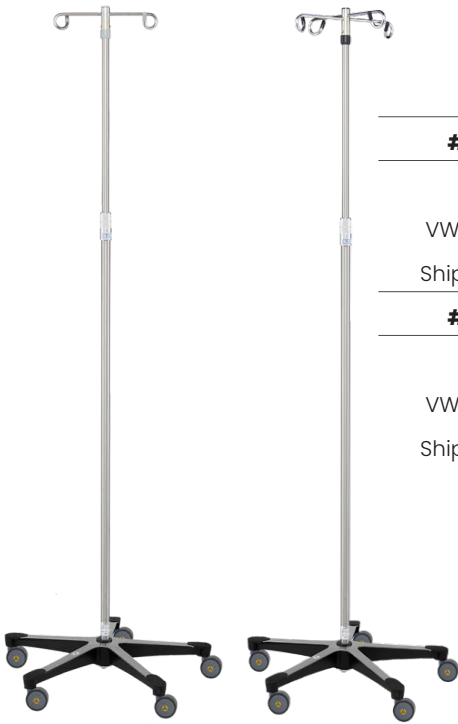

Name	MR Conditional IV Stands
# of Hooks	2 Hooks
Model #	7792MR-2
VWR Catalog #	76545-266
Shipping Weight	16 lbs.
# of Hooks	4 Hooks
Model #	7792MR-4
VWR Catalog #	76545-268
Shipping Weight	17 lbs.



IV Stands

Construction Specifications

Construction Specifications

Universal Dimensions	Height Range: 47" - 110" Base Diameter Range: 22" - 24"
Materials	304 stainless steel Chrome
Gauge	Stainless Steel IV Stands: Base Pole: 1" Diameter, 16GA; Extension Pole: ¾" Diameter, 18GA Exception: 1370/1371 series Base Pole: 1 ¼" Diameter, 16 GA; Extension Pole: 1" Diameter, 18GA
Welding	All welding by heliarc method and exposed welds are ground smooth
Finish	Exposed surfaces are free from pit marks, weld seams, scale and are polished with #4 polish to a satin finish

Shipping Information

All shipping weights are estimates and may be changed without prior notice. Specifications may be improved and dimensions changed without prior notice. Please contact customer service for estimated shipping costs at info@blickman.com.

IV Stand parts shipped separately to avoid damage and to lower costs.

Some assembly required.

THE BLICKMAN DIFFERENCE

IV Irrigation Safety Systems

- Our IV Irrigation Safety Systems (IVTs) offer superior safety, efficiency, and financial value for fluid intensive IV irrigation procedures, such as arthroscopic, hysteroscopic, endoscopic, TURP, etc.
- They maintain continuous fluid delivery by offering infinitely adjustable heights for up to eight 3 liter or 5 liter IV bags on four independent channels.
- Their easy, one handed operation allows nurses to attach IV bags at comfortable heights, starting at 5 ft., and then raise bags up to 8 ft. (or to 9 ft. with custom extender) simply by pushing down on "squeeze/release" handles, an ergonomic solution that meets OSHA standards and greatly reduces the threat of back injuries.

Key Product Features

Up to eight 3 liter or 5 liter IV bags on four independent channels


High quality, easy maintenance, light weight aluminum

Locking casters for safe mobility and 360° spin in tight spaces

Push down "Squeeze/Release" handles to raise IV bags

Indefinite adjustable heights from 5 ft. to 8 ft. (9 ft. with custom extender)

26 1/2" diameter, inverted, 5 legged base resistant to tipping

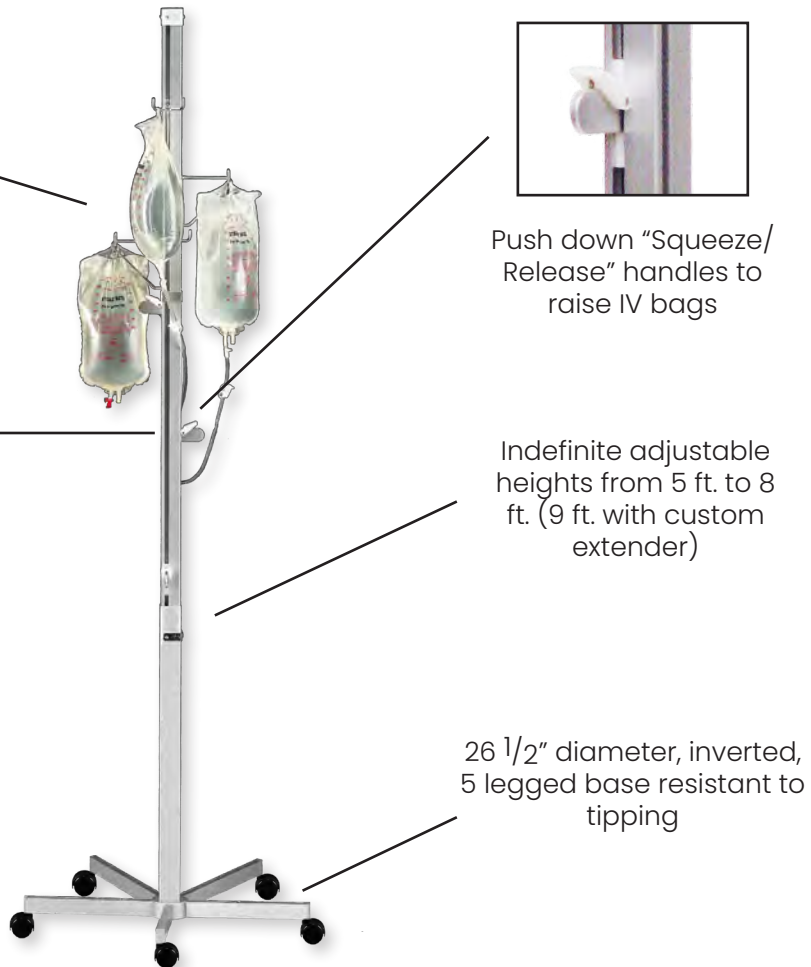

WHY CHOOSE OUR IVT?



Safe:

- This ergonomic solution meets OSHA standards and greatly reduces the threat of back injuries with its one-handed operation which allows nurses to attach IV bags at comfortable heights, simply by pushing down on “squeeze/release” handles, an
- Hard anodized aluminum construction with stain resistant finishes makes it so that no surfaces can chip, flake, or rust to introduce particle matter into a sterile or surgical environment.
- The lowest positions are cushioned by shock absorbers, thus sparing the common banging known to occur with the standard IV poles, which results in tearing of bag tabs and spilling remaining fluids onto the O.R. floor.
- The steady, five-legged bases resist potential tipping hazards, eliminating the need for “sand bagging”.

Efficient:

- There are no electrical aspects that require overnight charging, so they are always ready, whenever you are ready.
- The ability to stagger heights of up to eight 3 liter or 5 liter bags ensures continuous flow of liquids, as it allows you to change an IV bag, while the other bags are still irrigating.
- Easy to clean and maintain
- Compact footprint frees up valuable space on the O.R. floor.

Affordable:

- Our IV Irrigation systems are value leaders among a few comparable models, saving you up to **\$50/procedure** on IV tubes vs. those required with far more expensive, pump driven systems.

Irrigation Systems


IVTs

Description: We offer two types of IVTs, one with 2 hooks, which holds up to 4 IV bags and one with 4 hooks, which holds up to 8 IV bags. Hooks can raise and lower from 67" to 81". They travel on black plastic twin wheel hooded swivel casters. Both IVTs come with 26 1/2" diameter inverted, 5 leg base.

Our IVT stands have the same dimensions:

Width: 12 3/4" x Height: 89" x Depth: 12 3/4"

Note: these dimensions do not include the base

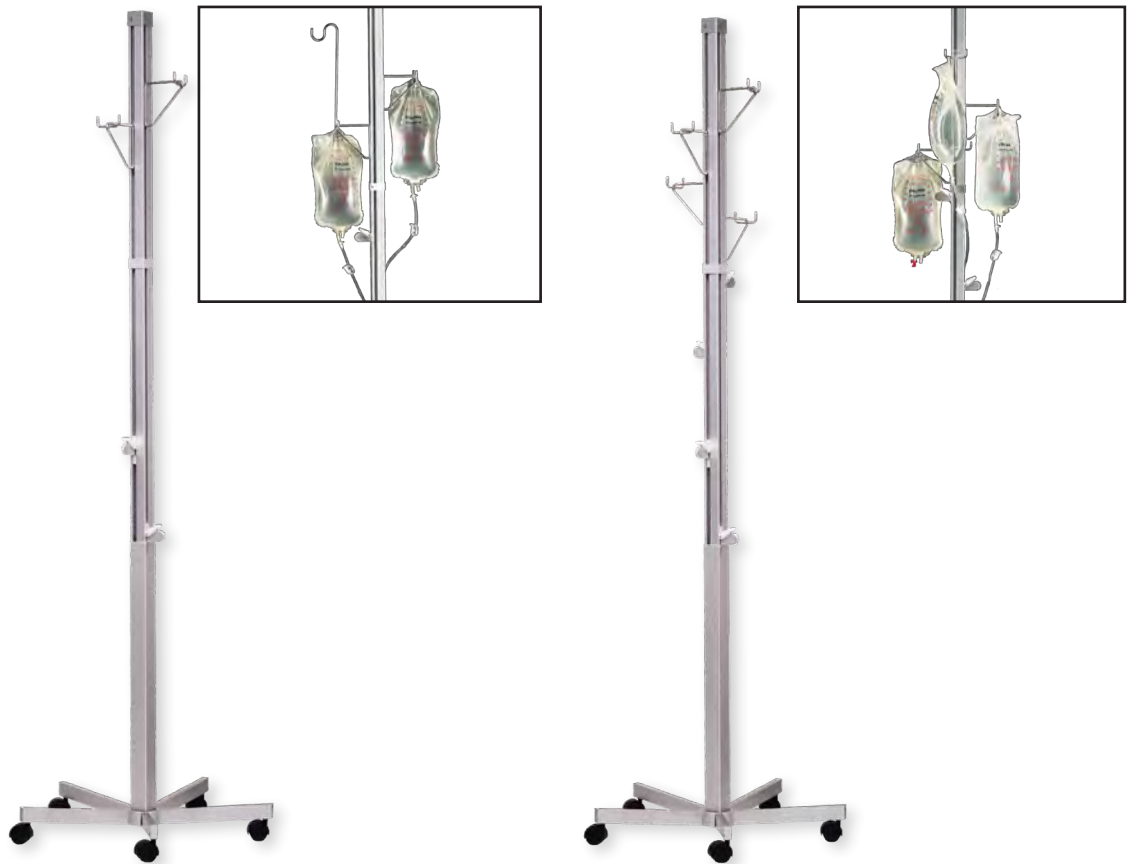

Name	IVT 2 Channels	IVT 4 Channels
Model #	IVT-2	IVT-4
VWR Catalog #	76545-280	76545-278
Shipping Weight	85 lbs.*	85 lbs.*

Shipping weight with an asterisk is the combined shipping weight of separate parts

Irrigation Systems

Construction Specifications

Construction Specifications

Universal Dimensions Tower: Width: 12 3/4" x Height: 89" x Depth: 12 3/4"
Base Diameter: 26 1/2"
Height Range of Hooks: 67" - 81"

Materials Hard anodized aluminum

Gauge N/A

Welding N/A

Finish Stain resistant finishes

Casters 2" black plastic, twin wheel, hooded swivel casters

Shipping Information

All shipping weights are estimates and may be changed without prior notice. Specifications may be improved and dimensions changed without prior notice. Please contact customer service for estimated shipping costs. IV Stand parts shipped separately to avoid damage and to lower costs. Some assembly required.

BLICKMAN

OUR PROMISE

For over 100 years, Blickman has been known throughout the healthcare industry for exceptional product quality. We pride ourselves on providing high-end and personalized service to ensure 100% customer satisfied. If you are not completely satisfied for any reason, we will make it right. All of our equipment is backed by our Limited Lifetime Warranty, so you can be certain that you have made the very best choice for your healthcare facility.

 **avantor**[™]
delivered by **VWR**[™]

VWR.COM

Prices, product, and/or services details are current when published and subject to change without notice. | Certain products or services may be limited by federal, state, provincial, or local regulations. | VWR, part of Avantor, makes no claims or warranties concerning sustainable/green products. Any claims concerning sustainable/green products are the sole claims of the manufacturer and not those of VWR International, LLC and/or Avantor, Inc. or affiliates. All prices are in US dollars unless otherwise noted. Offers valid in US and Canada unless otherwise noted, void where prohibited by law or company policy, while supplies last. | Trademarks are owned by Avantor, Inc. or its affiliates, unless otherwise noted. | Visit vwr.com to view our privacy policy, trademark owners, and additional disclaimers. © 2022 Avantor, Inc. All rights reserved.