

TECHNICAL DATA SHEET

violetFluor™ 450 Anti-Mouse CD45 (30-F11)

Catalog Number: 75-0451

PRODUCT INFORMATION

Contents: violetFluor™ 450 Anti-Mouse CD45 (30-F11)


Isotype: Rat IgG2b, kappa

Concentration: 0.2 mg/mL

Clone: 30-F11

Reactivity: Mouse

Formulation: 10 mM NaH₂PO₄, 150 mM NaCl, 0.09% NaN₃,
0.1% gelatin, pH7.2

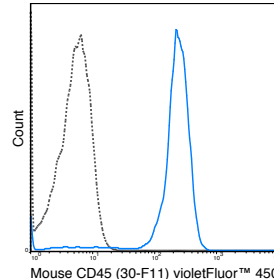

C57Bl/6 splenocytes were stained with 0.5 ug violetFluor™ 450 Anti-Mouse CD45 (75-0451) (solid line) or 0.5 ug violetFluor™ 450 Rat IgG2b isotype control (dashed line).

DESCRIPTION

The 30-F11 antibody reacts with mouse CD45, which is one of the most abundant hematopoietic markers and is expressed on all leukocytes (the Leukocyte Common Antigen, LCA). CD45 is a protein tyrosine phosphatase existing in several isoforms, each being generated and expressed in cell-specific patterns. With its broad cell distribution, CD45 is critical for many leukocyte functions, regulating signal transduction and cell activation associated with the T cell receptor, B cell receptor, and IL-2 receptor. Other forms of CD45, with restricted cellular expression, include CD45R (B220), CD45RA, CD45RB, CD45RO and others. The 30-F11 antibody is widely used as a leukocyte marker for B cells, T cell subsets and NK cell subsets.

PREPARATION & STORAGE

This monoclonal antibody was purified from tissue culture supernatant via affinity chromatography. The purified antibody was conjugated under optimal conditions, with unreacted dye removed from the preparation. It is recommended to store the product undiluted at 4°C, and protected from prolonged exposure to light. Do not freeze.

APPLICATION NOTES

This antibody preparation has been quality-tested for flow cytometry using mouse spleen cells, or an appropriate cell type (where indicated). The amount of antibody required for optimal staining of a cell sample should be determined empirically in your system.

violetFluor™ 450 dye is excited by the violet (405 nm) laser and has a peak emission of 450 nm. The most common band pass filters for this dye are 440/40 or 450/50. violetFluor™ 450 can be used as an alternative for Pacific Blue®, BD Horizon™ V450 or eFluor® 450.

REFERENCES

Panchal RG, Ulrich RL, Bradfute SB, Lane D, Ruthel G, Kenny TA, Iversen PL, Anderson AO, Gussio R, Raschke WC, and Bavari S. 2009. 284: 12874-12885. (immunoprecipitation)Cherpes TL, Busch JL, Sheridan BS, Harvey SAK, and Hendricks RL. 2008. J. Immunol. 181: 969-975. (Complement-mediated cell depletion)Cheng G, Zhang H, Yang X, Tzima E, Ewalt KL, Schimmel P, and Faber JE. 2008. Am. J. Physiol. Regul. Integr. Comp. Physiol. 295: R1138-R1146. (immunohistochemistry – paraffin embedded tissue)Nguyen JT, Evans DP, Galvan M, Pace KE, Leitenberg D, Bui TN, and Baum LG. 2001. J. Immunol. 167: 5697-5707. (immunofluorescence microscopy, immunoprecipitation)Czyzyk J, Leitenberg D, Taylor T, and Bottomly K. 2000. Mol. Cell. Biol. 20(23): 5740-5747. (western blot) Tsuboi S and Fukuda M. 1998. J. Biol. Chem. 273: 30680-30687. (western blot)

NOTE: Please choose the appropriate format for each application. Citations are provided as a convenience to you; please consult Materials and Methods sections for additional details about the use of any product in these publications.

For Research Use Only.

Not for use in diagnostic or therapeutic procedures. Not for resale. Not for distribution without written consent. Tonbo Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Tonbo Biosciences, Tonbo Biosciences Logo and all other trademarks are the property of Tonbo Biosciences Corporation. © 2013 Tonbo Biosciences.