


## 9 Pin Remote Control Cable (07595-47)

Contact Arrangement  
(Cable End)

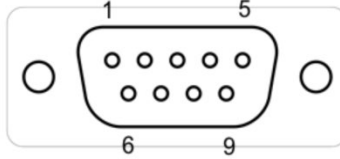

FUNCTION	(+) PIN (COLOR)	(-) PIN (COLOR)
INPUT; Current Speed Control (4-20 mA)	2 (GRY)	3 (RED)
INPUT; Voltage Speed Control (0-10V)	1 (BLK)	3 (RED)
INPUT; Start/Stop (must be used when using 0-10V/4-20 mA external inputs)	7 (ORG)	6 (BLU)
INPUT; Direction (CW/CCW)	8 (YEL)	6 (BLU)
INPUT; Local (keypad speed control)/Remote (external inputs speed control)	4 (GRN)	5 (BRN)
CHASSIS (Earth) GROUND		9 (VIO) DO NOT USE WITH INPUTS
SHELL		DRAIN WIRE