


# Normal Sera



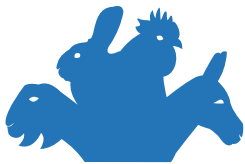

**ImmunoReagents, Inc.**

# Normal Sera

Normal non-immune serum has a variety of uses in immunoassays. It contains a full range of endogenous serum proteins, including all classes of immunoglobulins found in normal, healthy animals. Serum can be used as a positive or negative control in many immunoassays, or it is commonly used as a blocking agent in ELISA, IHC, and Western Blotting. A 5% (v/v) serum solution is sufficient to saturate and block nonspecific sites on fixed tissues or assay substrates.

## Specifications

**Serum:** Non-immune, delipidated

**Concentration:** 50-60 mg/ml

**Form:** Lyophilized or Bulk Liquid

**Unit Sizes:** 2 ml / 5 ml / 10 ml / Bulk Available

Name	Part # / VWR Cat. No.
Goat Serum	SP-004-VX10; 75806-252
Donkey Serum	SP-072-VX10; 75806-232
Guinea Pig Serum	SP-038-VX5; 10152-216
Horse Serum	SP-040-VX10; 10152-218
Human Serum	SP-001-VX10; 10152-222
Llama Serum	SP-132-VX10; 10152-226
Pig/Swine Serum	SP-122-VX10; 10152-234
Rabbit Serum	SP-003-VX10; 10152-238
Rat Serum	SP-006-VX5; 10152-244



800 932 5000 | [VWR.COM](http://VWR.COM)

Prices and product details are current when published and subject to change without notice. | Certain products may be limited by federal, state, provincial, or local regulations. | VWR, part of Avantor, makes no claims or warranties concerning sustainable/green products. Any claims concerning sustainable/green products are the sole claims of the manufacturer and not those of VWR International, LLC and/or Avantor, Inc. or affiliates. All prices are in US dollars unless otherwise noted. Offers valid in US and Canada, void where prohibited by law or company policy, while supplies last. | Trademarks are owned by Avantor, Inc. or its affiliates, unless otherwise noted. | Visit [vwr.com](http://vwr.com) to view our privacy policy, trademark owners, and additional disclaimers. © 2019 Avantor, Inc. All rights reserved.

0120 Lit. No. 090008W

\*Additional species available