

Material Safety Data Sheet

Creation Date 24-Feb-2010

Revision Date 19-Mar-2013

Revision Number 3

1. PRODUCT AND COMPANY IDENTIFICATION

Product Name Shandon Eosin-Y Alcoholic

Cat No. 6766007, 6766008, 9990001

Synonyms No information available.

Recommended Use Laboratory chemicals

Company
Richard Allan Scientific
A Subsidiary of Thermo Fisher Scientific
4481 Campus Drive
Kalamazoo, MI 49008
Tel: (800) 522-7270

Emergency Telephone Number
Chemtrec US: (800) 424-9300
Chemtrec EU: 001 (202) 483-7616

2. HAZARDS IDENTIFICATION

WARNING!

Emergency Overview

Flammable liquid and vapor. Harmful by inhalation, in contact with skin and if swallowed. Irritating to eyes and skin. May cause irritation of respiratory tract. May cause central nervous system effects. Aspiration hazard if swallowed - can enter lungs and cause damage. This substance has caused adverse reproductive and fetal effects in humans.

Appearance Light orange

Physical State Liquid

Odor Alcohol-like

Target Organs Eyes, Skin, Central nervous system (CNS), Liver, Kidney, spleen, Reproductive System

Potential Health Effects

Acute Effects

Principle Routes of Exposure

Eyes

Irritating to eyes.

Skin

Harmful in contact with skin. Irritating to skin.

Inhalation

Harmful by inhalation. Inhalation may cause central nervous system effects. May cause irritation of respiratory tract.

Ingestion

Harmful if swallowed. Aspiration hazard. May cause central nervous system effects. Ingestion may cause gastrointestinal irritation, nausea, vomiting and diarrhea.

Chronic Effects

This substance has caused adverse reproductive and fetal effects in humans. Tumorigenic effects have been reported in experimental animals.. May cause adverse liver effects. May cause adverse kidney effects.

See Section 11 for additional Toxicological information.

Aggravated Medical Conditions

Central nervous system disorders. Gastrointestinal tract. Preexisting eye disorders. Liver disorders. Skin disorders.

3. COMPOSITION/INFORMATION ON INGREDIENTS

Haz/Non-haz

Component	CAS-No	Weight %
Water	7732-18-5	37-40
Ethyl alcohol	64-17-5	54-56
Methyl alcohol	67-56-1	2-3
Isopropyl alcohol	67-63-0	2-3
Acetic acid	64-19-7	<1
Eosin-Y Dye	17372-87-1	<1

4. FIRST AID MEASURES

Eye Contact	Rinse immediately with plenty of water, also under the eyelids, for at least 15 minutes. Immediate medical attention is required.
Skin Contact	Wash off immediately with plenty of water for at least 15 minutes. Immediate medical attention is required.
Inhalation	Move to fresh air. If breathing is difficult, give oxygen. Do not use mouth-to-mouth resuscitation if victim ingested or inhaled the substance; induce artificial respiration with a respiratory medical device. Immediate medical attention is required.
Ingestion	Do not induce vomiting. Call a physician or Poison Control Center immediately.
Notes to Physician	Treat symptomatically.

5. FIRE-FIGHTING MEASURES

Flash Point	22.0°C / 71.6°F
Method -	No information available.
Autoignition Temperature	No information available.
Explosion Limits	
Upper	No data available
Lower	No data available
Suitable Extinguishing Media	Use water spray, alcohol-resistant foam, dry chemical or carbon dioxide.
Unsuitable Extinguishing Media	No information available.
Hazardous Combustion Products	No information available.
Sensitivity to mechanical impact	No information available.
Sensitivity to static discharge	No information available.

Specific Hazards Arising from the Chemical

Flammable. Risk of ignition. Vapors may form explosive mixtures with air. Vapors may travel to source of ignition and flash back. Containers may explode when heated.

Protective Equipment and Precautions for Firefighters

As in any fire, wear self-contained breathing apparatus pressure-demand, MSHA/NIOSH (approved or equivalent) and full protective gear. Thermal decomposition can lead to release of irritating gases and vapors.

NFPA **Health** 2 **Flammability** 3 **Instability** 0 **Physical hazards** N/A

6. ACCIDENTAL RELEASE MEASURES

Personal Precautions	Use personal protective equipment. Remove all sources of ignition. Take precautionary measures against static discharges. Do not get in eyes, on skin, or on clothing.
Environmental Precautions	Should not be released into the environment.
Methods for Containment and Clean Up	Remove all sources of ignition. Soak up with inert absorbent material. Take precautionary measures against static discharges. Keep in suitable, closed containers for disposal..

7. HANDLING AND STORAGE

Handling	Use only under a chemical fume hood. Wear personal protective equipment. Keep away from open flames, hot surfaces and sources of ignition. Take precautionary measures against static discharges. Do not breathe vapors or spray mist. Do not get in eyes, on skin, or on clothing. Do not ingest.
Storage	Keep containers tightly closed in a dry, cool and well-ventilated place. Keep away from heat and sources of ignition. Flammables area.

8. EXPOSURE CONTROLS/PERSONAL PROTECTION

Engineering Measures

Use only under a chemical fume hood. Use explosion-proof electrical/ventilating/lighting/equipment. Ensure that eyewash stations and safety showers are close to the workstation location.

Exposure Guidelines

Component	ACGIH TLV	OSHA PEL	NIOSH IDLH
Ethyl alcohol	STEL: 1000 ppm	(Vacated) TWA: 1000 ppm (Vacated) TWA: 1900 mg/m ³ TWA: 1000 ppm TWA: 1900 mg/m ³	IDLH: 3300 ppm TWA: 1000 ppm TWA: 1900 mg/m ³
Methyl alcohol	TWA: 200 ppm STEL: 250 ppm Skin	(Vacated) TWA: 200 ppm (Vacated) TWA: 260 mg/m ³ (Vacated) STEL: 250 ppm (Vacated) STEL: 325 mg/m ³ Skin TWA: 200 ppm TWA: 260 mg/m ³	IDLH: 6000 ppm TWA: 200 ppm TWA: 260 mg/m ³ STEL: 250 ppm STEL: 325 mg/m ³
Isopropyl alcohol	TWA: 200 ppm STEL: 400 ppm	(Vacated) TWA: 400 ppm (Vacated) TWA: 980 mg/m ³ (Vacated) STEL: 500 ppm (Vacated) STEL: 1225 mg/m ³ TWA: 400 ppm TWA: 980 mg/m ³	IDLH: 2000 ppm TWA: 400 ppm TWA: 980 mg/m ³ STEL: 500 ppm STEL: 1225 mg/m ³
Acetic acid	TWA: 10 ppm STEL: 15 ppm	(Vacated) TWA: 10 ppm (Vacated) TWA: 25 mg/m ³ TWA: 10 ppm TWA: 25 mg/m ³	IDLH: 50 ppm TWA: 10 ppm TWA: 25 mg/m ³ STEL: 15 ppm STEL: 37 mg/m ³

Component	Quebec	Mexico OEL (TWA)	Ontario TWAEV
Ethyl alcohol	TWA: 1000 ppm TWA: 1880 mg/m ³	TWA: 1000 ppm TWA: 1900 mg/m ³	STEL: 1000 ppm
Methyl alcohol	TWA: 200 ppm TWA: 262 mg/m ³ STEL: 250 ppm STEL: 328 mg/m ³ Skin	TWA: 200 ppm TWA: 260 mg/m ³ STEL: 250 ppm STEL: 310 mg/m ³	TWA: 200 ppm STEL: 250 ppm Skin
Isopropyl alcohol	TWA: 400 ppm TWA: 985 mg/m ³ STEL: 500 ppm STEL: 1230 mg/m ³	TWA: 400 ppm TWA: 980 mg/m ³ STEL: 500 ppm STEL: 1225 mg/m ³	TWA: 200 ppm STEL: 400 ppm
Acetic acid	TWA: 10 ppm TWA: 25 mg/m ³ STEL: 15 ppm STEL: 37 mg/m ³	TWA: 10 ppm TWA: 25 mg/m ³ STEL: 15 ppm STEL: 37 mg/m ³	TWA: 10 ppm STEL: 15 ppm

NIOSH IDLH: Immediately Dangerous to Life or Health

Personal Protective Equipment

Eye/face Protection

Wear appropriate protective eyeglasses or chemical safety goggles as described by OSHA's eye and face protection regulations in 29 CFR 1910.133 or European Standard EN166

Skin and body protection

Wear appropriate protective gloves and clothing to prevent skin exposure

Respiratory Protection

Follow the OSHA respirator regulations found in 29 CFR 1910.134 or European Standard EN 149. Use a NIOSH/MSHA or European Standard EN 149 approved respirator if exposure limits are exceeded or if irritation or other symptoms are experienced

9. PHYSICAL AND CHEMICAL PROPERTIES

Physical State

Liquid

9. PHYSICAL AND CHEMICAL PROPERTIES

Appearance	Light orange
Odor	Alcohol-like
Odor Threshold	No information available.
pH	Not applicable
Vapor Pressure	No information available.
Vapor Density	No information available.
Viscosity	No information available.
Boiling Point/Range	Not applicable
Melting Point/Range	No data available
Decomposition temperature	No information available.
Flash Point	22.0°C / 71.6°F
Evaporation Rate	No information available.
Specific Gravity	No information available.
Solubility	No information available.
log Pow	No data available
Molecular Formula	Solution

10. STABILITY AND REACTIVITY

Stability	Stable under normal conditions.
Conditions to Avoid	Incompatible products. Heat, flames and sparks.
Incompatible Materials	Strong oxidizing agents, Acids, Acid anhydrides, Acid chlorides
Hazardous Decomposition Products	Carbon monoxide (CO), Carbon dioxide (CO ₂)
Hazardous Polymerization	Hazardous polymerization does not occur.
Hazardous Reactions	None under normal processing.

11. TOXICOLOGICAL INFORMATION

Acute Toxicity

Product Information No acute toxicity information is available for this product

Component Information

Component	LD50 Oral	LD50 Dermal	LC50 Inhalation
Ethyl alcohol	7060 mg/kg (Rat)	Not listed	20000 ppm/10H (Rat)
Methyl alcohol	5628 mg/kg (Rat)	15800 mg/kg (Rabbit)	64000 ppm (Rat) 4 h 83.2 mg/L (Rat) 4 h
Isopropyl alcohol	5840 mg/kg (Rat)	13900 mg/kg (Rat) 12870 mg/kg (Rabbit)	72.6 mg/L (Rat) 4 h
Acetic acid	3310 mg/kg (Rat)	1060 mg/kg (Rabbit)	11.4 mg/L (Rat) 4 h

Irritation Irritating to eyes and skin

Toxicologically Synergistic Products No information available.

Chronic Toxicity

Carcinogenicity The table below indicates whether each agency has listed any ingredient as a carcinogen.

Component	ACGIH	IARC	NTP	OSHA	Mexico
Ethyl alcohol	A3	Group 1	Not listed	X	Not listed

IARC: (International Agency for Research on Cancer)

IARC: (International Agency for Research on Cancer)

Group 1 - Carcinogenic to Humans

Group 2A - Probably Carcinogenic to Humans

Group 2B - Possibly Carcinogenic to Humans

Sensitization	No information available.
Mutagenic Effects	No information available.
Reproductive Effects	Adverse reproductive effects have occurred in humans..
Developmental Effects	Substances known to cause developmental toxicity in humans.
Teratogenicity	Teratogenic effects have occurred in humans..
Other Adverse Effects	Tumorigenic effects have been reported in experimental animals.. See actual entry in RTECS for complete information.
Endocrine Disruptor Information	No information available

12. ECOLOGICAL INFORMATION

Ecotoxicity

Do not empty into drains.

Component	Freshwater Algae	Freshwater Fish	Microtox	Water Flea
Ethyl alcohol	Not listed	Leucidus idus: LC50 = 8.14 mg/L/48h	Photobacterium phosphoreum: EC50 = 34634 mg/L/30 min Photobacterium phosphoreum: EC50 = 35470 mg/L/5 min	EC50 = 9268 mg/L/48h EC50 = 10800 mg/L/24h
Methyl alcohol	Not listed	Pimephales promelas: LC50 > 10000 mg/L 96h	EC50 = 39000 mg/L 25 min EC50 = 40000 mg/L 15 min EC50 = 43000 mg/L 5 min	EC50 > 10000 mg/L 24h
Isopropyl alcohol	1000 mg/L EC50 > 72 h 1000 mg/L EC50 > 96 h	1400000 µg/L LC50 96 h 9640 mg/L LC50 96 h 11130 mg/L LC50 96 h	= 35390 mg/L EC50 Photobacterium phosphoreum 5 min	13299 mg/L EC50 = 48 h 9714 mg/L EC50 = 24 h
Acetic acid	Not listed	Pimephales promelas: LC50 = 88 mg/L/96h Lepomis macrochirus: LC50 = 75 mg/L/96h	Photobacterium phosphoreum: EC50 = 8.8 mg/L/15 min Photobacterium phosphoreum: EC50 = 8.8 mg/L/25 min Photobacterium phosphoreum: EC50 = 8.8 mg/L/5 min	EC50 = 95 mg/L/24h

Persistence and Degradability No information available**Bioaccumulation/ Accumulation** No information available**Mobility**

Component	log Pow
Ethyl alcohol	-0.32
Methyl alcohol	-0.74

Component	log Pow
Isopropyl alcohol	0.05
Acetic acid	-0.2

13. DISPOSAL CONSIDERATIONS

Waste Disposal Methods

Chemical waste generators must determine whether a discarded chemical is classified as a hazardous waste. Chemical waste generators must also consult local, regional, and national hazardous waste regulations to ensure complete and accurate classification

Component	RCRA - U Series Wastes	RCRA - P Series Wastes
Methyl alcohol - 67-56-1	U154	-

14. TRANSPORT INFORMATION

DOT

UN-No UN1987
Proper Shipping Name ALCOHOLS, N.O.S.
Proper technical name (Ethanol, Methanol)
Hazard Class 3
Packing Group II

TDG

UN-No UN1987
Proper Shipping Name ALCOHOLS, N.O.S.
Hazard Class 3
Packing Group II

IATA

UN-No UN1987
Proper Shipping Name ALCOHOLS, N.O.S.
Hazard Class 3
Packing Group II

IMDG/IMO

UN-No UN1987
Proper Shipping Name ALCOHOLS, N.O.S.
Hazard Class 3
Packing Group II

15. REGULATORY INFORMATION

International Inventories

Component	TSCA	DSL	NDSL	EINECS	ELINCS	NLP	PICCS	ENCS	AICS	CHINA	KECL
Water	X	X	-	231-791-2	-		X	-	X	X	X
Ethyl alcohol	X	X	-	200-578-6	-		X	X	X	X	X
Methyl alcohol	X	X	-	200-659-6	-		X	X	X	X	X

15. REGULATORY INFORMATION											
Isopropyl alcohol	X	X	-	200-661-7	-		X	X	X	X	X
Acetic acid	X	X	-	200-580-7	-		X	X	X	X	X
Eosin-Y Dye	X	X	-	241-409-6	-		X	X	X	X	X

Legend:

X - Listed

E - Indicates a substance that is the subject of a Section 5(e) Consent order under TSCA.

F - Indicates a substance that is the subject of a Section 5(f) Rule under TSCA.

N - Indicates a polymeric substance containing no free-radical initiator in its inventory name but is considered to cover the designated polymer made with any free-radical initiator regardless of the amount used.

P - Indicates a commenced PMN substance

R - Indicates a substance that is the subject of a Section 6 risk management rule under TSCA.

S - Indicates a substance that is identified in a proposed or final Significant New Use Rule

T - Indicates a substance that is the subject of a Section 4 test rule under TSCA.

XU - Indicates a substance exempt from reporting under the Inventory Update Rule, i.e. Partial Updating of the TSCA Inventory Data Base Production and Site Reports (40 CFR 710(B)).

Y1 - Indicates an exempt polymer that has a number-average molecular weight of 1,000 or greater.

Y2 - Indicates an exempt polymer that is a polyester and is made only from reactants included in a specified list of low concern reactants that comprises one of the eligibility criteria for the exemption rule.

U.S. Federal Regulations

TSCA 12(b) Not applicable

SARA 313

Component	CAS-No	Weight %	SARA 313 - Threshold Values %
Methyl alcohol	67-56-1	2-3	1.0
Isopropyl alcohol	67-63-0	2-3	1.0

SARA 311/312 Hazardous Categorization

Acute Health Hazard	Yes
Chronic Health Hazard	Yes
Fire Hazard	Yes
Sudden Release of Pressure Hazard	No
Reactive Hazard	No

Clean Water Act

Component	CWA - Hazardous Substances	CWA - Reportable Quantities	CWA - Toxic Pollutants	CWA - Priority Pollutants
Water	-	1 LB	-	-
Acetic acid	X	5000 lb	-	-

Clean Air Act

Component	HAPS Data	Class 1 Ozone Depletors	Class 2 Ozone Depletors
Methyl alcohol	X		-

OSHA

Not applicable

CERCLA

This material, as supplied, contains one or more substances regulated as a hazardous substance under the Comprehensive Environmental Response Compensation and Liability Act (CERCLA) (40 CFR 302)

Component	Hazardous Substances RQs	CERCLA EHS RQs
Methyl alcohol	5000 lb	-
Acetic acid	5000 lb	-

California Proposition 65

This product contains the following Proposition 65 chemicals: Ethyl alcohol is only a considered a Proposition 65 developmental hazard when it is ingested as an alcoholic beverage.

Component	CAS-No	California Prop. 65	Prop 65 NSRL
Ethyl alcohol	64-17-5	Developmental	-
Methyl alcohol	67-56-1	Methanol	-

State Right-to-Know

Component	Massachusetts	New Jersey	Pennsylvania	Illinois	Rhode Island
Ethyl alcohol	X	X	X	-	X
Methyl alcohol	X	X	X	X	X
Isopropyl alcohol	X	X	X	-	X
Acetic acid	X	X	X	-	X

U.S. Department of Transportation

Reportable Quantity (RQ): Y
 DOT Marine Pollutant N
 DOT Severe Marine Pollutant N

U.S. Department of Homeland Security

This product does not contain any DHS chemicals.

Other International Regulations


Mexico - Grade Moderate risk, Grade 2


Canada

This product has been classified in accordance with the hazard criteria of the Controlled Products Regulations (CPR) and the MSDS contains all the information required by the CPR.

WHMIS Hazard Class

B2 Flammable liquid
 D1B Toxic materials
 D2A Very toxic materials

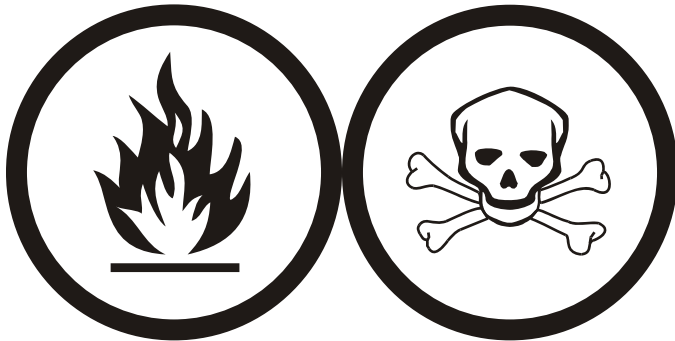

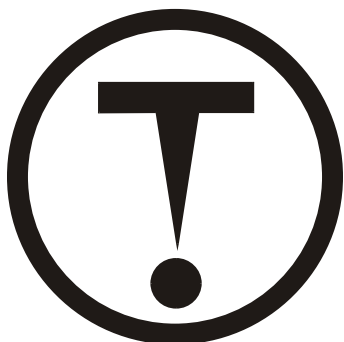

16. OTHER INFORMATION

Prepared By Regulatory Affairs
Thermo Fisher Scientific
Email: EMSDS.RA@thermofisher.com

Creation Date 24-Feb-2010

Print Date 19-Mar-2013

Revision Summary "****", and red text indicates revision

Disclaimer

The information provided on this Safety Data Sheet is correct to the best of our knowledge, information and belief at the date of its publication. The information given is designed only as a guide for safe handling, use, processing, storage, transportation, disposal and release and is not to be considered as a warranty or quality specification. The information relates only to the specific material designated and may not be valid for such material used in combination with any other material or in any process, unless specified in the text.

End of MSDS