

Biotix xTIP Pipette Tips

xPIPETTE™ and LTS™ Compatible

Perfect Ergonomic System

FlexFit®




Flexible Proximal End

- Requires less insertion force than competitors' to make a secure seal, reducing the risk of repetitive stress injury
- Tips don't fall off multichannel pipettes
- Completes the ergonomic system on LTS pipettes

StarStop™

Positive Stop Feature

- Gives consistently low ejection forces on LTS pipettes
- Aligns the tips for even liquid levels on multichannels

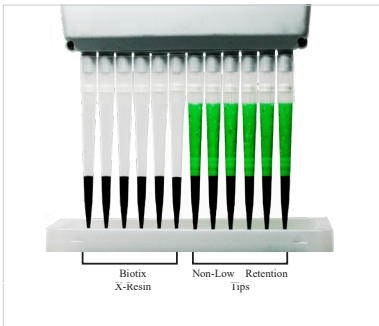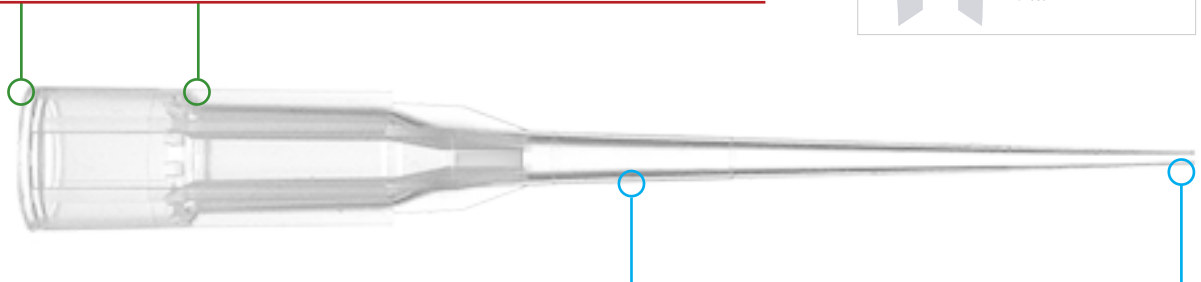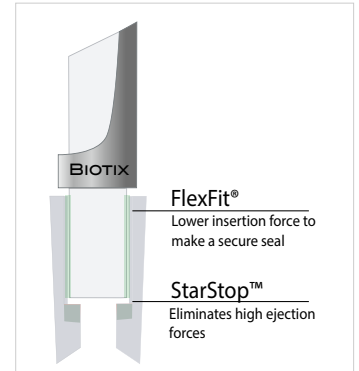

Clean Delivery

X-Resin®

Low Retention Surface

- Maximizes sample recovery
- Improves accuracy across your experiments

Blade®

Minimizes Surface Tension

- Eliminates hanging droplet formation
- Saves time by reducing the need for tip touch off
- Increases accuracy and precision

Volume μ L	Pre-Sterilized	Racked, Filtered	Racked	CleanPak Reload	Tip Eject Reload	Packaging
20	Yes	76314-472	76266-150		76266-128	960 tips/pack
	No		76266-098	76266-118		
200	Yes	76266-142				960 tips/pack
	No					
250	Yes		76266-152		76266-130	960 tips/pack
	No		76266-100	76266-120		
300	Yes	76266-144	76266-092			960 tips/pack
	No		76266-102	76266-122		
1000	Yes	76266-146	76266-094		76266-134	768 tips/pack
	No		76266-104	76266-124		
1200	Yes	76266-148	76266-096		76266-136	768 tips/pack
	No		76266-106	76266-126		

A More Sustainable Rack

The Biotix Rack provides top functionality and strength while using 1/3 less plastic than leading competitors

A more sustainable design with all the desired features:



- More rigid lid and base design for repeatable autoclave
- Locking clasp clicks to ensure rack is securely closed during transport
- Snap cards are easily removed with reloads by one hand
- Modified hole design prevents rack lift up during tip removal at any angle
- Lid opens past 180° and is easily removable for optimal multichannel pipetting
- Compatible with both CleanPak and Tip Eject Reloads

Reload Solutions

Biotix xTIP has two reload styles to reduce plastic waste in your lab. Simply reuse your empty racks.



xTIP Tip Eject™ Reload

Tip Eject is a self-aligning reload solution. Simply press on the top of the stack and the empty rack will be refilled at the bottom.

- Reduces waste by 30% when compared to racked product
- Sterilized after packaging for added protection
- Rack boxes are autoclaveable after loading at 121° C
- Certified free of DNase, RNase, endotoxins, pyrogens, ATP, and PCR inhibitors



xTIP CleanPak™ Reload

CleanPak reloads arrive hermetically sealed in a ziploc bag for superior tip purity. Available in both sterilized and non-sterile formats.

- Reduces waste by 50% when compared to racked product
- Color coded snapcards match the compatible pipette's plunger cap
- To refill an empty rack just reach into the bag and manually transfer a snapcard
- Rack boxes are autoclaveable after loading at 121° C
- Certified free of DNase, RNase, endotoxins, pyrogens, ATP, and PCR inhibitors