



PALL CORPORATION

avantor™  
delivered by VWR™

VWR.COM

## FOR A LIMITED TIME Get up to 25% off select Acrodisc® syringe filters

### wwPTFE Membrane


The water wettable hydrophilic polytetrafluoroethene (wwPTFE) membrane is an all-purpose, hydrophilic membrane for aqueous, acidic, basic, non-aggressive organic, and aggressive organic solutions. The wwPTFE membrane offers low protein binding and low levels of UV-absorbing extractables.

### wwPTFE Membrane Performance

The wwPTFE membrane retention efficiency can directly increase the life of your HPLC column. More retentive filters will also remove smaller non-dissolved particles that may rest on the head of your HPLC column and affect your analytical results as they dissolve during an analysis. The graph below demonstrates that the life of an HPLC column can be extended up to 52 times by filtering samples prior to injection.

### Column Life Extension

Effects of Filters on HPLC Column Life following injections of unfiltered and filtered 0.05% latex sphere suspensions. With unfiltered samples, the column failed due to plugging after 19 injections. Samples passed through Competitor filters plugged the columns after 500 injections. No increase in backpressure was observed after 1000 injections of samples filtered with the wwPTFE membrane.

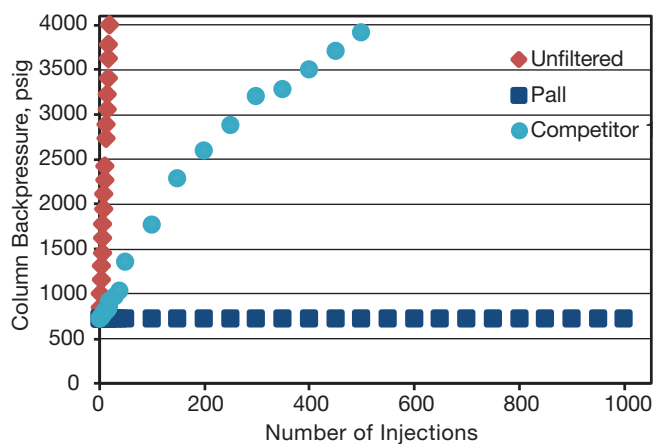

## Ordering Information

### Acrodisc Syringe Filters, 13 mm with Universal wwPTFE Membrane

VWR Cat. No.	Pall Part No.	Description	Pkg
MP76308-728	2400	0.2 µm wwPTFE, with minispikes, female luer lock/ male slip luer, non-sterile	300/pkg
MP76308-698	2401	0.2 µm wwPTFE, with minispikes, female luer lock/ male slip luer, non-sterile	1000/pkg
MP76308-700	2402	0.45 µm wwPTFE, with minispikes, female luer lock/ male slip luer, non-sterile	50/pkg
MP76308-702	2403	0.45 µm wwPTFE, with minispikes, female luer lock/ male slip luer, non-sterile	1000/pkg

### Acrodisc Syringe Filters, 25 mm with Universal wwPTFE Membrane

VWR Cat. No.	Pall Part No.	Description	Pkg
MP76308-634	4927	0.2 µm wwPTFE, non-sterile	200/pkg
MP76308-636	4929	0.2 µm wwPTFE, non-sterile	1000/pkg
MP76308-638	4914	0.45 µm wwPTFE, non-sterile	200/pkg
MP76308-640	4932	0.45 µm wwPTFE, non-sterile	1000/pkg

### Acrodisc One PSF Syringe Filters, 25 mm with wwPTFE Membrane

VWR Cat. No.	Pall Part No.	Description	Pkg
MP76308-662	AP-4910	0.2 µm wwPTFE, non-sterile	200/pkg
MP76308-666	AP-4912	0.2 µm wwPTFE, non-sterile	1000/pkg
MP76308-668	AP-4913	0.2 µm GXF/wwPTFE non-sterile	200/pkg
MP76308-672	AP-4915	0.2 µm GXF/wwPTFE non-sterile	1000/pkg
MP76308-674	AP-4916	0.45 µm wwPTFE non-sterile	200/pkg
MP76308-678	AP-4918	0.45 µm wwPTFE non-sterile	1000/pkg
MP76308-680	AP-4919	0.45 µm GXF/wwPTFE non-sterile	200/pkg
MP76308-684	AP-4921	0.45 µm GXF/wwPTFE non-sterile	1000/pkg


Receive up to 25% off of VWR published list price Pall Acrodisc Syringe Filters between December 3, 2021 to December 31, 2022. Products Ordered must be from list associated for this promotion. Promotion ends December 31, 2022. Visit [vwr.com/pall](http://vwr.com/pall) for information on qualifying products. For further details contact your VWR or Pall rep.

Promotional part numbers are available for a limited time, must be used at the time of purchase, and may be ordered via the "ORDER ENTRY" link on [vwr.com](http://vwr.com). (Promotional part numbers cannot be found via Search on [vwr.com](http://vwr.com).) For assistance with ordering, please contact VWR Customer Service at 1 800 932 5000.



Ordering:

[vwr.com](http://vwr.com)

© Copyright 2022, Pall Corporation. Pall, , Acrodisc One, AcroPrep, and Nanosep are trademarks of Pall Corporation. ® Indicates a trademark registered in the USA. Limit one per customer. Offer not valid where prohibited by law or institutional policy. Valid only in both US and CA. Allow up to 8 weeks for delivery of item. Offer expires December 31, 2022. GN21.1006

7/22



Prices, product, and/or services details are current when published and subject to change without notice. | Certain products or services may be limited by country, federal, state, provincial, or local regulations. | VWR, part of Avantor, makes no claims or warranties concerning sustainable/green products. Any claims concerning sustainable/green products are the sole claims of the manufacturer and not those of VWR International, LLC and/or Avantor, Inc. or affiliates. Offers valid in countries listed above, void where prohibited by law or company policy, while supplies last. | Trademarks are owned by Avantor, Inc. or its affiliates, unless otherwise noted. | Visit [vwr.com](http://vwr.com) to view our privacy policy, trademark owners, and additional disclaimers. © 2022 Avantor, Inc. All rights reserved.

0722 Lit. No. 160268WREV Promo Code 5813