

ACE-THRED™ FILTRATION APPARATUS

Using our famous Ace-Thred™ design, we are able to thread the filter funnel and filter flask together without the silicone stopper and troublesome clamping device. This design also eliminates leakage problems. Various porosities of filter discs are available. Fritted discs are easily removable, as it is not permanently sealed.

500mL Series Apparatus

- 500mL filter flask
- #7 Ace-Safe™ connector
- #25 to #25 Ace-Thred™ nylon adapter
- 25mm fritted disc
- 100mL graduated funnel

1000mL Series Apparatus

- 1000mL filter flask
- #7 Ace-Safe™ connector
- #25 to #50 Ace-Thred™ nylon adapter
- 47mm fritted disc
- 500mL graduated funnel

Instatherm Series Apparatus

- Same great design as 1000mL series
- Replaces the funnel with Instatherm® version
- Heats viscous materials for easier flow

Also Available...

Contact us for all your filtration needs - both, large and bench scale applications.



500mL Series



1000mL Series




FILTRATION APPARATUS continued...

2L Filter Flask Apparatus

- #15 Ace-Safe™ connector
- #50 to #80 Ace-Thred™ PTFE adapter
- 75mm fritted disc
- 1L graduated funnel

5L Bottle Apparatus

- #15 Ace-Safe™ connector
- #80 Ace-Thred™ to GLS80 PTFE adapter
- 75mm fritted disc
- 1L graduated funnel

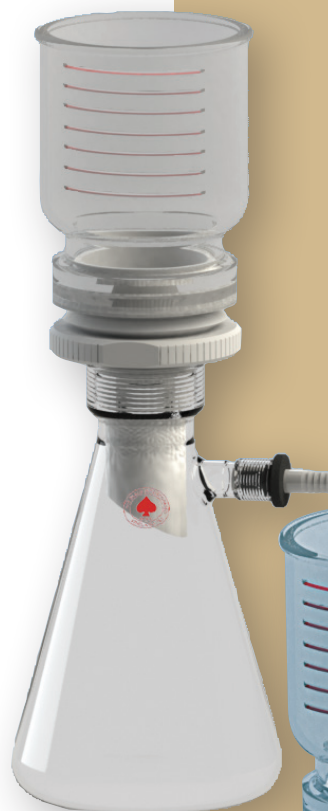

Jacketed Funnels

- For 75mm membranes or filter paper
- Ideal for solvent recycling and other filter applications where temperature control is required
- 1L & 2L sizes available



Contact us for...

Custom glass and repair services.



1.800.932.5000 | VWR.COM

Prices and product details are current when published and subject to change without notice. | Certain products may be limited by federal, state, provincial, or local regulations. | VWR, part of Avantor, makes no claims or warranties concerning sustainable/green products. Any claims concerning sustainable/green products are the sole claims of the manufacturer and not those of VWR International, LLC and/or Avantor, Inc. or affiliates. All prices are in US dollars unless otherwise noted. Offers valid in US, void where prohibited by law or company policy, while supplies last. | Trademarks are owned by Avantor, Inc. or its affiliates, unless otherwise noted. | Visit vwr.com to view our privacy policy, trademark owners, and additional disclaimers. © 2019 Avantor, Inc. All rights reserved.