



Moticam **4000BMH** HDMI & USB

FULL HD CAMERAS

The Moticam 4000 BMH is a C-Mount multitasking microscopy camera, that delivers an impressive 1080p Full HD live image with maximum colour fidelity on an integrated LCD screen.

It does not require a computer for standard operation, as it comes with an on-board software that will allow you to use it through its 11.6" Full HD integrated screen, with a mouse connected to the USB port of the camera.



<https://www.motic.com/en/>

FULL HD STANDALONE CAMERA





Moticam 4000BMH TECHNICAL SPECIFICATIONS

Sensor	IMX334
Sensor Size	1/1.8"(8.86mm)
Imaging Area	7.68 x 4.32mm
Capture Resolution	8MP
Scan mode	Progressive
Shutter mode	Rolling Shutter
Max. frames per second (fps*)	USB: 30fps@3840*2160 30fps@1920*1080 HDMI: 20fps@3840*2160
Image Capture	JPEG/TIFF Format with 4K(3840*2160) Resolution in SD Card
Video Record	4K(30fps@3840*2160) H264/H265 encoded MP4 file in SD Card
Screen	LCD 11.6" Full HD
Software	Motic Images Plus 3.1 for Windows/Mac OS
On-board software	Standard tools (Camera connected through HDMI to a HD screen and through USB to a mouse) Still image capture, Live image record, Zoom in & out, Magnifier, ROI, Mirror, Rotation, Split, Freeze, Cross line, Gallery.
	Image adjustments (AE/AWB) Auto/Manual Exposure mode, Gain, Target Brightness setting, Light Frequency, Auto/Manual White Balance mode.
	Advanced settings Saturation, Contrast, Gamma, Sharpness, Denoise, User settings, System's Time setting, Language (English, Chinese and Japanese), Format SD card.

*frames per second under optimal illumination conditions



*Design Change :The manufacturer reserves the right to make changes in instrument design in accordance with scientific & mechanical progress, without notice and without obligation.

ON-BOARD SOFTWARE & MOTIC IMAGES PLUS 3.1 SOFTWARE

THE [All in One BOX] CONCEPT / Everything that you may need to work with your Moticam is included in one box

On-board software

The Moticam 4000 BMH has its own built-in software, making it possible to view live images, make image adjustments or edit them. You can also capture images or record videos directly from a screen with the help of a mouse, storing them in an SD card.

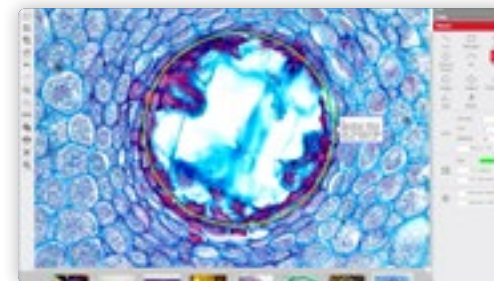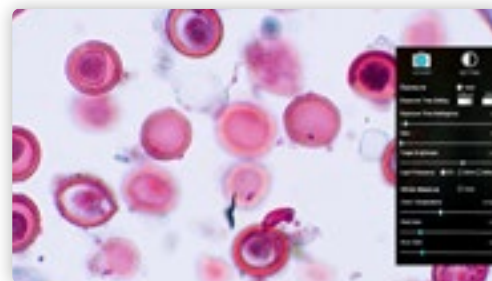

Motic Images Plus 3.1

By connecting your Moticam to your computer you will be able to work with our well-known software, that comes included in the package. View, capture, edit, measure, make reports... all its standard features and the new ones packed in a new user-friendly interface.