

Midas® III-Plus Automated Stainer


Automated Side Stainer for Hematology and Bacteriology

The Midas® III-Plus Automated Stainer (**VWR Cat. No. 15204-300**) gives you the ability to automate all hematology and bacteriology staining applications. Since the unit is fully programmable, you can duplicate your existing staining techniques and ensure consistent results for your peripheral blood and bone marrow smears as well as your bacteriological smears.

MilliporeSigma backs up the Midas® III-Plus Automated Stainer with a one year warranty and some of the best technical and customer service in the business. Don't forget we have a full line of stains (including the Harleco® product line) perfectly suited for use with the Midas® III-Plus Automated Stainer.

Key Features

- Ideally suited for Wright (**VWR Cat. No. 15204-228**) & Wright Giemsa (**VWR Cat. No. 15204-232**) staining of blood and bone marrow smears
- Ease of use for Gram and TB staining applications
- Attractive blue/light blue color for the instrument
- Optional clear colorless enclosure with door

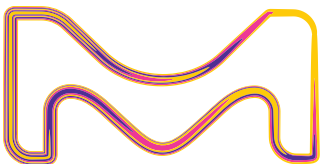

The Cornerstone of Your Laboratory

The Midas® III-Plus Automated Stainer is a fully automated, benchtop stainer with a small footprint. It is designed to perform routine hematology and bacteriology staining effortlessly. It utilizes the dip slide staining technique which ensures uniform staining for every sample. This reduces stain consumption up to 50%.

This versatile automated stainer combines the reliable electronics and mechanics of the original Midas® III Stainer, now with a new optional plastic hood enclosure to make your laboratory a safer environment. Ideally suited for staining up to 1000 slides per day, this stainer will become the new workhorse of your hematology/bacteriology lab.

Technical Specifications

Specification	Description
Staining Control	Microprocessor
Available programs	9 user-defined programs (up to 29 steps each), 1 test program
Application time	1 second to 99 minutes
Available stations	6 stations: 4 reagent, 1 drying, 1 rinse (Option to turn drying and rinse stations to reagent)
Stain Vessel	8 vessels & lids
Capacity	300 ml
Slide Carrier	1 20 position slide carrier
Capacity	20 or 30 slides
Rinse Station Flow	Program maintains 500 ml/min increments allowing for precise rinsing
Length	21" (53 cm)
Width	18" (45 cm)
Height	15" (38 cm)
Weight	28 lbs. (13 kg)
Shipping Weight	60 lbs.
UL/cUL approved	

Indication for Use:

For in vitro diagnostic use

IVD

Intended User:

For professional use only

ca.vwr.com/emdmillipore


Ordering Information

Product	VWR Cat. No.	Size	Package
Midas® III-Plus Slide Stainer	CA15204-300	1EA	Chipboard Box
Midas® Coverslip Carrier	25445-044	1EA	Chipboard Box
Midas® Slide Carrier 20 slides	CA64053-94	1EA	Chipboard Box
Midas® Slide Carrier 30 slides	CA64052E-94	1EA	Chipboard Box
Midas® Vessels and Lids White 8/pk	CAD64054W-95	1EA	Chipboard Box
Midas® Water Bath	CA64000WB-94	1EA	Chipboard Box
Midas® Water Bath Filters 2/pk	25445-016	1EA	Chipboard Box
Midas® Vessel Rack	CARD11953	1EA	Chipboard Box
Midas® Plastic hood enclosure	CA10183-814	1EA	
Midas® 110V power cord (US)	CARW918003	1EA	

We provide information and advice to our customers on application technologies and regulatory matters to the best of our knowledge and ability, but without obligation or liability. Existing laws and regulations are to be observed in all cases by our customers. This also applies in respect to any rights of third parties. Our information and advice do not relieve our customers of their own responsibility for checking the suitability of our products for the envisaged purpose.

© 2019 Merck KGaA, Darmstadt, Germany and/or its affiliates. All Rights Reserved. MilliporeSigma and the vibrant M are trademarks of Merck KGaA, Darmstadt, Germany or its affiliates. All other trademarks are the property of their respective owners. Detailed information on trademarks is available via publicly accessible resources.

2018-16704 12/2018

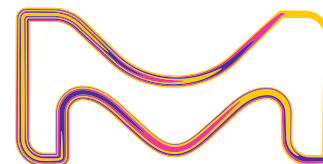

VWR
part of avantor

800 932 5000 | VWR.COM

Prices and product details are current when published and subject to change without notice. | Certain products may be limited by federal, state, provincial, or local regulations. | VWR, part of Avantor, makes no claims or warranties concerning sustainable/green products. Any claims concerning sustainable/green products are the sole claims of the manufacturer and not those of VWR International, LLC and/or Avantor, Inc. or affiliates. All prices are in Canadian dollars unless otherwise noted. Offers valid in Canada, void where prohibited by law or company policy, while supplies last. | Trademarks are owned by Avantor, Inc. or its affiliates, unless otherwise noted. | Visit vwr.com to view our privacy policy, trademark owners, and additional disclaimers. © 2019 Avantor, Inc. All rights reserved.