

MPR-722-PA | 76355-506  
MPR-722R-PA | 76355-508



MPR-722R-PA | 76355-508

# Pharmaceutical Refrigerator

High Performance Cooling  
with EPA SNAP Compliant Refrigerants

24.2 cu.ft. | 684 L  
23.7 cu.ft. | 671 L

The MPR-722-PA | 76355-506 (shelves) and MPR-722R-PA | 76355-508 (drawers) pharmaceutical refrigerators deliver a stable, reliable temperature control for clinical, pharmaceutical, biomedical storage and processing. Fast recovery following door openings assures safe storage during times of frequent access and high ambient conditions and to remove heat loads generated by powered instruments.

## Positive Airflow Assures Quick Recovery, Interior Uniformity

Vertical, forced air circulation with dual blower creates a uniform temperature at all shelf or drawer levels, top to bottom and front to back. Internal shelf or drawer clearances permit the positive airflow also restores temperature to setpoint following door openings, and safely accommodates busy, high traffic demand for access to stored contents in pharmacies, clinics and hospitals.



## Positive Airflow Prevents Temperature Stratification

The top-mounted air circulation blower creates a circular airflow within the chamber to provide uniform temperatures at all shelf levels to enhance and protect viability of stored product. Regardless of location, all products are subjected to the same safe temperature to assure repeatability of stored product viability.

## Precision Temperature Control Protects Inventory

Because many pharmaceuticals are degraded or destroyed if accidentally frozen, it is important that temperature be adjusted to desired setpoint in accordance with drug manufacturers' insert specifications. MPR Series pharmaceutical refrigerators operate over a temperature range of 2°C to 23°C, adjustable in 1.0°C increments. Actual interior air temperature is displayed to 0.1°C.



## Microprocessor Controls Simplify Operation

A microprocessor controller monitors all functions. Setpoint is factory set at 5°C. A soft key interface allows setpoint adjustment throughout the operating range. Audible and visual warnings indicate a door ajar or temperature deviation from setpoint. The large easy-to-read LED display confirms actual temperature.

## Choose Your Inventory Management Preference

The MPR Series pharmaceutical refrigerator offers two choices in inventory management. The MPR-722-PA | 76355-506 includes 4 open-wire shelves strong enough to hold bulk loads. The MPR-722R-PA | 76355-508 is equipped with heavy-duty, fully-extendable painted steel pull-out drawers for more convenient access to stored product such as large bottles or reagent kits.



## Adjustable Shelves Offer Flexibility in Storage Options

Open wire shelves (shown) are adjustable to allow for loads of varying heights such as bulk storage or process equipment such as fraction collectors. The Model MPR-722R-PA | 76355-508 with pull-out drawers (not shown) offers a practical solution with additional ergonomic advantages.

# Upright Pharmaceutical Refrigerator

MPR-722-PA | 76355-506  
MPR-722R-PA | 76355-508




Model MPR-722-PA | 76355-506

## Ideal for Storage or Process Instrumentation

The MPR-722-PA | 76355-506 Series pharmaceutical refrigerator is configured to provide close temperature control for both storage and process applications where heat is generated from instrumentation or apparatus. A powerful refrigeration system using high performance compressors and components delivers fast temperature pull-down, recovery and tolerance for high ambient temperature conditions. Energy efficiency helps improve LEED scores and contribute to global sustainability.


## Cycle Defrost Function

The cycle defrost and evaporator temperature sensor system ensures that defrost occurs only when necessary and automatically, so there is no need to turn off the power for defrosting. Irregular temperature increase during defrost is minimal with no temperature spikes. The evaporation heater also doubles as protection against drops in cabinet temperature caused by a low ambient temperature.

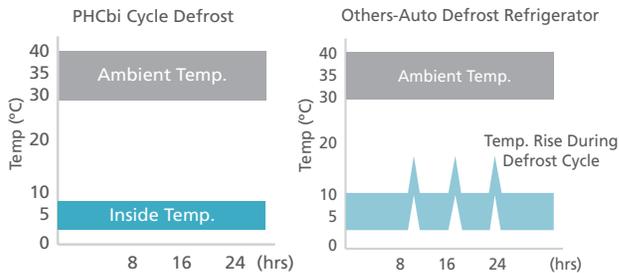

## Dimensions

Unit : inch | mm

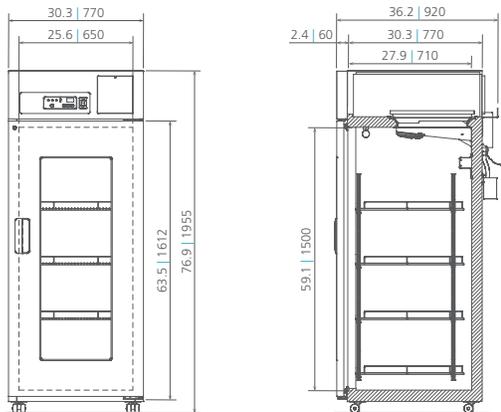

Model Number	MPR-722-PA   76355-506		MPR-722R-PA   76355-508	
External Dimensions (W x D x H) <sup>1)</sup>	inches   mm		30.3 x 36.2 x 76.9   770 x 920 x 1955	
Internal Dimensions (W x D x H)	inches   mm		25.6 x 27.9 x 59.1   650 x 710 x 1500	
Volume	cu.ft.   liters		24.2   684	
Net Weight	lbs   kg		334   174	
<b>Performance</b>				
Power Consumption	watts	400		
Warranty <sup>2)</sup>	2 year parts and labor, 3 year compressor parts			
Vaccine Storage Operating Temperature	°C	+2 to +8		
Factory Pre-Set Temperature	°C	+5		
Highest Ambient Temperature and Maintains Cabinet Temperature	°C	+2 to +23 (Ambient temperature 0 to +35, no load)		
Evaporator Prevents Vaccines from Freezing	°C	Operates above freezing at all times		
Temperature Control Range	°C	+2 to +23		
<b>Control</b>				
Microprocessor Controller, Adjustable	°C	Increments of 1.0		
Digital Temperature Display	°C	LED (displayed by 1, 0.1)		
Controller Security	Lockable with key pad			
Electronics Diagnostics	Total control system			
<b>Refrigeration</b>				
Cooling Method	Uniformity forced air			
Internal Airflow for Precise Temperature	Uniformity forced air			
Defrost Method Initiated Only as Needed	Electronically monitored evaporator			
Refrigeration System	Air cooled, CFC free			
Insulation	CFC free urethane			
<b>Construction</b>				
Outer Door	qty	1 (dual pane glass)		
Interior	Zinc galvanized steel, Acrylic finish			
Exterior	Zinc galvanized steel, Acrylic finish			
Outer Door Lock	Key			
Interior Light	Yes - with control panel switch			
Shelves	qty	4 (adjustable, wire)		—
Drawers	qty	—		5 (solid)
Casters	qty	4 (2 swivel, 2 fixed)		
Access Port	qty	3 (left, right, top)		
Access Port Diameter	inches   mm	1.1   30		
<b>Alarms</b> (V = Visual Alarm, B = Buzzer Alarm, R = Remote Alarm)				
Power Failure	R (V-B optional)			
Temperature (High and Low)	V-B-R			
Door Open	V-B			
Remote Alarm Contacts	Normally open, Normally closed, Common			
Remote Alarm Output	DC 24V 2A			
<b>Electrical and Noise Level</b>				
Power Supply			115V, 1Ø, 60Hz, NEMA 5-15P requires 5-15R receptacle	
Noise Level <sup>3)</sup>	dB(A)		48	
<b>Options</b>				
Blackout Panels to Cover Glass	qty	1		
<b>Optional Communication System</b>				
Wireless, Cloud-Based, Automatic Data Management	LabAlert® Monitoring System			

<sup>1)</sup> Exterior dimensions of main cabinet only, excluding handle and other external projections. Consult product sales rep for doorway entry instructions, less than 36.2".

<sup>2)</sup> Current warranty offered at time of printing and may be subject to change; US and Canada only

<sup>3)</sup> Actual value, background noise 20 dB(A)

Specifications are subject to change without notice.

Performance data herein is based on independent testing at time of publication.

Prices, product, and/or services details are current when published and subject to change without notice. | Certain products or services may be limited by federal, state, provincial, or local regulations. | VWR, part of Avantor, makes no claims or warranties concerning sustainable/green products. Any claims concerning sustainable/green products are the sole claims of the manufacturer and not those of VWR International, LLC and/or Avantor, Inc. or affiliates. All prices are in US dollars unless otherwise noted. Offers valid in US and Canada, void where prohibited by law or company policy, while supplies last. | Trademarks are owned by Avantor, Inc. or its affiliates, unless otherwise noted. | Visit [vwr.com](http://vwr.com) to view our privacy policy, trademark owners, and additional disclaimers. © 2020 Avantor, Inc. All rights reserved.