

BOEKEL PLATELET AGITATORS & INCUBATORS

Advancements in Blood Banking

Boekel Scientific offers two sizes of advanced platelet storage agitators and incubators. The systems were designed closely with blood banking professionals to incorporate useful and energy saving innovations.

- Patent pending Inventory Management System
- No refrigerant or compressor for service free heating and cooling
- Onboard thermal printer for keeping records of temperature data as well as a downloadable log
- Automated Hi-Lo Alarm Testing
- 7" industrial touch screen with intuitive icons
- Battery backup in the event of a power failure
- Two year warranty included



*Agitators and Incubators

Sold Separately



	SMALL UNIT	LARGE UNIT
Platelet Agitator	Model # 301200	Model # 301300
Platelet Incubator	Model # 301550	Model # 301650
Electrical	115V, 1 Phase, 50/60 H	115V, 1 Phase, 50/60 H
Bag Capacity	6 - 12 Apherisis, 18 Random (6 Shelves)	8 - 16 Apherisis, 48 Random (8Shelves)
Incubator Dimensions	21.25"W x 19.75"D x 26"H	25.3"W x 26.5"D x 30.0"H

PLATELET INCUBATORS & AGITATORS




Advancements in Blood Banking



On-Board Thermal Printer for quick and easy temperature record keeping



Incubator Thermo-electric heating and cooling eliminates the old style heater, compressor, and refrigerant

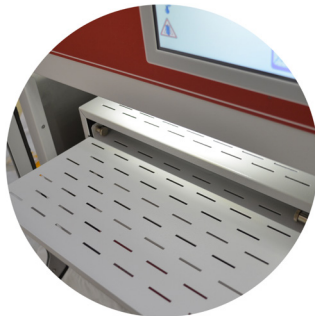

Easy sliding and removable agitator shelves



Patent Pending Inventory Management System



Giant color touch screen with USB downloadable event and temperature log