

Safety Data Sheet

SECTION 1: Identification

1.1. Product Identifier

Trade Name or Designation: Calcium Chloride, 10 g/L in Mixed Alcohol Solvent
for determination of Salts in Crude Oil

Product Number: 1718

Other Identifying Product Numbers: 1718-32

1.2. Recommended Use and Restrictions on Use

General Laboratory Reagent

1.3. Details of the Supplier of the Safety Data Sheet

Company: Ricca Chemical Company

Address: 448 West Fork Drive
Arlington, TX 76012 USA

Telephone: 888-467-4222

1.4. Emergency Telephone Number (24 hr)

CHEMTREC (USA) 800-424-9300
CHEMTREC (International) 1+ 703-527-3887

Safety Data Sheet

SECTION 2: Hazard(s) Identification

2.1. Classification of the Substance or Mixture (in accordance with OSHA HCS 29 CFR 1910.1200)

For the full text of the Hazard and Precautionary Statements listed below, see Section 16.

Hazard Class	Category	Hazard Statement	Precautionary Statements
Acute Toxicity - Oral	Category 3	H301	P264, P270, P301+P310, P321, P330, P405, P501
Acute Toxicity - Dermal	Category 4	H312	P280, P302+P352, P312, P321, P362, P501
Acute Toxicity - Inhalation	Category 4	H332	P261, P271, P304+P340, P312
Skin Corrosion / Irritation	Category 2	H315	P264, P280, P302+P352, P321, P332+P313, P362
Eye Damage / Irritation	Category 1	H318	P280, P305+P351+P338, P310
Reproductive Toxicity	Category 1	H360	P201, P202, P280, P308+P313, P405, P501
Specific Target Organs/Systemic Toxicity Following Single Exposure	Category 1	H370	P260, P264, P270, P307+P311, P321, P405, P501
Specific Target Organs/Systemic Toxicity Following Repeated Exposure	Category 1	H372	P260, P264, P270, P314, P501
Flammable Liquids	Category 2	H225	P210, P233, P240, P241, P242, P243, P280, P303+P361+P353, P370+P378, P403+P235, P501

2.2. GHS Label Elements

Pictograms:



Signal Word: **Danger**



Safety Data Sheet

Hazard Statements:

Hazard Number	Hazard Statement
H225	Highly flammable liquid and vapor.
H301	Toxic if swallowed.
H312	Harmful in contact with skin.
H315	Causes skin irritation.
H318	Causes serious eye damage.
H332	Harmful if inhaled.
H360	May damage fertility or the unborn child.
H370	Causes damage to organs.
H372	Causes damage to organs through prolonged or repeated exposure.



Safety Data Sheet

Precautionary Statements:

Precautionary Number	Precautionary Statement
P201	Obtain special instructions before use.
P202	Do not handle until all safety precautions have been read and understood.
P210	Keep away from heat, sparks and open flame. No smoking.
P233	Keep container tightly closed.
P240	Ground container and receiving equipment.
P241	Use explosion-proof equipment.
P242	Use only non-sparking tools.
P243	Take precautionary measures against static discharge.
P260	Do not breathe dust, fumes or mist.
P261	Avoid breathing dust, fumes or mist.
P264	Wash arms, hands and face thoroughly after handling.
P270	Do not eat, drink or smoke when using this product.
P271	Use only outdoors or in a well-ventilated area.
P280	Wear protective gloves and eye protection.
P301+P310	IF SWALLOWED: Immediately call a POISON CENTER or physician.
P302+P352	IF ON SKIN: Wash with plenty of soap and water.
P303+P361+P353	IF ON SKIN (or hair): Take off immediately all contaminated clothing. Rinse skin with water.
P304+P340	IF INHALED: Remove person to fresh air and keep comfortable for breathing.
P305+P351+P338	IF IN EYES: Rinse cautiously with water for several minutes. Remove contact lenses, if present and easy to do. Continue rinsing.
P307+P311	IF exposed: Call a POISON CENTER or physician.
P308+P313	IF exposed or concerned: Get medical attention.
P310	Immediately call a POISON CENTER or physician.
P312	Call a POISON CENTER or physician if you feel unwell.
P314	Get medical attention if you feel unwell.
P321	Specific treatment (Wash areas of contact with water).
P330	Rinse mouth.
P332+P313	If skin irritation occurs: Get medical attention.
P362	Take off contaminated clothing and wash it before reuse.
P405	Store locked up.
P501	Dispose of contents in accordance with local, state, federal and international regulations.

2.3. WHMIS Classification

WHMIS classification is not included based on the recommended option (Option 4) found in the Canada Gazette Part II, Vol. 149, No.3, page 458



Safety Data Sheet

2.4. Hazards not Otherwise Classified or Covered by GHS

Data not available.

SECTION 3: Composition / Information on Ingredients

3.1. Components of Substance or Mixture

Chemical Name	Formula	Molecular Weight	CAS Number	Weight%
1-Butanol	CH ₃ CH ₂ CH ₂ CH ₂ OH	74.12 g/mol	71-36-3	44.25%
Water	H ₂ O	18.01 g/mol	7732-18-5	29.00%
Methyl Alcohol	CH ₃ OH	32.04 g/mol	67-56-1	25.24%
Calcium Chloride Dihydrate	CaCl ₂ ·2H ₂ O	147.02 g/mol	10035-04-8	1.52%

SECTION 4: First-Aid Measures

4.1. General First Aid Information

Eye Contact: IF IN EYES: Rinse cautiously with water for several minutes. Remove contact lenses, if present and easy to do. Continue rinsing. May cause irritation with burning and stinging with possible damage to the cornea and conjunctiva.

Inhalation: IF INHALED: Remove person to fresh air and keep comfortable for breathing.

Skin Contact: IF ON SKIN (or hair): Take off immediately all contaminated clothing. Rinse skin with water. Results in drying and cracking which can lead to secondary infections and dermatitis. Dermal absorption causes many of the symptoms of inhalation.

Ingestion: IF SWALLOWED: Immediately call a POISON CENTER or physician. Dilute immediately with water or milk. Induce vomiting. Call a physician.

4.2. Most Important Symptoms and Effects, Acute and Delayed

Flammable liquid. Toxic by ingestion, inhalation and dermal absorption. Cannot be made non-poisonous. If ingested, give large quantity of water and induce vomiting. Call a physician. Contact may cause dryness and cracking of the skin. May cause irritation to the eyes. Wash areas of contact with water. May cause irritation of the respiratory system. EYE CONTACT: May cause irritation with burning and stinging with possible damage to the cornea and conjunctiva. SKIN CONTACT: Results in drying and cracking which can lead to secondary infections and dermatitis. Dermal absorption causes many of the symptoms of inhalation.

4.3. Medical Attention or Special Treatment Needed

Immediately call a POISON CENTER or physician. Specific treatment (Wash areas of contact with water).

SECTION 5: Fire-Fighting Measures

5.1. Extinguishing Media

Carbon dioxide, dry chemical, alcohol foam

5.2. Specific Hazards Arising from the Substance or Mixture

Highly flammable liquid and vapor. Moderate explosion hazard and dangerous fire hazard when exposed to heat, sparks and open flames. Above flash point, vapor-air mixtures are explosive within flammable limits noted above. Sensitive to static discharge.

Safety Data Sheet

5.3. Special Protective Equipment for Firefighters

Wear full protective clothing and NIOSH-approved self-contained breathing apparatus with full facepiece operated in the pressure demand or other positive pressure mode.

SECTION 6: Accidental Release Measures

6.1. Personal Precautions, Protective Equipment and Emergency Procedures

Ground container and receiving equipment. Use explosion-proof equipment. Use only non-sparking tools. Take precautionary measures against static discharge. Wear protective gloves and eye protection.

6.2. Cleanup and Containment Methods and Materials

Remove all sources of ignition. Contain spill. Do not flush to sewer. Absorb with suitable inert material (vermiculite, dry sand, etc) and place in a chemical waste container for proper disposal in an approved waste disposal facility. Ventilate area of spill. Have extinguishing agent available in case of fire. Use non-sparking tools and equipment. Dispose of in accordance with local regulations.

SECTION 7: Handling and Storage

7.1. Precautions for Safe Handling and Storage Conditions

Store locked up. As with all chemicals, wash hands thoroughly after handling. Avoid contact with eyes and skin. Protect from freezing and physical damage. Store in secure, flammable storage area away from all sources of ignition. Empty containers may be hazardous since they retain product residues.

SECTION 8: Exposure Controls / Personal Protection

8.1. Control Parameters

Chemical Name	Limit Type	Country	Exposure Limit	Information Source
1-Butanol (71-36-3)	TWA	USA	100 ppm TWA 300 mg/m ³ TWA	U.S. - OSHA - Final PELs - Time Weighted Averages (TWAs)
1-Butanol (71-36-3)	TLV-TWA	USA	20 ppm TWA	ACGIH - Threshold Limit Values - Time Weighted Averages (TLV-TWA)
Methyl Alcohol (67-56-1)	TWA	USA	200 ppm TWA 260 mg/m ³ TWA	U.S. - OSHA - Final PELs - Time Weighted Averages (TWAs)
Methyl Alcohol (67-56-1)	TLV-STEL	USA	250 ppm STEL	ACGIH - Threshold Limit Values - Short Term Exposure Limits (TLV-STEL)
Methyl Alcohol (67-56-1)	TLV-TWA	USA	200 ppm TWA	ACGIH - Threshold Limit Values - Time Weighted Averages (TLV-TWA)

8.2. Exposure Controls

Engineering Controls: Use only outdoors or in a well-ventilated area. A system of local and/or general exhaust is recommended to keep employee exposures below the Airborne Exposure Limit.

Respiratory Protection: Normal room ventilation is adequate. If the exposure limit is exceeded, a full facepiece respirator with organic vapor cartridge may be worn.

Skin Protection: Wear protective gloves and eye protection. Chemical resistant gloves.



Safety Data Sheet

Eye Protection: Wear protective gloves and eye protection. Safety glasses or goggles.

8.3. Personal Protective Equipment

Wear protective gloves and eye protection. Normal room ventilation is adequate. If the exposure limit is exceeded, a full facepiece respirator with organic vapor cartridge may be worn. Chemical resistant gloves. Safety glasses or goggles.

SECTION 9: Physical and Chemical Properties

9.1. Basic Physical and Chemical Properties

Appearance: Colorless liquid

Physical State: Liquid

Odor: Data not available.

Odor Threshold: Data not available.

pH: Data not available.

Melting/Freezing Point: Data not available.

Initial Boiling Point /Range: Approximately 118°C - Approximately 118°C

Flash Point: Approximately 29°C

Evaporation Rate: Data not available.

Flammability: Data not available.

Flammability/Explosive Limits: Data not available.

Vapor Pressure: Data not available.

Vapor Density: Data not available.

Relative Density: 0.9

Solubility: Miscible

Partition Coefficient (n-Octanol/Water): Data not available.

Auto-Ignition Temperature: Data not available.

Decomposition Temperature: Data not available.

Viscosity: Data not available.

Explosive Properties: Data not available.

Oxidizing Properties: Data not available.

SECTION 10: Stability and Reactivity

10.1. Reactivity and Chemical Stability

Stable under normal conditions of use and storage.

10.2. Possibility of Hazardous Reactions

Data not available.



Safety Data Sheet

10.3. Conditions to Avoid and Incompatible Materials

Keep away from heat, sparks and open flame. No smoking. Keep container tightly closed. Strong oxidizing agents such as Nitrates, Perchlorates or

10.4. Hazardous Decomposition Products

Will not occur.

SECTION 11: Toxicological Information

11.1. Information on Toxicological Effects

Acute Toxicity - Oral Exposure:

Toxic if swallowed. Wash arms, hands and face thoroughly after handling. Do not eat, drink or smoke when using this product. IF SWALLOWED: Immediately call a POISON CENTER or physician. Specific treatment (Wash areas of contact with water). Rinse mouth. Store locked up. Dispose of contents in accordance with local, state, federal and international regulations.

Acute Toxicity - Dermal Exposure:

Harmful in contact with skin. Wear protective gloves and eye protection. IF ON SKIN: Wash with plenty of soap and water. Call a POISON CENTER or physician if you feel unwell. Specific treatment (Wash areas of contact with water). Take off contaminated clothing and wash it before reuse. Dispose of contents in accordance with local, state, federal and international regulations.

Acute Toxicity - Inhalation Exposure:

Harmful if inhaled. Avoid breathing dust, fumes or mist. Use only outdoors or in a well-ventilated area. IF INHALED: Remove person to fresh air and keep comfortable for breathing. Call a POISON CENTER or physician if you feel unwell.

Acute Toxicity - Other Information:

LD50, Oral, Rat: (Methanol) 5628 mg/kg, details of toxic effects not reported other than lethal dose value; (n-Butanol) 790 mg/kg, liver, kidney, blood effects noted.

Skin Corrosion and Irritation:

Causes skin irritation. Wash arms, hands and face thoroughly after handling. Wear protective gloves and eye protection. IF ON SKIN: Wash with plenty of soap and water. Specific treatment (Wash areas of contact with water). If skin irritation occurs: Get medical attention. Take off contaminated clothing and wash it before reuse.

Serious Eye Damage and Irritation:

Causes serious eye damage. Wear protective gloves and eye protection. IF IN EYES: Rinse cautiously with water for several minutes. Remove contact lenses, if present and easy to do. Continue rinsing. Immediately call a POISON CENTER or physician.

Respiratory Sensitization:

Not applicable.

Skin Sensitization:

Not applicable.

Germ Cell Mutagenicity:

Not applicable.

Carcinogenicity:

Not applicable.



Safety Data Sheet

Reproductive Toxicity:

May damage fertility or the unborn child. Obtain special instructions before use. Do not handle until all safety precautions have been read and understood. Wear protective gloves and eye protection. IF exposed or concerned: Get medical attention. Store locked up. Dispose of contents in accordance with local, state, federal and international regulations.

Specific Target Organ Toxicity from Single Exposure:

Causes damage to organs. Do not breathe dust, fumes or mist. Wash arms, hands and face thoroughly after handling. Do not eat, drink or smoke when using this product. IF exposed: Call a POISON CENTER or physician. Specific treatment (Wash areas of contact with water). Store locked up. Dispose of contents in accordance with local, state, federal and international regulations.

Specific Target Organ Toxicity from Repeated Exposure:

Causes damage to organs through prolonged or repeated exposure. Do not breathe dust, fumes or mist. Wash arms, hands and face thoroughly after handling. Do not eat, drink or smoke when using this product. Get medical attention if you feel unwell. Dispose of contents in accordance with local, state, federal and international regulations.

Aspiration Hazard:

Not applicable.

Additional Toxicology Information:

Data not available.

SECTION 12: Ecological Information

12.1. Ecotoxicity

Not applicable.

12.2. Persistence and Degradability

Data not available.

12.3. Bioaccumulative Potential

Data not available.

12.4. Mobility in Soil

Data not available.

12.5. Other Adverse Ecological Effects

Data not available.

SECTION 13: Disposal Considerations

13.1. Waste Treatment Methods

Data not available.

Safety Data Sheet

SECTION 14: Transportation Information

14.1. Transportation by Land - Department of Transportation (DOT, United States of America)

UN Number: UN1987

Proper Shipping Name: Alcohols, n.o.s., (Butanol and Methanol)

Hazard Class: 3

Packing Group: II

Hazard Placard Labels:



14.2. Transportation by Air - International Air Transport Association (IATA)

UN Number: UN1987

Proper Shipping Name: Alcohols, n.o.s., (Butanol and Methanol)

Hazard Class: 3

Packing Group: II

Hazard Placard Labels:



SECTION 15: Regulatory Information

15.1. Occupational Safety and Health Administration (OSHA) Hazards

Not listed.

15.2. Superfund Amendments and Reauthorization Act (SARA) 302 Extremely Hazardous Substances

Not listed.

15.3. Superfund Amendments and Reauthorization Act (SARA) 311/312 Hazardous Chemicals

Methyl Alcohol (CAS # 67-56-1): 5000 lb final RQ; 2270 kg final RQ

1-Butanol (CAS # 71-36-3): 5000 lb final RQ; 2270 kg final RQ

15.4. Superfund Amendments and Reauthorization Act (SARA) 313 Toxic Release Inventory (TRI)

Methyl Alcohol (CAS # 67-56-1): 1.0 % de minimis concentration

1-Butanol (CAS # 71-36-3): 1.0 % de minimis concentration

15.5. Massachusetts Right-to-Know Substance List

Methyl Alcohol (CAS # 67-56-1): Present

1-Butanol (CAS # 71-36-3): Present

Safety Data Sheet

15.6. Pennsylvania Right-to-Know Hazardous Substances

Methyl Alcohol (CAS # 67-56-1): Environmental hazard
Methyl Alcohol (CAS # 67-56-1): Present
1-Butanol (CAS # 71-36-3): Present
Water (CAS # 7732-18-5): Present

15.7. New Jersey Worker and Community Right-to-Know Components

Methyl Alcohol (CAS # 67-56-1): flammable - third degree; teratogen
Methyl Alcohol (CAS # 67-56-1): sn 1222
Methyl Alcohol (CAS # 67-56-1): SN 1222 TPQ: 500 lb
1-Butanol (CAS # 71-36-3): flammable - third degree
1-Butanol (CAS # 71-36-3): sn 1330
1-Butanol (CAS # 71-36-3): SN 1330 TPQ: 500 lb

15.8. California Proposition 65

Methyl Alcohol (CAS # 67-56-1): developmental toxicity, initial date 3/16/12

15.9. Canada Domestic Substances List / Non-Domestic Substances List (DSL/NDSL)

Calcium Chloride Dihydrate (CAS # 10035-04-8): Present
Methyl Alcohol (CAS # 67-56-1): Present
1-Butanol (CAS # 71-36-3): Present
Water (CAS # 7732-18-5): Present

15.10. United States of America Toxic Substances Control Act (TSCA) List

Calcium Chloride Dihydrate (CAS # 10035-04-8): Present
Methyl Alcohol (CAS # 67-56-1): Present
1-Butanol (CAS # 71-36-3): Present
Water (CAS # 7732-18-5): Present

15.11. European Inventory of Existing Commercial Chemical Substances (EINECS), European List of Notified Chemical Substances (ELINCS), and No Longer Polymers (NLP)

Not listed.

Safety Data Sheet

SECTION 16: Other Information

16.1. Full Text of Hazard Statements and Precautionary Statements

Highly flammable liquid and vapor. Toxic if swallowed. Harmful in contact with skin. Causes skin irritation. Causes serious eye damage. Harmful if inhaled. May damage fertility or the unborn child. Causes damage to organs. Causes damage to organs through prolonged or repeated exposure.

Obtain special instructions before use. Do not handle until all safety precautions have been read and understood. Keep away from heat, sparks and open flame. No smoking. Keep container tightly closed. Ground container and receiving equipment. Use explosion-proof equipment. Use only non-sparking tools. Take precautionary measures against static discharge. Do not breathe dust, fumes or mist. Wash arms, hands and face thoroughly after handling. Do not eat, drink or smoke when using this product. Use only outdoors or in a well-ventilated area. Wear protective gloves and eye protection.

IF SWALLOWED: Immediately call a POISON CENTER or physician. IF ON SKIN: Wash with plenty of soap and water. IF ON SKIN (or hair): Take off immediately all contaminated clothing. Rinse skin with water. IF INHALED: Remove person to fresh air and keep comfortable for breathing. IF IN EYES: Rinse cautiously with water for several minutes. Remove contact lenses, if present and easy to do. Continue rinsing. IF exposed: Call a POISON CENTER or physician. Get medical attention if you feel unwell. Specific treatment (Wash areas of contact with water). Rinse mouth. If skin irritation occurs: Get medical attention. Take off contaminated clothing and wash it before reuse. In case of fire: Use dry chemical, foam or carbon dioxide to extinguish.

Store in a well-ventilated place. Keep cool. Store locked up.

Dispose of contents in accordance with local, state, federal and international regulations.

16.2. Miscellaneous Hazard Classes

Canadian Carcinogenicity Hazard Class: Not Applicable.


Physical Hazards Not Otherwise Classified (PHNOC): Not Applicable.

Health Hazards Not Otherwise Classified (HHNOC): Not Applicable.

Not Applicable.

16.3. National Fire Protection Association (NFPA) Rating

Health: 1
Flammability: 3
Reactivity: 0
Special Hazard:

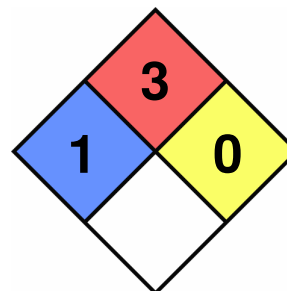

Safety Data Sheet

16.4. Document Revision

Last Revision Date: 5/1/2015

DISCLAIMER

When handled properly by qualified personnel, the product described herein does not present a significant health or safety hazard. Alteration of its characteristics by concentration, evaporation, addition of other substances, or other means may present hazards not specifically addressed herein and which must be evaluated by the user. The information furnished herein is believed to be accurate and represents the best data currently available to us. No warranty, expressed or implied, is made and RICCA CHEMICAL COMPANY assumes no legal responsibility or liability whatsoever resulting from its use.