



Safe Water?

Teacher Information

.....just add students™

Summary

An infant is admitted to the hospital with serious symptoms. The infant's family recently moved to a house with a private well, not a public water supply. Students test water from the family's well and identify health risks associated with the contaminated well water. The family moves to a house with a public water supply and must consider whether to use tap water or bottled water.

Core Concepts

- Well water should be tested regularly and protected from contamination from microorganisms and harmful chemicals.
- Public water supplies are tested and regulated by the Environmental Protection Agency (EPA) to provide clean and safe tap water.
- People should consider safety, costs, health consequences, and ecological impact of using tap water versus bottled water.

Time Required

Two 40-minute class periods

Kit contains

- Simulated samples of well water, tap water, and bottled water
- Nitrate Test Solution (simulated)
- Arsenic Test Solution (simulated)
- Coliform Bacteria Test solution (simulated)
- Labeled droppers
- *Well Water Testing* sheet
- *Bottled Water Testing* sheet
- *Tap Water Testing* sheet
- *Water Testing Color Charts*
- *Fact Sheet: Well Water Safety*
- *The Facts about Bottled Water*

Teacher Provides

- Safety goggles
- Paper towel for clean-up

Warning: Choking Hazard

This Science Take-Out kit contains small parts. Do not allow children under the age of seven to have access to any kit component.

Teacher Suggestions

- If students are working in teams, you may make copies of handouts for each student.
- The EPA publication **Water on Tap** provides detailed information on water use and water safety. <http://nepis.epa.gov/Exe/ZyPDF.cgi/P1008ZP0.PDF?Dockey=P1008ZP0.PDF>

Reusing *Safe Water?* Kits

Kits may be refilled and reused. Allow approximately 30 minutes for refilling 10 student kits. Teachers will need to instruct students on how to handle clean-up and return of the reusable kit materials. For example, teachers might provide the following information for students:

Rinse with water and dry with paper towel	Refill and Return to kit
<ul style="list-style-type: none">• Water Test Sheet	<ul style="list-style-type: none">• All tubes and droppers• <i>Well Water Test Sheet</i>• <i>Fact Sheet: Well Water Safety</i>• <i>The Facts About Bottled Water</i>




Note: It is not necessary to rinse or wash the droppers after use. Because the droppers are labeled, there is little chance for contamination. Washing the droppers may make the labels difficult to read. Simply ask students to squirt out any extra liquid from the droppers.

Consider laminating printed parts of the kits (such as information sheets) that will be reused.

Refills for the *Safe Water?* kits are available at www.sciencetakeout.com. The 10 Kit Refill Pack includes the following materials:

- Instructions and Quick Guide for refilling kits
- 5 mL simulated Well Water Sample
- 5 mL simulated Bottled Water Sample
- 5 mL simulated Tap Water Sample
- 5 mL simulated Coliform Bacteria Test Solution
- 5 mL simulated Nitrate Test Solution
- 5 mL simulated Arsenic Test Solution

Kit Contents Quick Guide

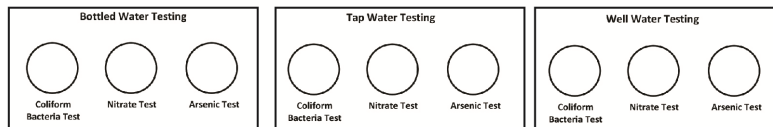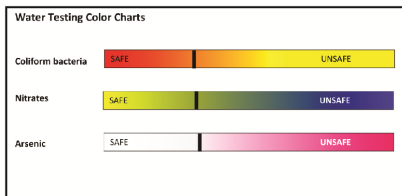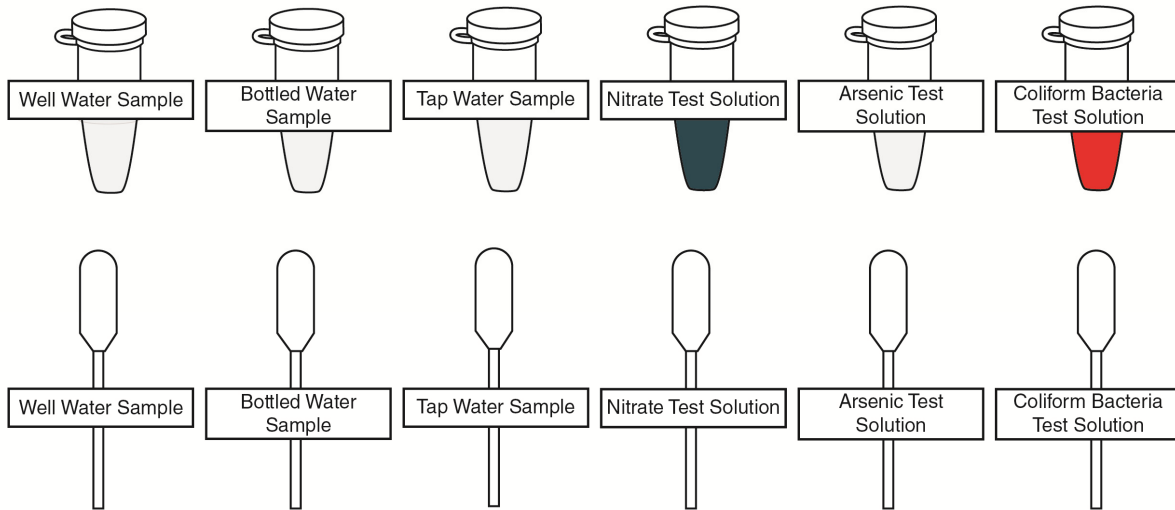

Fact Sheet: Well Water Safety

The Environmental Protection Agency (EPA) regulates public water systems but it does not have the authority to regulate private drinking water wells. Approximately 25 percent of Americans rely on their own private wells for drinking water and many do not have experts regularly checking the water's source and its quality. People who use private wells as their water supply must take special precautions to protect their well water from contamination.

If you use well water, you should have it tested yearly by a state certified water testing laboratory. You should test well water more frequently if someone in your house is pregnant or nursing, if there are unexplained illnesses in your family, if your neighbors find a contaminant in their well water, or if you note a change in water taste, odor, color, or clarity. If contamination is found, use bottled water for drinking and cooking until the well water supply is treated and tested to show that it is safe.

Sources of Well Water Contamination

Coliform Bacteria:
Coliform bacteria (such as *E. coli*) normally live in the intestines of people and animals. Most types of coliform bacteria are harmless. However, some types of coliform bacteria are harmful because they cause intestinal infections. Symptoms of intestinal infection caused by coliform bacteria include diarrhea, abdominal pain, and fever. More severe cases can lead to bloody diarrhea, dehydration, or kidney failure. Most intestinal bacterial infections are caused by contaminated food or water. Intestinal bacterial infections are usually treated with antibiotics and fluids to replace the fluids lost through vomiting and diarrhea.

Coliform bacteria are also found in human and animal wastes. Water flowing over or under land surfaces can pick up coliform bacteria pollutants from wildlife and soil, especially if flooding occurs. Bacterial pollution can be caused by drainage from septic tanks and by drainage from farm animals' manure, wild animals, landfills, and garbage dumps.

Nitrates:
Nitrates come from the breakdown of nitrogen compounds in the soil. Flowing ground water picks up nitrates from the soil. Nitrates are found in human wastes and animal wastes, so care must be taken to locate wells where they cannot be polluted by septic tanks and manure from farm animals. Fields that are fertilized with nitrate fertilizers, landfills, and garbage dumps are also sources of nitrate pollution.

High nitrate levels in water can cause methemoglobinemia, a serious condition found especially in infants less than six months of age. Nitrate converts hemoglobin (oxygen-carrying component of blood) into methemoglobin, which does not carry oxygen efficiently. This results in a reduced oxygen supply to all parts of the body. The most obvious symptom of methemoglobinemia is a bluish color of the skin. Other symptoms include headache, dizziness, weakness or difficulty in breathing. It is important to take babies with these symptoms to the hospital emergency room immediately. If methemoglobinemia is recognized in time, it can be treated with an injection of a chemical called methylene blue. If methemoglobinemia is not treated, it can result in brain damage and death.

Arsenic:
Arsenic interferes with enzymes essential for cellular respiration and decreases the production of ATP energy needed for life activities. Symptoms of arsenic poisoning begin with headaches, confusion, diarrhea, and drowsiness. Symptoms of arsenic poisoning often appear in children at doses of arsenic that are too low to have obvious effects on adults. When arsenic poisoning becomes more severe, patients may experience severe diarrhea, vomiting, blood in the urine, cramping muscles, hair loss, stomach pain, and convulsions. The final result of severe arsenic poisoning is coma and death. Treatment for arsenic poisoning involves the use of drugs that combine with arsenic and prevent it from interacting with cellular enzymes.

Arsenic is a naturally occurring element found throughout the earth's crust. Mining and construction can release large amounts of arsenic into nearby ground water sources. Some older fruit orchards may contain high levels of arsenic, since it was once used as a pesticide. These are common sources of well water arsenic contamination. Underground rocks and soils may contain arsenic; however, these natural sources usually do not contaminate wells at dangerous levels.

THE FACTS ABOUT BOTTLED WATER

Bottled water may be safer than tap water if your water supply is a private well or other water source that contains contaminants.

tested for <i>e.coli</i>	YES	NO
required to provide source	YES	NO
required to produce quality reports	YES	NO
can still be distributed when tap water standards aren't met	NO	YES

Bottling water produces 2.5 million tons of carbon dioxide yearly.

It takes 3x the amount of water to produce the bottles as it does to fill it.

40% of bottled water comes from public tap water supplies.

17 Million barrels of oil are used in the production of water bottles yearly.

\$10.00 per gallon vs \$0.0015 per gallon

Only 20% of bottles are recycled.

Bottles that are not recycled end up in landfills or end up as litter on land, in waterways, and oceans.

It takes 450 years for a plastic bottle to decompose.

Plastic bottles can leach potentially harmful chemicals such as BPA and phthalates into water in the bottle or into the environment.

Read these instructions before using Science Take-Out kits

Parental or Adult Supervision Required

This kit should be used only under the supervision of an adult who is committed to ensuring that the safety precautions below, and in the specific laboratory activity, are followed.

Safety Goggles and Gloves Strongly Recommended

We encourage students to adopt safe lab practices, and wear safety goggles and gloves when performing laboratory activities involving chemicals. Safety goggles and gloves are not provided in Science Take-Out kits. They may be purchased from a local hardware store or pharmacy.

Warning: Choking and Chemical Hazard

Science Take-Out kits contain small parts that could pose a choking hazard and chemicals that could be hazardous if ingested. Do not allow children under the age of seven to have access to any kit components. Safety Data Sheets (SDS) provide specific safety information regarding the chemical contents of the kits. SDS information for each kit is provided in the accompanying teacher instructions.

Chemicals Used in Science Take-Out Kits

Every effort has been made to reduce the use of hazardous chemicals in Science Take-Out kits. Most kits contain common household chemicals or chemicals that pose little or no risk.

General Safety Precautions

1. Work in a clean, uncluttered area. Cover the work area to protect the work surface.
2. Read and follow all instructions carefully.
3. Pay particular attention to following the specific safety precautions included in the kit activity instructions.
4. Goggles and gloves should be worn while performing experiments using chemicals.
5. Do not use the contents of this kit for any other purpose beyond those described in the kit instructions.
6. Do not leave experiment parts or kits where they could be used inappropriately by others.
7. Never taste or ingest any chemicals provided in the kit – they may be toxic.
8. Do not eat, drink, or apply make-up or contact lenses while performing experiments.
9. Wash your hands before and after performing experiments.
10. Chemicals used in Science Take-Out experiments may stain or damage skin, clothing or work surfaces. If spills occur, wash the area immediately and thoroughly.
11. At the end of the experiment, return ALL kit components to the kit plastic bag. Dispose of the plastic bag and contents in your regular household trash.

No blood or body fluids from humans or animals are used in Science Take-Out kits. Chemical mixtures are substituted as simulations of these substances.

Safe Water?

Part 1: Case of an Ailing Infant


Patient: Joey Nash

Age: 5 months

Symptoms: Vomiting, diarrhea, bluish skin, weakness, and increased breathing rate

Patient History: Joey's parents report that he was healthy until they moved several weeks ago into a home near a small farm. The water that the Nash family uses for drinking, cooking, and for preparing Joey's baby formula comes from their private well, not the town's public water supply.

Recommendation: Well water testing should be done because several families who live near the Nash family have had health problems due to contaminated well water.

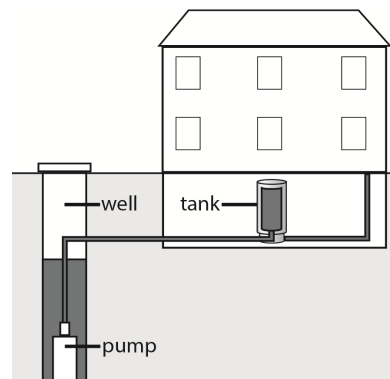

You will test the water from the Nash's well for three potential contaminants that could be causing Joey's symptoms—coliform bacteria, nitrates, and arsenic.

1. Place one drop of the well water sample into each of the circles on the **Well Water Testing** sheet.
2. Add one drop of each of the appropriate test solutions to the Well Water Testing circles:
 - a. Add one drop of Coliform Bacteria Test Solution to the Coliform Bacteria Test circle.
 - b. Add one drop of Nitrate Test Solution to the Nitrate Test circle.
 - c. Add one drop of Arsenic Test Solution to the Arsenic Test circle.

3. Use the **Water Testing Color Charts** in your kit to determine whether the concentration of each contaminant is safe or unsafe (too high). Record the results of the water tests in the data table below.

Potential Well Water Contaminant	Safe or Unsafe
Coliform Bacteria	
Nitrates	
Arsenic	

4. Based on the results of the well water tests, what conclusions can you draw about the safety of the water in the Nash's well?

Base your answers to questions 5 through 9 on the information in the **Fact Sheet: Well Water Safety**.

5. Which of Joey's symptoms may be caused by unsafe nitrate levels in the well water?

6. Name the condition caused by unsafe nitrate levels and the treatment used for this condition.

7. Which of Joey's symptoms may be caused by unsafe coliform bacteria levels in the well water?

8. Name the condition caused by unsafe coliform levels and the treatment used for this condition.

9. List at least three likely sources of the contaminants that made Joey sick? Be specific.

- ---
- ---
- ---

10. List at least three actions that Joey's family could take to prevent well water contamination from affecting their health in the future.

- ---
- ---
- ---

Part 2: Tap Water versus Bottled Water

Joey recovered from his illnesses, and for a while the Nash family used bottled water for all of their drinking and cooking. However, bottled water was expensive and their landlord refused to take action to treat their water supply. The Nash family moved to an apartment in town so that they would have tap water from the town's public water supply. The town water supply came from a local lake and was treated and tested before it was distributed to the town.

Town officials gave Mrs. Nash a copy of an official report that said the water had been tested and was safe. To be sure the town's water was safe, Mrs. Nash sent a sample of the town tap water to a state certified water testing laboratory. For comparison, she also sent a sample of "Pure Spring Water" sold by a local bottled water company.

You will test samples of the bottled water and the town tap water for three potential contaminants—coliform bacteria, nitrates, and arsenic.

1. Place one drop of the bottled water sample into each of the Bottled Water Testing circles on the **Bottled Water Testing** sheet.
2. Place one drop of the tap water sample into each of the Tap Water Testing circles on the **Tap Water Testing** sheet.
3. Add one drop of each of the appropriate test solutions to the Bottled Water Testing circles and the Tap Water Testing circles:
 - a. Add one drop of Coliform Bacteria Test Solution to the Coliform Bacteria Test circles.
 - b. Add one drop of Nitrate Test Solution to the Nitrate Test circles.
 - c. Add one drop of Arsenic Test Solution to the Arsenic Test circles.

4. Use the **Water Testing Color Charts** in your kit to determine whether the concentration of each contaminant is safe or unsafe (too high). Record the results of the water tests in the data table below.

Potential Water Contaminant	Bottled Water Safe or Unsafe	Tap Water Safe or Unsafe
Coliform Bacteria		
Nitrates		
Arsenic		

5. Based on the results of the bottled water and tap water tests, what conclusions can you draw about the safety of each type of water?

6. Mrs. Nash would like to continue using bottled water for drinking and cooking. List one possible reason why she might prefer to use bottled water instead of tap water.

7. Mr. Nash would prefer using tap water for drinking and cooking. List one possible reason why he might prefer to use tap water instead of bottled water.

Section 1 Chemical Product and Company Information

Science Take-Out 80 Office Park Way
Pittsford, NY 14534
(585)764-5400

**CHEMTREC 24 Hour Emergency
Phone Number (800) 424-9300**

For laboratory use only. Not for drug, food or household use

Product	Buffer Solution pH10
Synonyms	"Well Water Sample" (simulated)

Section 2 Hazards Identification

This substance or mixture has not been classified at this time according to the Globally Harmonized System (GHS) of Classification and Labeling of Chemicals.

Signal word: WARNING
Pictograms: None required
Target organs: None known

GHS Classification:
Skin irritation (Category 3)
Eye irritation (Category 2B)

GHS Label information: Hazard statement(s):
H316: Causes mild skin irritation.
H320: Causes eye irritation.

Precautionary statement(s):

P264: Wash hands thoroughly after handling.

P305+P351+P338: IF IN EYES: Rinse cautiously with water for several minutes. Remove contact lenses, if present and easy to do. Continue rinsing.

P332+P313: If skin irritation occurs: Get medical attention.

P337+P313: If eye irritation persists: Get medical attention.

Ca Prop 65 - This product does not contain any chemicals known to the State of California to cause cancer, birth defects, or any other reproductive harm.

Section 3 Composition / Information on Ingredients

Chemical Name	CAS #	%	EINECS
Water	7732-18-5	99.08%	231-791-2
Potassium chloride	7447-40-7	0.40%	231-211-8
Boric acid	10043-35-3	0.33%	233-139-2
Sodium hydroxide	1310-73-2	0.19%	215-185-5

Section 4 First Aid Measures

INGESTION: Call physician or Poison Control Center immediately. Induce vomiting only if advised by appropriate medical personnel. Never give anything by mouth to an unconscious person.

INHALATION: Remove to fresh air. If not breathing, give artificial respiration. If breathing is difficult, give oxygen. Get medical attention.

EYE CONTACT: Check for and remove contact lenses. Flush thoroughly with water for at least 15 minutes, lifting upper and lower eyelids occasionally. Get immediate medical attention.

SKIN ABSORPTION: Remove contaminated clothing. Flush thoroughly with mild soap and water. If irritation occurs, get medical attention.

Section 5 Fire Fighting Measures

Suitable Extinguishing Media: Use any media suitable for extinguishing supporting fire.

Protective Actions for Fire-fighters: In fire conditions, wear a NIOSH/MSHA-approved self-contained breathing apparatus and full protective gear. Use water spray to keep fire-exposed containers cool.

Specific Hazards: During a fire, irritating and highly toxic gases may be generated by thermal decomposition or combustion.

Section 6 Accidental Release Measures

Personal Precautions: Evacuate personnel to safe area. Use proper personal protective equipment as indicated in Section 8. Provide adequate ventilation.

Environmental Precautions: Avoid runoff into storm sewers and ditches which lead to waterways.

Containment and Cleanup: Absorb with inert dry material, sweep or vacuum up and place in a suitable container for proper disposal. Wash spill area with soap and water.

Section 7 Handling and Storage

Precautions for Safe Handling: Read label on container before using. Do not wear contact lenses when working with chemicals. Keep out of reach of children. Avoid contact with eyes, skin and clothing. Do not inhale vapors, spray or mist. Use with adequate ventilation. Avoid ingestion. Wash thoroughly after handling. Remove and wash clothing before reuse.

Conditions for Safe Storage: Store in a cool, well-ventilated area away from incompatible substances.

Section 8 Exposure controls / Personal Protection

Exposure Limits:	Chemical Name	ACGIH (TLV)	OSHA (PEL)	NIOSH (REL)
	Potassium chloride	None established	None established	None established

Engineering controls: Facilities storing or utilizing this material should be equipped with an eyewash facility and a safety shower and fire extinguishing material. Personnel should wear safety glasses, goggles, or faceshield, lab coat or apron, appropriate protective gloves. Use adequate ventilation to keep airborne concentrations low.

Respiratory protection: None should be needed in normal laboratory handling at room temperatures. If misty conditions prevail, work in fume hood or wear a NIOSH/MSHA approved respirator.

Section 9 Physical and Chemical Properties

Appearance: Clear, colorless liquid. Odor: No odor. Odor threshold: Data not available. pH: 10.0 Melting/Freezing point: Approx. 0°C (32°F) (water) Boiling point: Approx. 100°C (212°F) (water) Flash point: Data not available	Evaporation rate (Water = 1): <1 Flammability (solid/gas): Data not available. Explosion limits: Lower/Upper: Data not available Vapor pressure (mm Hg): 14 (water) Vapor density (Air = 1): 0.7 (water) Relative density (Specific gravity): Approx. 1.0 (water) Solubility(ies): Complete in water.	Partition coefficient: Data not available Auto-ignition temp.: Data not available Decomposition temp.: Data not available Viscosity: Data not available. Molecular formula: Mixture Molecular weight: Mixture
---	--	--

Section 10 Stability and Reactivity

Chemical stability: Stable

Hazardous polymerization: Will not occur.

Conditions to avoid: Excessive temperatures which cause evaporation.

Incompatibilities with other materials: Acids, alkalis, and air will change the buffer's ability.

Hazardous decomposition products: Boron oxide and chlorine gas.

Section 11 Toxicological Information

Acute toxicity: Data not available

Serious eye damage/irritation: Data not available

Germ cell mutagenicity: Data not available

Skin corrosion/irritation: Data not available

Respiratory or skin sensitization: Data not available

Carcinogenicity: Data not available

NTP: No component of this product present at levels greater than or equal to 0.1% is identified as a known or anticipated carcinogen by NTP.

IARC: No component of this product present at levels greater than or equal to 0.1% is identified as probable, possible or confirmed human carcinogen by IARC.

OSHA: No component of this product present at levels greater than or equal to 0.1% is identified as a carcinogen or potential carcinogen by OSHA.

Reproductive toxicity: Data not available

STOT-single exposure: Data not available

Aspiration hazard: Data not available

STOT-repeated exposure: Data not available

Potential health effects:

Inhalation: May be harmful if inhaled.

Ingestion: May be harmful if swallowed.

Skin: May cause mild irritation.

Eyes: May cause mild irritation.

Signs and symptoms of exposure: To the best of our knowledge the chemical, physical and toxicological properties have not been thoroughly investigated. Specific data is not available. Exercise appropriate procedures to minimize potential hazards.

Additional information: RTECS #: Data not available

Section 12 Ecological Information

Toxicity to fish: No data available

Toxicity to daphnia and other aquatic invertebrates: No data available

Toxicity to algae: No data available

Persistence and degradability: No data available

Bioaccumulative potential: No data available

Mobility in soil: No data available

PBT and vPvB assessment: No data available

Other adverse effects: An environmental hazard cannot be excluded in the event of unprofessional handling or disposal.

Section 13 Disposal Considerations

These disposal guidelines are intended for the disposal of catalog-size quantities only. Federal regulations may apply to empty container. State and/or local regulations may be different. Dispose of in accordance with all local, state and federal regulations or contract with a licensed chemical disposal agency.

Section 14 Transport Information

UN/NA number: Not applicable

Shipping name: Not Regulated

Hazard class: Not applicable

Packing group: Not applicable

Reportable Quantity: No

Marine pollutant: No

Exceptions: Not applicable

2012 ERG Guide # Not applicable

Section 15 Regulatory Information

A chemical is considered to be listed if the CAS number for the anhydrous form is on the Inventory list.

Component	TSCA	CERLCA (RQ)	RCRA code	DSL	NDSL	WHMIS Classification
Potassium Chloride	Listed	Not Listed	Not Listed	Listed	Not Listed	Uncontrolled Product
Sodium hydroxide	Listed	1,000 lbs (454 kg)	D002	Listed	Not Listed	E

Section 16 Additional Information

The information contained herein is furnished without warranty of any kind. Employers should use this information only as a supplement to other information gathered by them and must make independent determinations of suitability and completeness of information from all sources to assure proper use of these materials and the safety and health of employees.

NTP: National Toxicology Program, IARC: International Agency for Research on Cancer, OSHA: Occupational Safety and Health Administration, STOT: Specific Target Organ Toxicity, SE: Single Exposure, RE: Repeated Exposure, ERG: Emergency Response Guidebook.

Section 1 Chemical Product and Company Information

Science Take-Out 80 Office Park Way
Pittsford, NY 14534
(585)764-5400

**CHEMTREC 24 Hour Emergency
Phone Number (800) 424-9300**
For laboratory use only. Not for drug, food or household use

Product	Buffer Solution pH3
Synonyms	"Bottled Water Sample" (simulated), "Tap Water Sample" (simulated)

Section 2 Hazards Identification

This substance or mixture has not been classified at this time according to the Globally Harmonized System (GHS) of Classification and Labeling of Chemicals.

Signal word: WARNING
Pictograms: None required
Target organs: None known

GHS Classification:
Skin Irritation (Category 3)
Eye irritation (Category 2B)

GHS Label information: Hazard statement(s):
H316: Causes minor skin irritation.
H320: Causes eye irritation.

Precautionary statement(s):

P264: Wash hands thoroughly after handling.

P305+P351+P338: IF IN EYES: Rinse cautiously with water for several minutes. Remove contact lenses, if present and easy to do. Continue rinsing.

P332+P313: If skin irritation occurs: Get medical attention.

P337+P313: If eye irritation persists: Get medical attention.

Ca Prop 65 - This product does not contain any chemicals known to the State of California to cause cancer, birth defects, or any other reproductive harm.

Section 3 Composition / Information on Ingredients

Chemical Name	CAS #	%	EINECS
Water	7732-18-5	99.74%	231-791-2
Potassium biphthalate	877-24-7	1.08%	212-889-4
Hydrochloric acid	7647-01-0	0.18%	231-595-7

Section 4 First Aid Measures

INGESTION: Call physician or Poison Control Center immediately. Induce vomiting only if advised by appropriate medical personnel. Never give anything by mouth to an unconscious person.

INHALATION: Remove to fresh air. If not breathing, give artificial respiration. If breathing is difficult, give oxygen. Get medical attention.

EYE CONTACT: Check for and remove contact lenses. Flush thoroughly with water for at least 15 minutes, lifting upper and lower eyelids occasionally. Get immediate medical attention.

SKIN ABSORPTION: Remove contaminated clothing. Flush thoroughly with mild soap and water. If irritation occurs, get medical attention.

Section 5 Fire Fighting Measures

Suitable Extinguishing Media: Use any media suitable for extinguishing supporting fire.

Protective Actions for Fire-fighters: In fire conditions, wear a NIOSH/MSHA-approved self-contained breathing apparatus and full protective gear. Use water spray to keep fire-exposed containers cool.

Specific Hazards: During a fire, irritating and highly toxic gases may be generated by thermal decomposition or combustion.

Section 6 Accidental Release Measures

Personal Precautions: Evacuate personnel to safe area. Use proper personal protective equipment as indicated in Section 8. Provide adequate ventilation.

Environmental Precautions: Avoid runoff into storm sewers and ditches which lead to waterways.

Containment and Cleanup: Absorb with inert dry material, sweep or vacuum up and place in a suitable container for proper disposal. Wash spill area with soap and water.

Section 7 Handling and Storage

Precautions for Safe Handling: Read label on container before using. Do not wear contact lenses when working with chemicals. Keep out of reach of children. Avoid contact with eyes, skin and clothing. Do not inhale vapors, spray or mist. Use with adequate ventilation. Avoid ingestion. Wash thoroughly after handling. Remove and wash clothing before reuse.

Conditions for Safe Storage: Store in a cool, well-ventilated area away from incompatible substances.

Section 8 Exposure controls / Personal Protection

Exposure Limits:	Chemical Name	ACGIH (TLV)	OSHA (PEL)	NIOSH (REL)
	None	None established	None established	None established

Engineering controls: Facilities storing or utilizing this material should be equipped with an eyewash facility and a safety shower and fire extinguishing material. Personnel should wear safety glasses, goggles, or faceshield, lab coat or apron, appropriate protective gloves. Use adequate ventilation to keep airborne concentrations low.

Respiratory protection: None should be needed in normal laboratory handling at room temperatures. If misty conditions prevail, work in fume hood or wear a NIOSH/MSHA approved respirator.

Section 9 Physical and Chemical Properties

Appearance: Clear, colorless liquid. Odor: No odor. Odor threshold: Data not available. pH: 3.0 Melting/Freezing point: Approx. 0°C (32°F) (water) Boiling point: Approx. 100°C (212°F) (water) Flash point: Data not available	Evaporation rate (Water = 1): <1 Flammability (solid/gas): Data not available. Explosion limits: Lower/Upper: Data not available Vapor pressure (mm Hg): 14 (water) Vapor density (Air = 1): 0.7 (water) Relative density (Specific gravity): Approx. 1.0 (water) Solubility(ies): Complete in water.	Partition coefficient: Data not available Auto-ignition temp.: Data not available Decomposition temp.: Data not available Viscosity: Data not available. Molecular formula: Mixture Molecular weight: Mixture
--	--	--

Section 10 Stability and Reactivity

Chemical stability: Stable

Hazardous polymerization: Will not occur.

Conditions to avoid: Excessive temperatures which cause evaporation.

Incompatibilities with other materials: Acids, alkalis, and air will change the buffer's ability.

Hazardous decomposition products: Carbon oxides and chlorine gas.

Section 11 Toxicological Information

Acute toxicity: Data not available

Serious eye damage/irritation: Data not available

Germ cell mutagenicity: Data not available

Skin corrosion/irritation: Data not available

Respiratory or skin sensitization: Data not available

Carcinogenicity: Data not available

NTP: No component of this product present at levels greater than or equal to 0.1% is identified as a known or anticipated carcinogen by NTP.

IARC: No component of this product present at levels greater than or equal to 0.1% is identified as probable, possible or confirmed human carcinogen by IARC.

OSHA: No component of this product present at levels greater than or equal to 0.1% is identified as a carcinogen or potential carcinogen by OSHA.

Reproductive toxicity: Data not available

Aspiration hazard: Data not available

STOT-single exposure: Data not available

STOT-repeated exposure: Data not available

Potential health effects:

Inhalation: May be harmful if inhaled.

Ingestion: May be harmful if swallowed.

Skin: May cause mild irritation.

Eyes: May cause mild irritation.

Signs and symptoms of exposure: To the best of our knowledge the chemical, physical and toxicological properties have not been thoroughly investigated. Specific data is not available. Exercise appropriate procedures to minimize potential hazards.

Additional information: RTECS #: Data not available

Section 12 Ecological Information

Toxicity to fish: No data available

Toxicity to daphnia and other aquatic invertebrates: No data available

Toxicity to algae: No data available

Persistence and degradability: No data available

Bioaccumulative potential: No data available

Mobility in soil: No data available

PBT and vPvB assessment: No data available

Other adverse effects: An environmental hazard cannot be excluded in the event of unprofessional handling or disposal.

Section 13 Disposal Considerations

These disposal guidelines are intended for the disposal of catalog-size quantities only. Federal regulations may apply to empty container. State and/or local regulations may be different. Dispose of in accordance with all local, state and federal regulations or contract with a licensed chemical disposal agency.

Section 14 Transport Information

UN/NA number: Not applicable

Hazard class: Not applicable

Exceptions: Not applicable

Shipping name: Not Regulated

Packing group: Not applicable

2012 ERG Guide # Not applicable

Reportable Quantity: No

Marine pollutant: No

Section 15 Regulatory Information

A chemical is considered to be listed if the CAS number for the anhydrous form is on the Inventory list.

Component	TSCA	CERLCA (RQ)	RCRA code	DSL	NDSL	WHMIS Classification
Potassium bipthalate	Listed	Not Listed	Not Listed	Listed	Not Listed	Not listed
Hydrochloric acid	Listed	5,000 lbs (2270 kg)	Not Listed	Listed	Not Listed	E

Section 16 Additional Information

The information contained herein is furnished without warranty of any kind. Employers should use this information only as a supplement to other information gathered by them and must make independent determinations of suitability and completeness of information from all sources to assure proper use of these materials and the safety and health of employees.

NTP: National Toxicology Program, IARC: International Agency for Research on Cancer, OSHA: Occupational Safety and Health Administration, STOT: Specific Target Organ Toxicity, SE: Single Exposure, RE: Repeated Exposure, ERG: Emergency Response Guidebook.

Section 1 Chemical Product and Company Information

Science Take-Out
80 Office Park Way
Pittsford, NY 14534
(585)764-5400

**CHEMTREC 24 Hour Emergency
Phone Number (800) 424-9300**
For laboratory use only. Not for drug, food or household use

Product	Bromothymol blue, 0.1% aqueous solution
Synonyms	"Nitrate Test Solution" (simulated)

Section 2 Hazards Identification

This substance or mixture has not been classified at this time according to the Globally Harmonized System (GHS) of Classification and Labeling of Chemicals.

Signal word: Not classified
Pictograms: Not classified
Target organs: None known

GHS Classification: Not classified
GHS Label information: Not classified
Precautionary Statement: Not classified

Supplementary information:

Do not breathe vapors, spray or mist. Do not get in eyes, on skin, or on clothing. Wear protective gloves/protective clothing/eye protection/face protection. Wash hands thoroughly after handling. Get medical attention if you feel unwell.

Ca Prop 65 - This product does not contain any chemicals known to the State of California to cause cancer, birth defects, or any other reproductive harm.

Section 3 Composition / Information on Ingredients

Chemical Name	CAS #	%	EINECS
Water	7732-18-5	99.9%	231-791-2
Bromothymol blue, sodium salt	34722-90-2	0.1%	252-169-7

Section 4 First Aid Measures

INGESTION: MAY BE HARMFUL IF SWALLOWED. Call physician or Poison Control Center immediately. Induce vomiting only if advised by appropriate medical personnel. Never give anything by mouth to an unconscious person.

INHALATION: MAY BE HARMFUL IF INHALED. MAY CAUSE RESPIRATORY TRACT IRRITATION. Remove to fresh air. If not breathing, give artificial respiration. If breathing is difficult, give oxygen. Get medical attention.

EYE CONTACT: MAY CAUSE EYE IRRITATION. Check for and remove contact lenses. Flush thoroughly with water for at least 15 minutes, lifting upper and lower eyelids occasionally. Get immediate medical attention.

SKIN ABSORPTION: MAY CAUSE SKIN IRRITATION. Remove contaminated clothing. Flush thoroughly with mild soap and water. If irritation occurs, get medical attention.

Section 5 Fire Fighting Measures

Suitable Extinguishing Media: Use any media suitable for extinguishing supporting fire.

Protective Actions for Fire-fighters: In fire conditions, wear a NIOSH/MSHA-approved self-contained breathing apparatus and full protective gear. Use water spray to keep fire-exposed containers cool.

Specific Hazards: In fire conditions, water may evaporate from this solution which may cause hazardous decomposition products to be formed as dust or fume.

Section 6 Accidental Release Measures

Personal Precautions: Evacuate personnel to safe area. Use proper personal protective equipment as indicated in Section 8. Provide adequate ventilation.

Environmental Precautions: Avoid runoff into storm sewers and ditches which lead to waterways.

Containment and Cleanup: Absorb with inert dry material, sweep or vacuum up and place in a suitable container for proper disposal. Wash spill area with soap and water.

Section 7 Handling and Storage

Precautions for Safe Handling: Read label on container before using. Do not wear contact lenses when working with chemicals. Keep out of reach of children. Avoid contact with eyes, skin and clothing. Do not inhale vapors, spray or mist. Use with adequate ventilation. Avoid ingestion. Wash thoroughly after handling. Remove and wash clothing before reuse.

Conditions for Safe Storage: Store in a cool, well-ventilated area away from incompatible substances. Protect from light.

Section 8 Exposure controls / Personal Protection

Exposure Limits:	Chemical Name	ACGIH (TLV)	OSHA (PEL)	NIOSH (REL)
	Bromothymol blue	None established	None established	None established

Engineering controls: Facilities storing or utilizing this material should be equipped with an eyewash facility and a safety shower and fire extinguishing material. Personnel should wear safety glasses, goggles, or faceshield, lab coat or apron, appropriate protective gloves. Use adequate ventilation to keep airborne concentrations low.

Respiratory protection: None should be needed in normal laboratory handling at room temperatures. If misty conditions prevail, work in fume hood or wear a NIOSH/MSHA approved respirator.

Section 9 Physical and Chemical Properties

Appearance: Liquid, clear, blue-green. Odor: No odor. Odor threshold: Data not available. pH: Data not available Melting/Freezing point: Approx. 0°C (32°F) (water) Boiling point: Approx. 100°C (212°F) (water) Flash point: Data not available	Evaporation rate (Water = 1): <1 Flammability (solid/gas): Data not available. Explosion limits: Lower/Upper: Data not available Vapor pressure (mm Hg): 14 (water) Vapor density (Air = 1): 0.7 (water) Relative density (Specific gravity): Approx. 1.0 (water) Solubility(ies): Complete in water.	Partition coefficient: Data not available Auto-ignition temp.: Data not available Decomposition temp.: Data not available Viscosity: Data not available. Molecular formula: Mixture Molecular weight: Mixture
---	--	--

Section 10 Stability and Reactivity

Chemical stability: Stable

Hazardous polymerization: Will not occur.

Conditions to avoid: Excessive temperatures which cause evaporation. Protect from light.

Incompatibilities with other materials: Strong oxidizers.

Hazardous decomposition products: Carbon oxides, sulfur oxides and bromine gas.

Section 11 Toxicological Information

Acute toxicity: Data not available

Serious eye damage/irritation: Data not available

Germ cell mutagenicity: Data not available

Skin corrosion/irritation: Data not available

Respiratory or skin sensitization: Data not available

Carcinogenicity: Data not available

NTP: No component of this product present at levels greater than or equal to 0.1% is identified as a known or anticipated carcinogen by NTP.

IARC: No component of this product present at levels greater than or equal to 0.1% is identified as probable, possible or confirmed human carcinogen by IARC.

OSHA: No component of this product present at levels greater than or equal to 0.1% is identified as a carcinogen or potential carcinogen by OSHA.

Reproductive toxicity: Data not available

STOT-single exposure: Data not available

Aspiration hazard: Data not available

STOT-repeated exposure: Data not available

Potential health effects:

Inhalation: May be harmful if inhaled.

Ingestion: May be harmful if swallowed.

Skin: May cause irritation.

Eyes: May cause irritation.

Signs and symptoms of exposure: To the best of our knowledge the chemical, physical and toxicological properties have not been thoroughly investigated. Specific data is not available. Exercise appropriate procedures to minimize potential hazards.

Additional information: RTECS #: Data not available

Section 12 Ecological Information

Toxicity to fish: No data available

Toxicity to daphnia and other aquatic invertebrates: No data available

Toxicity to algae: No data available

Persistence and degradability: No data available

Bioaccumulative potential: No data available

Mobility in soil: No data available

PBT and vPvB assessment: No data available

Other adverse effects: An environmental hazard cannot be excluded in the event of unprofessional handling or disposal.

Section 13 Disposal Considerations

These disposal guidelines are intended for the disposal of catalog-size quantities only. Federal regulations may apply to empty container. State and/or local regulations may be different. Dispose of in accordance with all local, state and federal regulations or contract with a licensed chemical disposal agency.

Section 14 Transport Information

UN/NA number: Not applicable

Shipping name: Not Regulated

Hazard class: Not applicable

Packing group: Not applicable

Reportable Quantity: No

Marine pollutant: No

Exceptions: Not applicable

2012 ERG Guide # Not applicable

Section 15 Regulatory Information

A chemical is considered to be listed if the CAS number for the anhydrous form is on the Inventory list.

Component	TSCA	CERLCA (RQ)	RCRA code	DSL	NDSL	WHMIS Classification
Bromothymol blue, sodium salt	Listed	Not Listed	Not Listed	Listed	Not Listed	Not Listed

Section 16 Additional Information

The information contained herein is furnished without warranty of any kind. Employers should use this information only as a supplement to other information gathered by them and must make independent determinations of suitability and completeness of information from all sources to assure proper use of these materials and the safety and health of employees.

NTP: National Toxicology Program, IARC: International Agency for Research on Cancer, OSHA: Occupational Safety and Health Administration, STOT: Specific Target Organ Toxicity, SE: Single Exposure, RE: Repeated Exposure, ERG: Emergency Response Guidebook.

Section 1 Chemical Product and Company Information

Science Take-Out
80 Office Park Way
Pittsford, NY 14534
(585)764-5400

**CHEMTREC 24 Hour Emergency
Phone Number (800) 424-9300**
For laboratory use only. Not for drug, food or household use

Product	Water
Synonyms	"Arsenic Test Solution" (simulated)

Section 2 Hazards Identification

This substance or mixture has not been classified at this time according to the Globally Harmonized System (GHS) of Classification and Labeling of Chemicals.

Signal word: Not Classified
Pictograms: None required
Target organs: None known

GHS Classification: Not classified

GHS Label information: Hazard statement(s): Not classified

Ca Prop 65 - This product does not contain any chemicals known to the State of California to cause cancer, birth defects, or any other reproductive harm.

Section 3 Composition / Information on Ingredients

Chemical Name	CAS #	%	EINECS
Water	7732-18-15	100%	231-791-2

Section 4 First Aid Measures

INGESTION:

INHALATION:

EYE CONTACT:

SKIN ABSORPTION:

Section 5 Fire Fighting Measures

Suitable Extinguishing Media: Use any media suitable for extinguishing supporting fire.

Protective Actions for Fire-fighters: In fire conditions, wear a NIOSH/MSHA-approved self-contained breathing apparatus and full protective gear. Use water spray to keep fire-exposed containers cool.

Specific Hazards: Data not available.

Section 6 Accidental Release Measures

Personal Precautions: Evacuate personnel to safe area. Use proper personal protective equipment as indicated in Section 8. Provide adequate ventilation.

Environmental Precautions: Avoid runoff into storm sewers and ditches which lead to waterways.

Containment and Cleanup: Absorb with inert dry material, sweep or vacuum up and place in a suitable container for proper disposal.

Section 7 Handling and Storage

Precautions for Safe Handling: Data not available.

Conditions for Safe Storage: Data not available.

Section 1 Chemical Product and Company Information

Science Take-Out
80 Office Park Way
Pittsford, NY 14534
(585)764-5400

**CHEMTREC 24 Hour Emergency
Phone Number (800) 424-9300**
For laboratory use only. Not for drug, food or household use

Product	Methyl red, 0.05% aqueous solution
Synonyms	"Coliform Bacteria Test Solution" (simulated)

Section 2 Hazards Identification

This substance or mixture has not been classified at this time according to the Globally Harmonized System (GHS) of Classification and Labeling of Chemicals.

Signal word: Not classified
Pictograms: Not classified
Target organs: None known

GHS Classification: Not classified
GHS Label information: Not classified
Precautionary Statement: Not classified

Supplementary information:

Do not breathe vapors, spray or mist. Do not get in eyes, on skin, or on clothing. Wear protective gloves/protective clothing/eye protection/face protection. Wash hands thoroughly after handling. Get medical attention if you feel unwell.

Ca Prop 65 - This product does not contain any chemicals known to the State of California to cause cancer, birth defects, or any other reproductive harm.

Section 3 Composition / Information on Ingredients

Chemical Name	CAS #	%	EINECS
Water	7732-18-5	99.95%	231-791-2
Methyl red, sodium salt	845-10-3	0.05%	212-682-9

Section 4 First Aid Measures

INGESTION: MAY BE HARMFUL IF SWALLOWED. Call physician or Poison Control Center immediately. Induce vomiting only if advised by appropriate medical personnel. Never give anything by mouth to an unconscious person.

INHALATION: MAY BE HARMFUL IF INHALED. MAY CAUSE RESPIRATORY TRACT IRRITATION. Remove to fresh air. If not breathing, give artificial respiration. If breathing is difficult, give oxygen. Get medical attention.

EYE CONTACT: MAY CAUSE EYE IRRITATION. Check for and remove contact lenses. Flush thoroughly with water for at least 15 minutes, lifting upper and lower eyelids occasionally. Get immediate medical attention.

SKIN ABSORPTION: MAY CAUSE SKIN IRRITATION. Remove contaminated clothing. Flush thoroughly with mild soap and water. If irritation occurs, get medical attention.

Section 5 Fire Fighting Measures

Suitable Extinguishing Media: Carbon dioxide, dry chemical, dry sand, alcohol foam.

Protective Actions for Fire-fighters: In fire conditions, wear a NIOSH/MSHA-approved self-contained breathing apparatus and full protective gear. Use water spray to keep fire-exposed containers cool.

Specific Hazards: In fire conditions, water may evaporate from this solution which may cause hazardous decomposition products to be formed as dust or fume.

Section 6 Accidental Release Measures

Personal Precautions: Evacuate personnel to safe area. Use proper personal protective equipment as indicated in Section 8. Provide adequate ventilation.

Environmental Precautions: Avoid runoff into storm sewers and ditches which lead to waterways.

Containment and Cleanup: Absorb with inert dry material, sweep or vacuum up and place in a suitable container for proper disposal. Wash spill area with soap and water.

Section 7 Handling and Storage

Precautions for Safe Handling: Read label on container before using. Do not wear contact lenses when working with chemicals. Keep out of reach of children. Avoid contact with eyes, skin and clothing. Do not inhale vapors, spray or mist. Use with adequate ventilation. Avoid ingestion. Wash thoroughly after handling. Remove and wash clothing before reuse.

Conditions for Safe Storage: Store in a cool, well-ventilated area away from incompatible substances. Protect from light.

Section 8 Exposure controls / Personal Protection

Exposure Limits:	Chemical Name	ACGIH (TLV)	OSHA (PEL)	NIOSH (REL)
	Methyl red	None established	None established	None established

Engineering controls: Facilities storing or utilizing this material should be equipped with an eyewash facility and a safety shower and fire extinguishing material. Personnel should wear safety glasses, goggles, or faceshield, lab coat or apron, appropriate protective gloves. Use adequate ventilation to keep airborne concentrations low.

Respiratory protection: None should be needed in normal laboratory handling at room temperatures. If misty conditions prevail, work in fume hood or wear a NIOSH/MSHA approved respirator.

Section 9 Physical and Chemical Properties

Appearance: Liquid, clear, red. Odor: No odor. Odor threshold: Data not available. pH: Data not available Melting/Freezing point: Approx. 0°C (32°F) (water) Boiling point: Approx. 100°C (212°F) (water) Flash point: Data not available	Evaporation rate (Water = 1): <1 Flammability (solid/gas): Data not available. Explosion limits: Lower/Upper: Data not available Vapor pressure (mm Hg): 14 (water) Vapor density (Air = 1): 0.7 (water) Relative density (Specific gravity): Approx. 1.0 (water) Solubility(ies): Complete in water.	Partition coefficient: Data not available Auto-ignition temp.: Data not available Decomposition temp.: Data not available Viscosity: Data not available. Molecular formula: Mixture Molecular weight: Mixture
--	--	--

Section 10 Stability and Reactivity

Chemical stability: Stable

Hazardous polymerization: Will not occur.

Conditions to avoid: Excessive temperatures which cause evaporation. Protect from light.

Incompatibilities with other materials: Strong oxidizers, reducing agents.

Hazardous decomposition products: Carbon oxides, nitrogen oxides and sodium oxides.

Section 11 Toxicological Information

Acute toxicity: Oral-rat TDL: 12000 mg/kg [Methyl red]

Serious eye damage/irritation: Data not available

Germ cell mutagenicity: Data not available

Skin corrosion/irritation: Data not available

Respiratory or skin sensitization: Data not available

Carcinogenicity: Data not available

NTP: No component of this product present at levels greater than or equal to 0.1% is identified as a known or anticipated carcinogen by NTP.

IARC: No component of this product present at levels greater than or equal to 0.1% is identified as probable, possible or confirmed human carcinogen by IARC.

OSHA: No component of this product present at levels greater than or equal to 0.1% is identified as a carcinogen or potential carcinogen by OSHA.

Reproductive toxicity: Data not available

STOT-single exposure: Data not available

Aspiration hazard: Data not available

STOT-repeated exposure: Data not available

Potential health effects:

Inhalation: May be harmful if inhaled.

Ingestion: May be harmful if swallowed.

Skin: May cause irritation.

Eyes: May cause irritation.

Signs and symptoms of exposure: To the best of our knowledge the chemical, physical and toxicological properties have not been thoroughly investigated. Specific data is not available. Exercise appropriate procedures to minimize potential hazards.

Additional information: RTECS #: DG8960000 [Methyl red]

Section 12 Ecological Information

Toxicity to fish: No data available

Toxicity to daphnia and other aquatic invertebrates: No data available

Toxicity to algae: No data available

Persistence and degradability: No data available

Bioaccumulative potential: No data available

Mobility in soil: No data available

PBT and vPvB assessment: No data available

Other adverse effects: An environmental hazard cannot be excluded in the event of unprofessional handling or disposal.

Section 13 Disposal Considerations

These disposal guidelines are intended for the disposal of catalog-size quantities only. Federal regulations may apply to empty container. State and/or local regulations may be different. Dispose of in accordance with all local, state and federal regulations or contract with a licensed chemical disposal agency.

Section 14 Transport Information

UN/NA number: Not applicable

Shipping name: Not Regulated

Hazard class: Not applicable

Packing group: Not applicable

Reportable Quantity: No

Marine pollutant: No

Exceptions: Not applicable

2012 ERG Guide # Not applicable

Section 15 Regulatory Information

A chemical is considered to be listed if the CAS number for the anhydrous form is on the Inventory list.

Component	TSCA	CERLCA (RQ)	RCRA code	DSL	NDSL	WHMIS Classification
Methyl red, sodium salt	Listed	Not Listed	Not Listed	Listed	Not Listed	Not Listed

Section 16 Additional Information

The information contained herein is furnished without warranty of any kind. Employers should use this information only as a supplement to other information gathered by them and must make independent determinations of suitability and completeness of information from all sources to assure proper use of these materials and the safety and health of employees.

NTP: National Toxicology Program, IARC: International Agency for Research on Cancer, OSHA: Occupational Safety and Health Administration, STOT: Specific Target Organ Toxicity, SE: Single Exposure, RE: Repeated Exposure, ERG: Emergency Response Guidebook.

Section 1 L'information de produit chimique et de compagnie

Science Take-Out
80 Office Park Way
Pittsford, NY 14534
(585)764-5400

**CHEMTREC 24 Numéros De Téléphone
De Secours D'Heure (800) 424-9300**
Pour l'usage de laboratoire seulement.
Pas pour l'usage de drogue, de nourriture ou de ménage.

Produit	Solution de tampon pH10
Synonymes	"Well Water Sample" (simulated)

Section 2 Identification De Risques

Cette substance ou un mélange n'a pas été classé à ce moment selon le Système général harmonisé (SGH) de classification et d'étiquetage des produits chimiques.

Mention d'avertissement: AVERTISSEMENT

Pictogrammes: Aucune requise

Les organes cibles: Aucun connu.

Classification par le GHS:

Skin irritation (Category 3)

Eye irritation (Category 2B)

Renseignements sur l'étiquette GHS: Mention de danger(s):

H316: Provoque une légère irritation cutanée.

H320: Provoque une irritation des yeux.

Déclarations de précaution(s):

P264: Se laver les mains soigneusement après manipulation.

P305+P351+P338: EN CAS DE CONTACT AVEC LES YEUX: Rincer avec précaution à l'eau pendant plusieurs minutes. Enlever les lentilles de contact si la victime en porte et si elles peuvent être facilement enlevées. Continuer à rincer.

P332+P313: En cas d'irritation cutanée: Obtenir des soins médicaux.

P337+P313: Si l'irritation oculaire persiste: Obtenir des soins médicaux.

CA Prop 65 - Ce produit ne contient pas de produits chimiques connus à l'État de Californie pour causer le cancer, des malformations congénitales, ou toute autre atteinte à la reproduction.

Section 3 Composition / Information Sur Des Ingrédients

Nommé Chimique	CAS #	%	EINECS
L'eau	7732-18-5	99.08%	231-791-2
Chlorure de potassium	7447-40-7	0.40%	231-211-8
Acide borique	10043-35-3	0.33%	233-139-2
Hydroxide de sodium	1310-73-2	0.19%	215-185-5

Section 4 Mesures De Premiers Soins

INGESTION: Appeler un médecin ou un centre antipoison immédiatement. Provoquer le vomissement seulement si elle est informée par le personnel compétent médicaux. Ne jamais rien donner par la bouche à une personne inconsciente.

INHALATION: Sortir au grand air. Si elle ne respire pas, pratiquer la respiration artificielle. Si la respiration est difficile, donner de l'oxygène. Obtenir des soins médicaux.

CONTACT AVEC LES YEUX: Vérifier et enlever les lentilles de contact. Rincer abondamment à l'eau pendant au moins 15 minutes, en soulevant les paupières inférieures et supérieures de temps en temps. Obtenez une attention médicale immédiate.

ABSORPTION PAR LA PEAU: Enlever les vêtements contaminés. Rincer soigneusement avec du savon doux et d'eau. En cas d'irritation, consulter un médecin.

Section 5 Mesures De Lutte Contre l'Incendie

Moyens d'extinction: Utilisez des supports adaptés pour éteindre le feu à l'appui.

Actions de protection pour les sapeurs-pompiers: En cas d'incendie, porter un appareil respiratoire NIOSH / MSHA approuvé autonome et un équipement complet de protection. Utiliser un jet d'eau pour maintenir incendie refroidir les conteneurs exposés.

Dangers spécifiques: En cas d'incendie, des gaz irritants et très toxiques peuvent être générés par la décomposition thermique ou la combustion.

Section 6 Mesures De Déchargement Accidentel

Précautions personnelles: Évacuer le personnel vers la zone sûre. Utiliser un équipement de protection personnelle comme indiqué dans la Section 8. Assurer une ventilation adéquate.

Précautions environnementales: Éviter tout ruissellement vers les égouts pluviaux et les fossés qui aboutissent aux voies navigables.

Confinement et de nettoyage: Absorbent avec le matériel sec inerte, balayez ou nettoyez à l'aspirateur vers le haut et placez dans un récipient approprié pour la disposition appropriée. Laver la zone de déversement avec du savon et de l'eau.

Section 7 Manipulation Et Stockage

Précautions pour la manutention en toute sécurité: Lire l'étiquette sur le contenant avant d'utiliser. Ne pas porter de lentilles cornéennes lorsque vous travaillez avec des produits chimiques. Tenir hors de portée des enfants. Éviter tout contact avec les yeux, la peau et les vêtements. Ne pas inhaler les vapeurs, les embruns ou le brouillard. Utiliser avec une ventilation adéquate. Éviter l'ingestion. Bien se laver après la manipulation. Retirer et laver les vêtements avant de les réutiliser.

Conditions de stockage: Stocker dans un endroit frais et bien aéré, loin des substances incompatibles.

Section 8 Commandes D'Exposition / Protection Personnelle

Limites d'exposition:	Nommé Chimique	ACGIH (TLV)	OSHA (PEL)	NIOSH (REL)
	Chlorure de potassium	Aucun établi.	Aucun établi.	Aucun établi.

Contrôles d'ingénierie: Les installations d'entreposage ou d'utilisation de ce matériel doit être équipé d'une douche oculaire et une douche de sécurité et le matériel d'extinction d'incendie. Le personnel doit porter des lunettes de sécurité, des lunettes, ou un écran facial, une blouse de laboratoire ou tablier, des gants protecteurs appropriés. Utiliser une ventilation adéquate pour maintenir les concentrations atmosphériques faible.

Protection respiratoire: Aucun ne devrait être nécessaire dans le laboratoire normal manipulant aux températures ambiantes. Si les conditions brumeuses prévaloir, travailler dans la hotte ou de porter un masque respiratoire approuvé NIOSH / MSHA.

Section 9 Propriétés Physiques Et Chimiques

Apparence: Clair, liquide incolore. Odeur: Aucun odeur. Seuil de l'odeur: Données non disponibles. pH: 10.0 Point de fusion / congélation: Environ 0°C (32°F) (eau) Point d'ébullition: Environ 100°C (212°F) (eau) Point d'éclair: Données non disponibles	Taux d'évaporation (Eau = 1): <1 Inflammabilité (solide / gaz): Données non disponibles. Limites d'explosivité: Bas / Max: Données non disponibles Pression de vapeur (mm Hg): 14 (eau) Densité de vapeur (Air = 1): 0.7 (eau) Densité relative (gravité spécifique): Environ 1.0 (eau) Solubilité (s): Complet dans l'eau.	Coefficient de partage: Données non disponibles Auto-inflammation: Données non disponibles Température de décomposition: Données non disponibles. Viscosité: Données non disponibles. Formule moléculaire: Mélange Poids moléculaire: Mélange
--	--	--

Section 10 Stabilité Et Réactivité

Stabilité chimique: Stable **Polymérisation dangereuse:** N'aura pas lieu.
Conditions à éviter: Les températures excessives qui causent l'évaporation.
Incompatibilités avec d'autres matériaux: Acides, alcalis, et l'air changera la capacité de l'amortisseur.
Produits de décomposition dangereux: Oxyde de bore et gaz de chlore.

Section 11 L'Information Toxicologique

Toxicité aiguë: Données non disponibles
La corrosion de la peau et l'irritation: Données non disponibles
Des lésions oculaires graves / irritation: Données non disponibles
Respiratoire ou sensibilisation de la peau: Données non disponibles
Mutagénicité des cellules germinales: Données non disponibles
Cancérogène: Données non disponibles
NTP: Aucun composant de ce produit présent à des niveaux supérieurs ou égaux à 0,1% n'a été identifié comme cancérigène reconnu ou présumé par NTP.
IARC: Aucun composant de ce produit présent à des niveaux supérieurs ou égaux à 0,1% n'a été identifié comme cancérigène probable, possible ou confirmé par IARC.
OSHA: Aucun composant de ce produit présent à des niveaux supérieurs ou égaux à 0,1% n'a été identifié comme cancérigène ni comme cancérigène possible par OSHA.
Reproductive toxicity: Données non disponibles
STOT-exposition unique: Données non disponibles.
STOT-une exposition répétée: Données non disponibles
Risque d'aspiration: Données non disponibles
Effets d'une surexposition:
Inhalation: Peut être nocif en cas d'inhalation.
Ingestion: Peut être nocif en cas d'ingestion.
Peau: Peut causer une légère irritation.
Yeux: Peut causer une légère irritation.
Les signes et les symptômes de l'exposition: Pour le meilleur de notre connaissance les propriétés chimiques, physiques et toxicologiques n'ont pas été étudiées à fond.
Les données spécifiques n'est pas disponible. Exercice des procédures appropriées afin de minimiser les dangers potentiels.
Informations complémentaires: RTECS #: Données non disponibles

Section 12 L'Information Écologique

Toxicité pour les poissons: Pas de données disponible
Toxicité pour les daphnies et autres invertébrés aquatiques: Pas de données disponible
Toxicité pour les algues: Pas de données disponible
Persistance et dégradabilité: Pas de données disponible **Potentiel de bioaccumulation:** Pas de données disponible
Mobilité dans le sol: Pas de données disponibles **Évaluation PBT et vPvB:** Pas de données disponibles
Autres effets indésirables: Un danger pour l'environnement ne peut pas être exclu dans l'éventualité d'une manipulation ou d'élimination.

Section 13 Considérations De Disposition

Ces lignes directrices sont destinées à l'élimination de la disposition d'un catalogue de taille seules les quantités. Les règlements fédéraux peuvent s'appliquer aux contenants vides. Des réglementations nationales et / ou local peut être différent. Éliminer conformément à toutes les réglementations locales, provinciales et fédérales ou d'un contrat avec une agence élimination des produits chimiques sous licence..

Section 14 L'Information De Transport

Numéro UN / NA: Non applicable **Nom d'expédition:** Non réglé
Classe de danger: Non applicable **Groupe d'emballage:** Non applicable **Quantité à déclarer:** Non **Polluant marin:** Non
Exceptions: Non applicable **2012 ERG Guide #:** Non applicable

Section 15 L'Information De Normalisation

Un produit chimique est considéré comme inscrit si le numéro CAS pour la forme anhydre est sur la liste d'inventaire.

Composant	TSCA	CERLCA (RQ)	RCRA code	DSL	NDSL	Classification SIMDUT
Phosphate de potassium	Listed	Not Listed	Not Listed	Listed	Not Listed	Produit non contrôlé

Section 16 L'Information Additionnelle

Les informations contenues dans ce document sont fournis sans garantie d'aucune sorte. Les employeurs devraient considérer cette information seulement comme complément à d'autres informations recueillies par eux et doivent prendre des décisions indépendantes de la pertinence et l'exhaustivité de l'information de toutes les sources afin d'assurer une utilisation correcte de ces matériaux et de la sécurité et la santé des employés. NTP: National Toxicology Program, IARC: International Agency for Research on Cancer, OSHA: Occupational Safety and Health Administration, STOT: Specific Target Organ Toxicity, SE: Single Exposure, RE: Repeated Exposure.

Section 1 L'information de produit chimique et de compagnie

Science Take-Out
80 Office Park Way
Pittsford, NY 14534
(585)764-5400

**CHEMTREC 24 Numéros De Téléphone
De Secours D'Heure (800) 424-9300**
Pour l'usage de laboratoire seulement.
Pas pour l'usage de drogue, de nourriture ou de ménage.

Produit	Solution de tampon pH3
Synonymes	"Bottled Water Sample" (simulated), "Tap Water Sample" (simulated)

Section 2 Identification De Risques

Cette substance ou un mélange n'a pas été classé à ce moment selon le Système général harmonisé (SGH) de classification et d'étiquetage des produits chimiques.

Mention d'avertissement: AVERTISSEMENT

Pictogrammes: Aucune requise

Les organes cibles: Aucun connu.

Classification par le GHS:

Skin irritation (Category 3)

Eye irritation (Category 2B)

Renseignements sur l'étiquette GHS: Mention de danger(s):

H316: Provoque une légère irritation cutanée.

H320: Provoque une irritation des yeux.

Déclarations de précaution(s):

P264: Se laver les mains soigneusement après manipulation.

P305+P351+P338: EN CAS DE CONTACT AVEC LES YEUX: Rincer avec précaution à l'eau pendant plusieurs minutes. Enlever les lentilles de contact si la victime en porte et si elles peuvent être facilement enlevées. Continuer à rincer.

P332+P313: En cas d'irritation cutanée: Obtenir des soins médicaux.

P337+P313: Si l'irritation oculaire persiste: Obtenir des soins médicaux.

CA Prop 65 - Ce produit ne contient pas de produits chimiques connus à l'État de Californie pour causer le cancer, des malformations congénitales, ou toute autre atteinte à la reproduction.

Section 3 Composition / Information Sur Des Ingrédients

Nommé Chimique	CAS #	%	EINECS
L'eau	7732-18-5	99.74%	231-791-2
Biphthalate de potassium	877-24-7	1.08%	212-889-4
Acide chlorhydrique	7647-01-0	0.18%	231-595-7

Section 4 Mesures De Premiers Soins

INGESTION: Appeler un médecin ou un centre antipoison immédiatement. Provoquer le vomissement seulement si elle est informée par le personnel compétent médicaux. Ne jamais rien donner par la bouche à une personne inconsciente.

INHALATION: Sortir au grand air. Si elle ne respire pas, pratiquer la respiration artificielle. Si la respiration est difficile, donner de l'oxygène. Obtenir des soins médicaux.

CONTACT AVEC LES YEUX: Vérifier et enlever les lentilles de contact. Rincer abondamment à l'eau pendant au moins 15 minutes, en soulevant les paupières inférieures et supérieures de temps en temps. Obtenez une attention médicale immédiate.

ABSORPTION PAR LA PEAU: Enlever les vêtements contaminés. Rincer soigneusement avec du savon doux et d'eau. En cas d'irritation, consulter un médecin.

Section 5 Mesures De Lutte Contre l'Incendie

Moyens d'extinction: Utilisez des supports adaptés pour éteindre le feu à l'appui.

Actions de protection pour les sapeurs-pompiers: En cas d'incendie, porter un appareil respiratoire NIOSH / MSHA approuvé autonome et un équipement complet de protection. Utiliser un jet d'eau pour maintenir incendie refroidir les conteneurs exposés.

Dangers spécifiques: En cas d'incendie, des gaz irritants et très toxiques peuvent être générés par la décomposition thermique ou la combustion.

Section 6 Mesures De Déchargement Accidentel

Précautions personnelles: Évacuer le personnel vers la zone sûre. Utiliser un équipement de protection personnelle comme indiqué dans la Section 8. Assurer une ventilation adéquate.

Précautions environnementales: Éviter tout ruissellement vers les égouts pluviaux et les fossés qui aboutissent aux voies navigables.

Confinement et de nettoyage: Absorbent avec le matériel sec inerte, balayez ou nettoyez à l'aspirateur vers le haut et placez dans un récipient approprié pour la disposition appropriée. Laver la zone de déversement avec du savon et de l'eau.

Section 7 Manipulation Et Stockage

Précautions pour la manutention en toute sécurité: Lire l'étiquette sur le contenant avant d'utiliser. Ne pas porter de lentilles cornéennes lorsque vous travaillez avec des produits chimiques. Tenir hors de portée des enfants. Éviter tout contact avec les yeux, la peau et les vêtements. Ne pas inhaler les vapeurs, les embruns ou le brouillard. Utiliser avec une ventilation adéquate. Éviter l'ingestion. Bien se laver après la manipulation. Retirer et laver les vêtements avant de les réutiliser.

Conditions de stockage: Stocker dans un endroit frais et bien aéré, loin des substances incompatibles.

Section 8 Commandes D'Exposition / Protection Personnelle

Limites d'exposition:	Nommé Chimique	ACGIH (TLV)	OSHA (PEL)	NIOSH (REL)
	Aucun	Aucun établi.	Aucun établi.	Aucun établi.

Contrôles d'ingénierie: Les installations d'entreposage ou d'utilisation de ce matériel doit être équipé d'une douche oculaire et une douche de sécurité et le matériel d'extinction d'incendie. Le personnel doit porter des lunettes de sécurité, des lunettes, ou un écran facial, une blouse de laboratoire ou tablier, des gants protecteurs appropriés. Utiliser une ventilation adéquate pour maintenir les concentrations atmosphériques faible.

Protection respiratoire: Aucun ne devrait être nécessaire dans le laboratoire normal manipulant aux températures ambiantes. Si les conditions brumeuses prévaloir, travailler dans la hotte ou de porter un masque respiratoire approuvé NIOSH / MSHA.

Section 9 Propriétés Physiques Et Chimiques

Apparence: Clair, liquide incolore. Odeur: Aucun odeur. Seuil de l'odeur: Données non disponibles. pH: 3.0 Point de fusion / congélation: Environ 0°C (32°F) (eau) Point d'ébullition: Environ 100°C (212°F) (eau) Point d'éclair: Données non disponibles	Taux d'évaporation (Eau = 1): <1 Inflammabilité (solide / gaz): Données non disponibles. Limites d'explosivité: Bas / Max: Données non disponibles Pression de vapeur (mm Hg): 14 (eau) Densité de vapeur (Air = 1): 0.7 (eau) Densité relative (gravité spécifique): Environ 1.0 (eau) Solubilité (s): Complet dans l'eau.	Coefficient de partage: Données non disponibles Auto-inflammation: Données non disponibles Température de décomposition: Données non disponibles. Viscosité: Données non disponibles. Formule moléculaire: Mélange Poids moléculaire: Mélange
---	--	--

Section 10 Stabilité Et Réactivité

Stabilité chimique: Stable **Polymérisation dangereuse:** N'aura pas lieu.
Conditions à éviter: Les températures excessives qui causent l'évaporation.
Incompatibilités avec d'autres matériaux: Acides, alcalis, et l'air changera la capacité de l'amortisseur.
Produits de décomposition dangereux: Oxyde de carbone et gaz de chlore.

Section 11 L'Information Toxicologique

Toxicité aiguë: Données non disponibles
La corrosion de la peau et l'irritation: Données non disponibles
Des lésions oculaires graves / irritation: Données non disponibles
Respiratoire ou sensibilisation de la peau: Données non disponibles
Mutagénicité des cellules germinales: Données non disponibles
Cancérogène: Données non disponibles
NTP: Aucun composant de ce produit présent à des niveaux supérieurs ou égaux à 0,1% n'a été identifié comme cancérigène reconnu ou présumé par NTP.
IARC: Aucun composant de ce produit présent à des niveaux supérieurs ou égaux à 0,1% n'a été identifié comme cancérigène probable, possible ou confirmé par IARC.
OSHA: Aucun composant de ce produit présent à des niveaux supérieurs ou égaux à 0,1% n'a été identifié comme cancérigène ni comme cancérigène possible par OSHA.
Reproductive toxicity: Données non disponibles
STOT-exposition unique: Données non disponibles.
STOT-une exposition répétée: Données non disponibles
Risque d'aspiration: Données non disponibles
Effets d'une surexposition:
Inhalation: Peut être nocif en cas d'inhalation.
Ingestion: Peut être nocif en cas d'ingestion.
Peau: Peut causer une légère irritation.
Yeux: Peut causer une légère irritation.
Les signes et les symptômes de l'exposition: Pour le meilleur de notre connaissance les propriétés chimiques, physiques et toxicologiques n'ont pas été étudiées à fond.
Les données spécifiques n'est pas disponible. Exercice des procédures appropriées afin de minimiser les dangers potentiels.
Informations complémentaires: RTECS #: Données non disponibles

Section 12 L'Information Écologique

Toxicité pour les poissons: Pas de données disponible
Toxicité pour les daphnies et autres invertébrés aquatiques: Pas de données disponible
Toxicité pour les algues: Pas de données disponible
Persistance et dégradabilité: Pas de données disponible **Potentiel de bioaccumulation:** Pas de données disponible
Mobilité dans le sol: Pas de données disponibles **Évaluation PBT et vPvB:** Pas de données disponibles
Autres effets indésirables: Un danger pour l'environnement ne peut pas être exclu dans l'éventualité d'une manipulation ou d'élimination.

Section 13 Considérations De Disposition

Ces lignes directrices sont destinées à l'élimination de la disposition d'un catalogue de taille seules les quantités. Les règlements fédéraux peuvent s'appliquer aux contenants vides. Des réglementations nationales et / ou local peut être différent. Éliminer conformément à toutes les réglementations locales, provinciales et fédérales ou d'un contrat avec une agence élimination des produits chimiques sous licence..

Section 14 L'Information De Transport

Numéro UN / NA: Non applicable **Nom d'expédition:** Non réglé
Classe de danger: Non applicable **Groupe d'emballage:** Non applicable **Quantité à déclarer:** Non **Polluant marin:** Non
Exceptions: Non applicable **2012 ERG Guide #:** Non applicable

Section 15 L'Information De Normalisation

Un produit chimique est considéré comme inscrit si le numéro CAS pour la forme anhydre est sur la liste d'inventaire.

Composant	TSCA	CERLCA (RQ)	RCRA code	DSL	NDSL	Classification SIMDUT
Biphtalate de potassium	Listed	Not Listed	Not Listed	Listed	Not Listed	Produit non contrôlé
Acide chlorhydrique	Listed	5,000 lbs (2270 kg)	Not Listed	Listed	Not Listed	E

Section 16 L'Information Additionnelle

Les informations contenues dans ce document sont fournis sans garantie d'aucune sorte. Les employeurs devraient considérer cette information seulement comme complément à d'autres informations recueillies par eux et doivent prendre des décisions indépendantes de la pertinence et l'exhaustivité de l'information de toutes les sources afin d'assurer une utilisation correcte de ces matériaux et de la sécurité et la santé des employés. NTP: National Toxicology Program, IARC: International Agency for Research on Cancer, OSHA: Occupational Safety and Health Administration, STOT: Specific Target Organ Toxicity, SE: Single Exposure, RE: Repeated Exposure.

Section 1 L'information de produit chimique et de compagnie

Science Take-Out
80 Office Park Way
Pittsford, NY 14534
(585)764-5400

**CHEMTREC 24 Numéros De Téléphone
De Secours D'Heure (800) 424-9300**
Pour l'usage de laboratoire seulement.
Pas pour l'usage de drogue, de nourriture ou de ménage.

Produit	Blue de bromothymol, 0.1% solute
Synonymes	"Nitrate Test Solution" (simulated)

Section 2 Identification De Risques

Cette substance ou un mélange n'a pas été classé à ce moment selon le Système général harmonisé (SGH) de classification et d'étiquetage des produits chimiques.

Mention d'avertissement: Non classé

Pictogrammes: Non classé

Les organes cibles: Aucun connu.

Classification par le GHS: Non classé

Renseignements sur l'étiquette GHS: Mention de danger(s): Non classé

Déclarations de précaution(s): Non classé

Des renseignements supplémentaires:

Ne pas respirer les vapeurs, les embruns ou le brouillard. Éviter tout contact avec les yeux, la peau ou les vêtements. Porter des gants de protection / des vêtements de protection / un équipement de protection des yeux / du visage. Se laver les mains soigneusement après manipulation. Consulter un médecin en cas de malaise.

CA Prop 65 - Ce produit ne contient pas de produits chimiques connus à l'État de Californie pour causer le cancer, des malformations congénitales, ou toute autre atteinte à la reproduction.

Section 3 Composition / Information Sur Des Ingrédients

Nommé Chimique	CAS #	%	EINECS
L'eau	7732-18-5	99.9%	231-791-2
Bleu de bromothymol, sel sodique	34722-90-2	0.1%	252-169-7

Section 4 Mesures De Premiers Soins

INGESTION: Appeler un médecin ou un centre antipoison immédiatement. Provoquer le vomissement seulement si elle est informée par le personnel compétent médicaux. Ne jamais rien donner par la bouche à une personne inconsciente.

INHALATION: Sortir au grand air. Si elle ne respire pas, pratiquer la respiration artificielle. Si la respiration est difficile, donner de l'oxygène. Obtenir des soins médicaux.

CONTACT AVEC LES YEUX: Vérifier et enlever les lentilles de contact. Rincer abondamment à l'eau pendant au moins 15 minutes, en soulevant les paupières inférieures et supérieures de temps en temps. Obtenez une attention médicale immédiate.

ABSORPTION PAR LA PEAU: Enlever les vêtements contaminés. Rincer soigneusement avec du savon doux et d'eau. En cas d'irritation, consulter un médecin.

Section 5 Mesures De Lutte Contre l'Incendie

Moyens d'extinction: Utilisez des supports adaptés pour éteindre le feu à l'appui.

Actions de protection pour les sapeurs-pompiers: En cas d'incendie, porter un appareil respiratoire NIOSH / MSHA approuvé autonome et un équipement complet de protection. Utiliser un jet d'eau pour maintenir incendie refroidir les conteneurs exposés.

Dangers spécifiques: En cas de feu, de l'eau peut s'évaporer à partir de cette solution, qui peut causer les produits dangereux de décomposition à être formée comme poussière ou vapeur.

Section 6 Mesures De Déchargement Accidentel

Précautions personnelles: Évacuer le personnel vers la zone sûre. Utiliser un équipement de protection personnelle comme indiqué dans la Section 8. Assurer une ventilation adéquate.

Précautions environnementales: Éviter tout ruissellement vers les égouts pluviaux et les fossés qui aboutissent aux voies navigables.

Confinement et de nettoyage: Absorbent avec le matériel sec inerte, balayez ou nettoyez à l'aspirateur vers le haut et placez dans un récipient approprié pour la disposition appropriée. Laver la zone de déversement avec du savon et de l'eau.

Section 7 Manipulation Et Stockage

Précautions pour la manutention en toute sécurité: Lire l'étiquette sur le contenant avant d'utiliser. Ne pas porter de lentilles cornéennes lorsque vous travaillez avec des produits chimiques. Tenir hors de portée des enfants. Éviter tout contact avec les yeux, la peau et les vêtements. Ne pas inhaler les vapeurs, les embruns ou le brouillard. Utiliser avec une ventilation adéquate. Éviter l'ingestion. Bien se laver après la manipulation. Retirer et laver les vêtements avant de les réutiliser.

Conditions de stockage: Stocker dans un endroit frais et bien aéré, loin des substances incompatibles.

Section 8 Commandes D'Exposition / Protection Personnelle

Limites d'exposition:	Nommé Chimique	ACGIH (TLV)	OSHA (PEL)	NIOSH (REL)
	Bleu de bromothymol	Aucun établi.	Aucun établi.	Aucun établi.

Contrôles d'ingénierie: Les installations d'entreposage ou d'utilisation de ce matériel doit être équipé d'une douche oculaire et une douche de sécurité et le matériel d'extinction d'incendie. Le personnel doit porter des lunettes de sécurité, des lunettes, ou un écran facial, une blouse de laboratoire ou tablier, des gants protecteurs appropriés. Utiliser une ventilation adéquate pour maintenir les concentrations atmosphériques faible.

Protection respiratoire: Aucun ne devrait être nécessaire dans le laboratoire normal manipulant aux températures ambiantes. Si les conditions brumeuses prévaloir, travailler dans la hotte ou de porter un masque respiratoire approuvé NIOSH / MSHA.

Section 9 Propriétés Physiques Et Chimiques

Apparence: Liquide, clair, bleu-vert. Odeur: Aucun odeur. Seuil de l'odeur: Données non disponibles. pH: Données non disponibles. Point de fusion / congélation: Environ 0°C (32°F) (eau) Point d'ébullition: Environ 100°C (212°F) (eau) Point d'éclair: Données non disponibles	Taux d'évaporation (Eau = 1): <1 Inflammabilité (solide / gaz): Données non disponibles. Limites d'explosivité: Bas / Max: Données non disponibles Pression de vapeur (mm Hg): 14 (eau) Densité de vapeur (Air = 1): 0.7 (eau) Densité relative (gravité spécifique): Environ 1.0 (eau) Solubilité (s): Complet dans l'eau.	Coefficient de partage: Données non disponibles Auto-inflammation: Données non disponibles Température de décomposition: Données non disponibles. Viscosité: Données non disponibles. Formule moléculaire: Mélange Poids moléculaire: Mélange
--	--	--

Section 10 Stabilité Et Réactivité

Stabilité chimique: Stable **Polymérisation dangereuse:** N'aura pas lieu.
Conditions à éviter: Les températures excessives qui causent l'évaporation. Protéger de la lumière.
Incompatibilités avec d'autres matériaux: Combustible fortes.
Produits de décomposition dangereux: Oxydes de carbone, l'oxydes de soufre et brome gazeux.

Section 11 L'Information Toxicologique

Toxicité aiguë: Données non disponibles
La corrosion de la peau et l'irritation: Données non disponibles
Des lésions oculaires graves / irritation: Données non disponibles
Respiratoire ou sensibilisation de la peau: Données non disponibles
Mutagénicité des cellules germinales: Données non disponibles
Cancérogène: Données non disponibles
NTP: Aucun composant de ce produit présent à des niveaux supérieurs ou égaux à 0,1% n'a été identifié comme cancérigène reconnu ou présumé par NTP.
IARC: Aucun composant de ce produit présent à des niveaux supérieurs ou égaux à 0,1% n'a été identifié comme cancérigène probable, possible ou confirmé par IARC.
OSHA: Aucun composant de ce produit présent à des niveaux supérieurs ou égaux à 0,1% n'a été identifié comme cancérigène ni comme cancérigène possible par OSHA.
Reproductive toxicity: Données non disponibles
STOT-exposition unique: Données non disponibles.
STOT-une exposition répétée: Données non disponibles
Risque d'aspiration: Données non disponibles
Effets d'une surexposition:
Inhalation: Peut être nocif en cas d'inhalation.
Ingestion: Peut être nocif en cas d'ingestion.
Peau: Peut causer une légère irritation.
Yeux: Peut causer une légère irritation.
Les signes et les symptômes de l'exposition: Pour le meilleur de notre connaissance les propriétés chimiques, physiques et toxicologiques n'ont pas été étudiées à fond.
Les données spécifiques n'est pas disponible. Exercice des procédures appropriées afin de minimiser les dangers potentiels.
Informations complémentaires: RTECS #: Données non disponibles

Section 12 L'Information Écologique

Toxicité pour les poissons: Pas de données disponible
Toxicité pour les daphnies et autres invertébrés aquatiques: Pas de données disponible
Toxicité pour les algues: Pas de données disponible
Persistance et dégradabilité: Pas de données disponible **Potentiel de bioaccumulation:** Pas de données disponible
Mobilité dans le sol: Pas de données disponibles **Évaluation PBT et vPvB:** Pas de données disponibles
Autres effets indésirables: Un danger pour l'environnement ne peut pas être exclu dans l'éventualité d'une manipulation ou d'élimination.

Section 13 Considérations De Disposition

Ces lignes directrices sont destinées à l'élimination de la disposition d'un catalogue de taille seules les quantités. Les règlements fédéraux peuvent s'appliquer aux contenants vides. Des réglementations nationales et / ou local peut être différent. Éliminer conformément à toutes les réglementations locales, provinciales et fédérales ou d'un contrat avec une agence élimination des produits chimiques sous licence..

Section 14 L'Information De Transport

Numéro UN / NA: Non applicable **Nom d'expédition:** Non réglé
Classe de danger: Non applicable **Groupe d'emballage:** Non applicable **Quantité à déclarer:** Non **Polluant marin:** Non
Exceptions: Non applicable **2012 ERG Guide #:** Non applicable

Section 15 L'Information De Normalisation

Un produit chimique est considéré comme inscrit si le numéro CAS pour la forme anhydre est sur la liste d'inventaire.

Composant	TSCA	CERLCA (RQ)	RCRA code	DSL	NDSL	Classification SIMDUT
Bleu de bromothymol, sel sodique	Listed	Not Listed	Not Listed	Listed	Not Listed	Produit non contrôlé

Section 16 L'Information Additionnelle

Les informations contenues dans ce document sont fournis sans garantie d'aucune sorte. Les employeurs devraient considérer cette information seulement comme complément à d'autres informations recueillies par eux et doivent prendre des décisions indépendantes de la pertinence et l'exhaustivité de l'information de toutes les sources afin d'assurer une utilisation correcte de ces matériaux et de la sécurité et la santé des employés. NTP: National Toxicology Program, IARC: International Agency for Research on Cancer, OSHA: Occupational Safety and Health Administration, STOT: Specific Target Organ Toxicity, SE: Single Exposure, RE: Repeated Exposure.

Section 1 L'information de produit chimique et de compagnie

Science Take-Out
80 Office Park Way
Pittsford, NY 14534
(585)764-5400

**CHEMTREC 24 Numéros De Téléphone
De Secours D'Heure (800) 424-9300**
Pour l'usage de laboratoire seulement.
Pas pour l'usage de drogue, de nourriture ou de ménage.

Produit	L'eau
Synonymes	"Arsenic Test Solution" (simulated)

Section 2 Identification De Risques

Cette substance ou un mélange n'a pas été classé à ce moment selon le Système général harmonisé (SGH) de classification et d'étiquetage des produits chimiques.

Mention d'avertissement: Non classé

Pictogrammes: Non classé

Les organes cibles: Aucun connu.

Classification par le GHS: Non classé

Renseignements sur l'étiquette GHS: Mention de danger(s): Non classé

Déclarations de précaution(s): Non classé

CA Prop 65 - Ce produit ne contient pas de produits chimiques connus à l'État de Californie pour causer le cancer, des malformations congénitales, ou toute autre atteinte à la reproduction.

Section 3 Composition / Information Sur Des Ingrédients

Nommé Chimique	CAS #	%	EINECS
L'eau	7732-18-5	100%	231-791-2

Section 4 Mesures De Premiers Soins

INGESTION:

INHALATION:

CONTACT AVEC LES YEUX:

ABSORPTION PAR LA PEAU:

Section 5 Mesures De Lutte Contre l'Incendie

Moyens d'extinction: Utilisez des supports adaptés pour éteindre le feu à l'appui.

Actions de protection pour les sapeurs-pompiers: En cas d'incendie, porter un appareil respiratoire NIOSH / MSHA approuvé autonome et un équipement complet de protection. Utiliser un jet d'eau pour maintenir incendie refroidir les conteneurs exposés.

Dangers spécifiques: Données non disponibles

Section 6 Mesures De Déchargement Accidentel

Précautions personnelles: Évacuer le personnel vers la zone sûre. Utiliser un équipement de protection personnelle comme indiqué dans la Section 8. Assurer une ventilation adéquate.

Précautions environnementales: Éviter tout ruissellement vers les égouts pluviaux et les fossés qui aboutissent aux voies navigables.

Confinement et de nettoyage: Absorbent avec le matériel sec inerte, balayez ou nettoyez à l'aspirateur vers le haut et placez dans un récipient approprié pour la disposition appropriée. Laver la zone de déversement avec du savon et de l'eau.

Section 7 Manipulation Et Stockage

Précautions pour la manutention en toute sécurité: Lire l'étiquette sur le contenant avant d'utiliser. Ne pas porter de lentilles cornéennes lorsque vous travaillez avec des produits chimiques. Tenir hors de portée des enfants. Éviter tout contact avec les yeux, la peau et les vêtements. Ne pas inhaler les vapeurs, les embruns ou le brouillard. Utiliser avec une ventilation adéquate. Éviter l'ingestion. Bien se laver après la manipulation. Retirer et laver les vêtements avant de les réutiliser.

Conditions de stockage: Stocker dans un endroit frais et bien aéré, loin des substances incompatibles.

Section 8 Commandes D'Exposition / Protection Personnelle

Limites d'exposition:	Nommé Chimique	ACGIH (TLV)	OSHA (PEL)	NIOSH (REL)
	L'eau	Aucun établi.	Aucun établi.	Aucun établi.

Contrôles d'ingénierie: Les installations d'entreposage ou d'utilisation de ce matériel doit être équipé d'une douche oculaire et une douche de sécurité et le matériel d'extinction d'incendie. Le personnel doit porter des lunettes de sécurité, des lunettes, ou un écran facial, une blouse de laboratoire ou tablier, des gants protecteurs appropriés. Utiliser une ventilation adéquate pour maintenir les concentrations atmosphériques faible.

Protection respiratoire: Aucun ne devrait être nécessaire dans le laboratoire normal manipulant aux températures ambiantes. Si les conditions brumeuses prévaloir, travailler dans la hotte ou de porter un masque respiratoire approuvé NIOSH / MSHA.

Section 9 Propriétés Physiques Et Chimiques

Apparence: Liquide, jaune clair. Odeur: Aucun odeur. Seuil de l'odeur: Données non disponibles. pH: Données non disponibles. Point de fusion / congélation: Environ 0°C (32°F) (eau) Point d'ébullition: 100oC @760mmHg (212oF) Point d'éclair: Non combustible	Taux d'évaporation (Eau = 1): 1 Inflammabilité (solide / gaz): Données non disponibles. Limites d'explosivité: Bas / Max: Données non disponibles Pression de vapeur (mm Hg): 14 Densité de vapeur (Air = 1): 0.7 Densité relative (gravité spécifique): 1.0 Solubilité (s): non applicable	Coefficient de partage: Données non disponibles Auto-inflammation: Données non disponibles Température de décomposition: Données non disponibles. Viscosité: Données non disponibles. Formule moléculaire: H ₂ O Poids moléculaire: 18.02
--	--	---

Section 10 Stabilité Et Réactivité

Stabilité chimique: Stable **Polymérisation dangereuse:** N'aura pas lieu.
Conditions à éviter: Données non disponibles.
Incompatibilités avec d'autres matériaux: Données non disponibles.
Produits de décomposition dangereux: Données non disponibles.

Section 11 L'Information Toxicologique

Toxicité aiguë: Oral-rat LD50: >90 mL/kg
La corrosion de la peau et l'irritation: Données non disponibles
Des lésions oculaires graves / irritation: Données non disponibles
Respiratoire ou sensibilisation de la peau: Données non disponibles
Mutagénicité des cellules germinales: Données non disponibles
Cancérogène: Données non disponibles
NTP: Aucun composant de ce produit présent à des niveaux supérieurs ou égaux à 0,1% n'a été identifié comme cancérogène reconnu ou présumé par NTP.
IARC: Aucun composant de ce produit présent à des niveaux supérieurs ou égaux à 0,1% n'a été identifié comme cancérogène probable, possible ou confirmé par IARC.
OSHA: Aucun composant de ce produit présent à des niveaux supérieurs ou égaux à 0,1% n'a été identifié comme cancérogène ni comme cancérogène possible par OSHA.
Reproductive toxicity: Données non disponibles
STOT-exposition unique: Données non disponibles.
STOT-une exposition répétée: Données non disponibles
Risque d'aspiration: Données non disponibles
Effets d'une surexposition:
Inhalation: Peut être nocif en cas d'inhalation.
Ingestion: Peut être nocif en cas d'ingestion.
Peau: Peut causer une légère irritation.
Yeux: Peut causer une légère irritation.
Les signes et les symptômes de l'exposition: Pour le meilleur de notre connaissance les propriétés chimiques, physiques et toxicologiques n'ont pas été étudiées à fond.
Les données spécifiques n'est pas disponible. Exercice des procédures appropriées afin de minimiser les dangers potentiels.
Informations complémentaires: RTECS #: Données non disponibles

Section 12 L'Information Écologique

Toxicité pour les poissons: Pas de données disponible
Toxicité pour les daphnies et autres invertébrés aquatiques: Pas de données disponible
Toxicité pour les algues: Pas de données disponible
Persistence et dégradabilité: Pas de données disponible **Potentiel de bioaccumulation:** Pas de données disponible
Mobilité dans le sol: Pas de données disponibles **Évaluation PBT et vPvB:** Pas de données disponibles
Autres effets indésirables: Un danger pour l'environnement ne peut pas être exclu dans l'éventualité d'une manipulation ou d'élimination.

Section 13 Considérations De Disposition

Ces lignes directrices sont destinées à l'élimination de la disposition d'un catalogue de taille seules les quantités. Les règlements fédéraux peuvent s'appliquer aux contenants vides. Des réglementations nationales et / ou local peut être différent. Éliminer conformément à toutes les réglementations locales, provinciales et fédérales ou d'un contrat avec une agence élimination des produits chimiques sous licence..

Section 14 L'Information De Transport

Numéro UN / NA: Non applicable **Nom d'expédition:** Non réglé
Classe de danger: Non applicable **Groupe d'emballage:** Non applicable **Quantité à déclarer:** Non **Polluant marin:** Non
Exceptions: Non applicable **2012 ERG Guide #:** Non applicable

Section 15 L'Information De Normalisation

Un produit chimique est considéré comme inscrit si le numéro CAS pour la forme anhydre est sur la liste d'inventaire.

Composant	TSCA	CERLCA (RQ)	RCRA code	DSL	NDSL	Classification SIMDUT
L'eau	Listed	Not Listed	Not Listed	Listed	Not Listed	Produit non contrôlé

Section 16 L'Information Additionnelle

Les informations contenues dans ce document sont fournis sans garantie d'aucune sorte. Les employeurs devraient considérer cette information seulement comme complément à d'autres informations recueillies par eux et doivent prendre des décisions indépendantes de la pertinence et l'exhaustivité de l'information de toutes les sources afin d'assurer une utilisation correcte de ces matériaux et de la sécurité et la santé des employés. NTP: National Toxicology Program, IARC: International Agency for Research on Cancer, OSHA: Occupational Safety and Health Administration, STOT: Specific Target Organ Toxicity, SE: Single Exposure, RE: Repeated Exposure.

Section 1 L'information de produit chimique et de compagnie

Science Take-Out
80 Office Park Way
Pittsford, NY 14534
(585)764-5400

**CHEMTREC 24 Numéros De Téléphone
De Secours D'Heure (800) 424-9300**
Pour l'usage de laboratoire seulement.
Pas pour l'usage de drogue, de nourriture ou de ménage.

Produit	Rouge de méthyl, 0.05% solute
Synonymes	"Coliform Bacteria Test Solution" (simulated)

Section 2 Identification De Risques

Cette substance ou un mélange n'a pas été classé à ce moment selon le Système général harmonisé (SGH) de classification et d'étiquetage des produits chimiques.

Mention d'avertissement: Non classé

Pictogrammes: Non classé

Les organes cibles: Aucun connu.

Classification par le GHS: Non classé

Renseignements sur l'étiquette GHS: Mention de danger(s): Non classé

Déclarations de précaution(s): Non classé

Des renseignements supplémentaires:

Ne pas respirer les vapeurs, les embruns ou le brouillard. Éviter tout contact avec les yeux, la peau ou les vêtements. Porter des gants de protection / des vêtements de protection / un équipement de protection des yeux / du visage. Se laver les mains soigneusement après manipulation. Consulter un médecin en cas de malaise.

CA Prop 65 - Ce produit ne contient pas de produits chimiques connus à l'État de Californie pour causer le cancer, des malformations congénitales, ou toute autre atteinte à la reproduction.

Section 3 Composition / Information Sur Des Ingrédients

Nommé Chimique	CAS #	%	EINECS
L'eau	7732-18-5	99.95%	231-791-2
Rouge de méthyle, sel sodique	845-10-3	0.05%	212-682-9

Section 4 Mesures De Premiers Soins

INGESTION: PEUT ÊTRE NOCIF EN CAS D'INGESTION. Appeler un médecin ou un centre antipoison immédiatement. Provoquer le vomissement seulement si elle est informée par le personnel compétent médicaux. Ne jamais rien donner par la bouche à une personne inconsciente.

INHALATION: PEUT ÊTRE NOCIF EN CAS D'INHALATION. PEUT IRRITER LES VOIES RESPIRATOIRES. Sortir au grand air. Si elle ne respire pas, pratiquer la respiration artificielle. Si la respiration est difficile, donner de l'oxygène. Obtenir des soins médicaux.

CONTACT AVEC LES YEUX: PEUT CAUSER L'IRRITATION DES YEUX. Vérifier et enlever les lentilles de contact. Rincer abondamment à l'eau pendant au moins 15 minutes, en soulevant les paupières inférieures et supérieures de temps en temps. Obtenez une attention médicale immédiate.

ABSORPTION PAR LA PEAU: PEUT CAUSER UNE IRRITATION DE LA PEAU. Enlever les vêtements contaminés. Rincer soigneusement avec du savon doux et d'eau. En cas d'irritation, consulter un médecin.

Section 5 Mesures De Lutte Contre l'Incendie

Moyens d'extinction: Dioxyde de carbone, produit chimique sec, du sable sec, mousse anti-alcool.

Actions de protection pour les sapeurs-pompiers: En cas d'incendie, porter un appareil respiratoire NIOSH / MSHA approuvé autonome et un équipement complet de protection. Utiliser un jet d'eau pour maintenir incendie refroidir les conteneurs exposés.

Dangers spécifiques: En cas de feu, de l'eau peut s'évaporer à partir de cette solution, qui peut causer les produits dangereux de décomposition à être formée comme poussière ou vapeur.

Section 6 Mesures De Déchargement Accidentel

Précautions personnelles: Évacuer le personnel vers la zone sûre. Utiliser un équipement de protection personnelle comme indiqué dans la Section 8. Assurer une ventilation adéquate.

Précautions environnementales: Éviter tout ruissellement vers les égouts pluviaux et les fossés qui aboutissent aux voies navigables.

Confinement et de nettoyage: Absorbent avec le matériel sec inerte, balayez ou nettoyez à l'aspirateur vers le haut et placez dans un récipient approprié pour la disposition appropriée. Laver la zone de déversement avec du savon et de l'eau.

Section 7 Manipulation Et Stockage

Précautions pour la manutention en toute sécurité: Lire l'étiquette sur le contenant avant d'utiliser. Ne pas porter de lentilles cornéennes lorsque vous travaillez avec des produits chimiques. Tenir hors de portée des enfants. Éviter tout contact avec les yeux, la peau et les vêtements. Ne pas inhaler les vapeurs, les embruns ou le brouillard. Utiliser avec une ventilation adéquate. Éviter l'ingestion. Bien se laver après la manipulation. Retirer et laver les vêtements avant de les réutiliser.

Conditions de stockage: Stocker dans un endroit frais et bien aéré, loin des substances incompatibles.

Section 8 Commandes D'Exposition / Protection Personnelle

Limites d'exposition:	Nommé Chimique	ACGIH (TLV)	OSHA (PEL)	NIOSH (REL)
	Rouge de méthyle	Aucun établi.	Aucun établi.	Aucun établi.

Contrôles d'ingénierie: Les installations d'entreposage ou d'utilisation de ce matériel doit être équipé d'une douche oculaire et une douche de sécurité et le matériel d'extinction d'incendie. Le personnel doit porter des lunettes de sécurité, des lunettes, ou un écran facial, une blouse de laboratoire ou tablier, des gants protecteurs appropriés. Utiliser une ventilation adéquate pour maintenir les concentrations atmosphériques faible.

Protection respiratoire: Aucun ne devrait être nécessaire dans le laboratoire normal manipulant aux températures ambiantes. Si les conditions brumeuses prévaloir, travailler dans la hotte ou de porter un masque respiratoire approuvé NIOSH / MSHA.

Section 9 Propriétés Physiques Et Chimiques

Apparence: Liquide, clair, incolore rouge. Odeur: Aucun odeur. Seuil de l'odeur: Données non disponibles. pH: Données non disponibles. Point de fusion / congélation: Environ 0°C (32°F) (eau) Point d'ébullition: Environ 100°C (212°F) (eau) Point d'éclair: Données non disponibles	Taux d'évaporation (Eau = 1): <1 Inflammabilité (solide / gaz): Données non disponibles. Limites d'explosivité: Bas / Max: Données non disponibles Pression de vapeur (mm Hg): 14 (eau) Densité de vapeur (Air = 1): 0.7 (eau) Densité relative (gravité spécifique): Environ 1.0 (eau) Solubilité (s): Complet dans l'eau.	Coefficient de partage: Données non disponibles Auto-inflammation: Données non disponibles Température de décomposition: Données non disponibles. Viscosité: Données non disponibles. Formule moléculaire: Mélange Poids moléculaire: Mélange
---	--	--

Section 10 Stabilité Et Réactivité

Stabilité chimique: Stable **Polymérisation dangereuse:** N'aura pas lieu.
Conditions à éviter: Les températures excessives qui causent l'évaporation. Protéger de la lumière.
Incompatibilités avec d'autres matériaux: Combustible fortes, l'agents reducteurs.
Produits de décomposition dangereux: Oxydes de carbone, l'oxydes d'azote et l'oxydes de sodium.

Section 11 L'Information Toxicologique

Toxicité aiguë: Oral-rat TDLo: 12000 mg/kg [Rouge de méthyle]
La corrosion de la peau et l'irritation: Données non disponibles
Des lésions oculaires graves / irritation: Données non disponibles
Respiratoire ou sensibilisation de la peau: Données non disponibles
Mutagénicité des cellules germinales: Données non disponibles
Cancérogène: Données non disponibles
NTP: Aucun composant de ce produit présent à des niveaux supérieurs ou égaux à 0,1% n'a été identifié comme cancérigène reconnu ou présumé par NTP.
IARC: Aucun composant de ce produit présent à des niveaux supérieurs ou égaux à 0,1% n'a été identifié comme cancérigène probable, possible ou confirmé par IARC.
OSHA: Aucun composant de ce produit présent à des niveaux supérieurs ou égaux à 0,1% n'a été identifié comme cancérigène ni comme cancérigène possible par OSHA.
Reproductive toxicity: Données non disponibles
STOT-exposition unique: Données non disponibles.
STOT-une exposition répétée: Données non disponibles
Risque d'aspiration: Données non disponibles
Effets d'une surexposition:
Inhalation: Peut être nocif en cas d'inhalation.
Ingestion: Peut être nocif en cas d'ingestion.
Peau: Peut causer une légère irritation.
Yeux: Peut causer une légère irritation.
Les signes et les symptômes de l'exposition: Pour le meilleur de notre connaissance les propriétés chimiques, physiques et toxicologiques n'ont pas été étudiées à fond.
Les données spécifiques n'est pas disponible. Exercice des procédures appropriées afin de minimiser les dangers potentiels.
Informations complémentaires: RTECS #: DG8960000 [Rouge de méthyle]

Section 12 L'Information Écologique

Toxicité pour les poissons: Pas de données disponible
Toxicité pour les daphnies et autres invertébrés aquatiques: Pas de données disponible
Toxicité pour les algues: Pas de données disponible
Persistance et dégradabilité: Pas de données disponible **Potentiel de bioaccumulation:** Pas de données disponible
Mobilité dans le sol: Pas de données disponibles **Évaluation PBT et vPvB:** Pas de données disponibles
Autres effets indésirables: Un danger pour l'environnement ne peut pas être exclu dans l'éventualité d'une manipulation ou d'élimination.

Section 13 Considérations De Disposition

Ces lignes directrices sont destinées à l'élimination de la disposition d'un catalogue de taille seules les quantités. Les règlements fédéraux peuvent s'appliquer aux contenants vides. Des réglementations nationales et / ou local peut être différent. Éliminer conformément à toutes les réglementations locales, provinciales et fédérales ou d'un contrat avec une agence élimination des produits chimiques sous licence..

Section 14 L'Information De Transport

Numéro UN / NA: Non applicable **Nom d'expédition:** Non réglé
Classe de danger: Non applicable **Groupe d'emballage:** Non applicable **Quantité à déclarer:** Non **Polluant marin:** Non
Exceptions: Non applicable **2012 ERG Guide #:** Non applicable

Section 15 L'Information De Normalisation

Un produit chimique est considéré comme inscrit si le numéro CAS pour la forme anhydre est sur la liste d'inventaire.

Composant	TSCA	CERLCA (RQ)	RCRA code	DSL	NDSL	Classification SIMDUT
Rouge de méthyle, sel sodique	Listed	Not Listed	Not Listed	Listed	Not Listed	Not Listed

Section 16 L'Information Additionnelle

Les informations contenues dans ce document sont fournis sans garantie d'aucune sorte. Les employeurs devraient considérer cette information seulement comme complément à d'autres informations recueillies par eux et doivent prendre des décisions indépendantes de la pertinence et l'exhaustivité de l'information de toutes les sources afin d'assurer une utilisation correcte de ces matériaux et de la sécurité et la santé des employés. NTP: National Toxicology Program, IARC: International Agency for Research on Cancer, OSHA: Occupational Safety and Health Administration, STOT: Specific Target Organ Toxicity, SE: Single Exposure, RE: Repeated Exposure.