



Anti-c-myc IgG conjugated to Dylight® 650

Product Number D10-1714
Amount 100 µg total protein or 1 mg total protein
Clone 9E10 (sequence EQKLISEEDL)
Lot Number DY650012-12-001
Conc 0.909 mg/ml
Store at 2-8°C

Form/Shipping & Storage

Supplied liquid. Upon receipt, store at 2-8 °C. Fluorescent dyes are sensitive to light. Please store in the dark and avoid exposure to sunlight.

Handling

We recommend that the investigator determine the appropriate working concentration for their specific application.

Buffer

0.1M sodium phosphate, 0.15M NaCl, pH7.4, 1% BSA, 0.02% NaN₃

Stability


Product should be stored at 2-8°C in the dark and be used within 1 year. If further dilution of the conjugate is required, use diluted material within one month.

Note

For research use only, not for diagnostic or therapeutic use.

Spectral Characteristic

Visible absorption maxima 652 nm
Emission maximum 672 nm

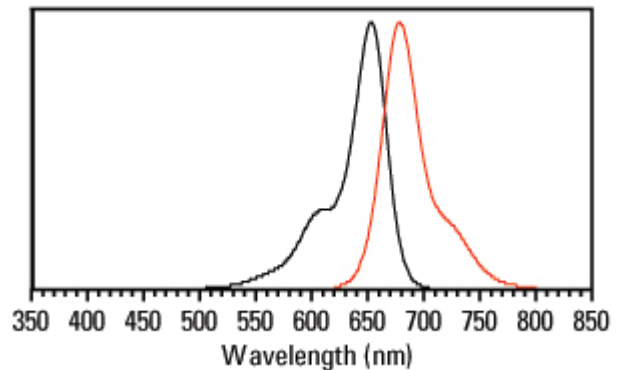

References:

Gazitt Y, He YJ, Erdos GW, Chang L, Ashktorab H, Cohen RJ. Development of a two color immunofluorescence stain and immunolocalization method for N-myc and c-myc oncoproteins with a newly generated mouse IgM anti N-myc antibody. J Immunol Methods. 1992 Apr 8;148(1-2):159-69.

Kieke MC, Cho BK, Boder ET, Kranz DM, Wittrup KD. Isolation of anti-T cell receptor scFv mutants by yeast surface display. Protein Eng. 1997 Nov;10(11):1303-10.

Lincoln ST, Bauer KD. Limitations in the measurement of c-myc oncoprotein and other nuclear antigens by flow cytometry. Cytometry. 1989 Jul;10(4):456-62.



COLUMBIA BIOSCIENCES