

# CENTRA

CENTRALIZED WATER




## CENTRA-RDS

The CENTRA has revolutionised the way that large volumes of pure water are produced, stored and distributed. Instead of a central lab system, comprising of an untidy grouping of components and pipework, you can now have just one integrated system box. The CENTRA-RDS is a storage and recirculating distribution system designed to work with other CENTRA or PURELAB systems or existing purification installations. In addition to storing and distributing pure water, the CENTRA-RDS maintains and improves water purity.

- The compact design provides more flexible installation options for new buildings and refurbishments. The CENTRA's small footprint means that it can be placed closer to the laboratory, avoiding the negative cost and design implications of long pipe work loops
- Reliable, continual supply of pure water using unique access controls, leak detection systems and full AV alarms with optional building management system (BMS) connectivity
- Optimized inorganic water quality through use of in-line purification technologies. Recirculated water is UV treated, filtered and (if fitted) can be improved through deionization
- 350 l reservoir of purified water available at up to 30 l/min from a distribution loop. A wide range of water purities are possible from RO permeate up to 18.2 MΩ-cm Type I purity
- Low microbial counts achieved through vent filtration, spray-ball water distribution and smooth internal surface in the reservoir combined with UV oxidation and 0.2 µm filtration in the loop

### Process Flow CENTRA-RDS

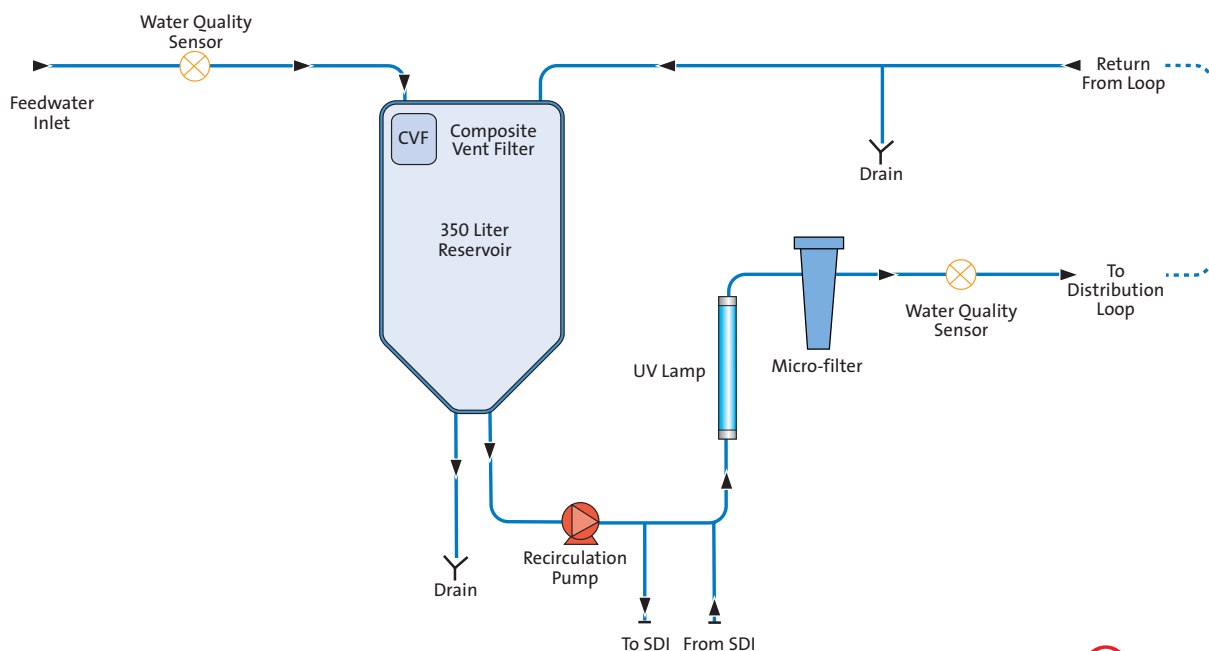

# CENTRA

## CENTRALIZED WATER

### Treated Water Specifications

Model	Voltage/ Freq	RDS	RDS HFV/HFR	
Flowrate (exit from unit) - l/min (USG/min)	230/50	16 (4.2)	30 (8)	
	230/50/60	50	n/a	30 (8)
		60	n/a	38 (10)
	115/60	18 (4.8)	38 (10)	
Outlet pressure (exit from unit) - bar (psi)		3 (45)	3 (45)	
Inorganics (resistivity @ 25°C)		Up to 18.2MΩ-cm <sup>1</sup>	Up to 18.2MΩ-cm <sup>1</sup>	
Bacteria - typical		<5 CFU/ml <sup>2</sup>	<5 CFU/ml <sup>2</sup>	
Particles		0.2 µm filtration	0.2 µm filtration	

<sup>1</sup> Ion-exchange cylinder (SDI) installed

<sup>2</sup> Dependent on feedwater quality and system being correctly installed and regularly sanitized.

### Dimensions and weights

Dimensions	Height 1820mm (71.7in), Width 730mm (28.75in), Depth 890mm (35in),
Shipping weight	180kg (396lb)
Operational weight	530kg (1,168lb)

### Feedwater Requirements

Source	Potable tap water, preferably reverse osmosis or de-ionized	
Fouling Index (maximum)	5	5
Feedwater flow – l/min (maximum)	20	40
Feedwater pressure	10 bar (150 psi) maximum, 0.1 bar (1.5 psi) minimum	
Drain requirements l/min (gravity fall with air gap)	20	70

### Electrical Requirements

Mains input	230V ac, 50Hz, 115V ac, 60Hz
System control voltage	24V dc
Power consumption during recirculation	1200VA
Electrical protection rating	20 amps
Reservoir level connection for control of additional ELGA products	Jack Plug 3.5mm DIN plug 6 way
Noise level during recirculation	<65dBA

#### ELGA LabWater

Tel: +44 (0) 1494 887500 Fax: +44 (0) 1494 887505

Email: [info@elgalabwater.com](mailto:info@elgalabwater.com) Website: [www.elgalabwater.com](http://www.elgalabwater.com)

ELGA is the global laboratory water brand name of Veolia Water Solutions and Technologies. The information contained in this document is the property of VWS (UK) Ltd, trading as ELGA LabWater, and is supplied without liability for errors or omissions.

© VWS (UK) Ltd. 2010 – All rights reserved. ELGA® and CENTRA® are registered trademarks of VWS (UK) Ltd.

LITR 38762-03