


VWR Irradiated Non-Skid Shoe Covers

- Gamma Irradiated to a Sterility Assurance Level (SAL) of 10^{-4}
- Non-Skid Sole Provides Maximum Traction on Wet and Smooth Surfaces
- Superior Coefficient of Friction
- Impervious Material: >99.9% Bacterial Filtration Efficiency
- Clean, Low-Linting Fabric

VWR Irradiated Non-Skid Shoe Covers have been engineered to optimize skid resistance, barrier protection, cleanliness, and comfort. These shoe covers are gamma irradiated to address the concern of bio-burden contamination. Manufactured and packaged in an ISO Class 8 (FED-STD-209E Class 100,000/M6.5) cleanroom, shoe covers are bulk packed in sealed polyethylene bags contained within a larger carton liner. Each case is supplied with a certificate of irradiation.

Shoe covers are rigorously tested and manufactured in an ISO Certified facility under stringent process controls to ensure that each product meets exacting quality standards and performs to specification. Our products are validated through independent lab testing.

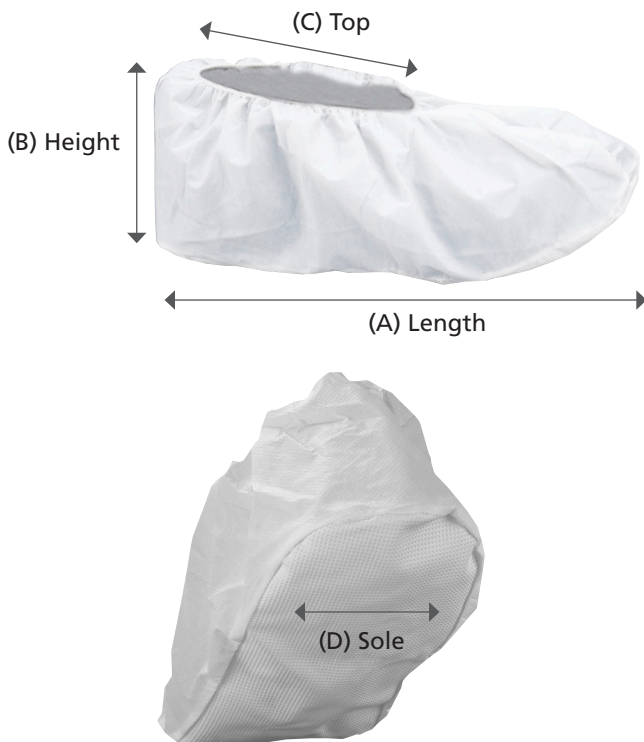

Shoe Cover Dimensions

Size	M	L	Universal	XL
Length (A)	11½"	12¾"	14"	15"
Height (B)	5"	5"	5"	5¾"
Top, Relaxed (C)	4"	4¾"	5"	5"
Top, Stretched (C)	8¾"	9¼"	10½"	11¾"
Sole Width (D)	4¼"	4½"	4½"	4¾"

VWR Irradiated Non-Skid Shoe Covers

Size	Cat. No.	Case of
White		
Medium	414004-478	200
Large	414004-479	200
Universal	414004-476	200
X-Large	414004-477	200

To order, visit vwr.com or call your VWR Sales Representative today.

Material Properties for VWR® Irradiated Non-Skid Shoe Covers

	Test Item	Result	Test Standard
PHYSICAL PROPERTIES	Particle Shedding (Helmke Drum)	Level I	IEST-RP-CC003.3
	Weight (g/m ²)	58.58	ASTM D3776
	Wyzenbeek Abrasion Mass Loss	6.17%	200 Rubs Using Zero Grade Emery Paper
	Mullen Bursting Strength (Avg. PSI)	18.0	ASTM D3786
	Abrasion Resistance of Shoe Cover Sole (Avg. cycles at which fabric shows noticeable wear)	1000	ASTM D3884 (H-18 Wheel, 500g Load)
	Tensile Strength (Avg. lbs./in.)		
	Warp	24.4	ASTM D5034
	Filling	20.3	
	Trapezoid Tearing Strength (Avg. lbs./in.)		
	Lengthwise Yarns	7.7	ASTM D4533
	Widthwise Yarns	13.7	
BARRIER PROPERTIES	Bacterial Filtration Efficiency (28.3L/min, 1CFM)	>99.9%	ASTM F2101
	Water Resistance	Pass	AATCC Method 42
SKID-RESISTANCE PROPERTIES	Coefficient of Friction (COF)		
	Kinetic COF	3.54	ASTM D1894
	Static COF	2.88	

All test results provided by independent third-party testing laboratories located in USA.

WARNING: These garments and associated materials are not suitable for use in some environments containing chemicals and/or hazardous agents. It is the responsibility of the user to determine the level of risk in a particular environment and the proper personal protection equipment needed. Garments manufactured from synthetic non-woven material may generate static electricity. Garments that contain an anti-static treatment are not intended to be used as a safety feature. These garments are not recommended to be used in a flammable and/or explosive environment. Contact VWR International for garment/fabric safety data. The application of these products is out of VWR International's control. Therefore, VWR International, LLC makes no warranties, expressed or implied, and assumes no liability as to the performance of these products for a particular use. Caution: Avoid heat and/or open flame.



Prices, product appearance and specifications are current at the time of printing, subject to change without notice. Availability for certain products may be limited by federal, state, provincial or local licensing requirements. VWR makes no claims or warranties concerning sustainable/green products. Any claims concerning sustainable/green products are the sole claims of the manufacturer and not those of VWR International, LLC. All prices are in U.S. dollars unless otherwise noted. Offers valid in USA and Canada, void where prohibited by law or company policy, while supplies last. Visit vwr.com to view our privacy policy and additional disclaimers.

VWR, forms of VWR, and the VWR logo and/or design are either registered trademarks ®, trademarks ™, or service marks ™ of VWR International, LLC in the United States and/or other countries. All other marks referenced herein are registered trademarks, trademarks or service marks of their respective owner(s). For a complete list of trademark owners please visit vwr.com.

©2010 VWR International, LLC. All rights reserved.

1010 REV2